



# **Pharmacy Information System (PhIS) and Clinic Pharmacy System (CPS)**

---

## **Pharmacy Based User Manual Radiopharmaceuticals**

<b>Version</b>	<b>: 13<sup>th</sup> Edition</b>
<b>Document ID</b>	<b>: PB_U.MANUAL_RADIOPHARMACEUTICALS</b>



© 2011-2023 Pharmacy Information System & Clinic Pharmacy System (PhIS& CPS) Project

**CONFIDENTIAL COPYRIGHTED MATERIAL** – *The information includes all concepts, comments, recommendations, and material, contained herein shall remain the property of Pharmacy Information System & Clinic Pharmacy System (PhIS & CPS) Project. No portion of this document shall be disclosed, duplicated or used in whole or in part of any purpose other than the purpose of the Pharmacy Information System & Clinic Pharmacy System (PhIS & CPS) Project execution only.*

*Reference ID : PB\_U. MANUAL\_RADIOPHARMACEUTICALS-13<sup>th</sup> E*

*Application reference: PhIS & CPS v2.6*



## Table of Contents

1.0	Introduction .....	1
1.1	Overview of PhIS.....	1
1.2	Purpose and Objectives .....	1
1.3	Organised Sections .....	1
2.0	Application Standard Features.....	2
2.1	PhIS Legend.....	2
3.0	Radiopharmaceuticals .....	4
	Overview .....	4
	User Group.....	4
	Functional Diagram .....	4
	Functional Description.....	4
3.1	Radiopharmaceuticals.....	5
3.1.1	Create New Appointment .....	5
3.2	Search Registered Patient.....	12
3.3	Radiopharmaceuticals.....	15
3.4	Radiopharmaceuticals Verification .....	19
3.4.1	Intervention.....	21
3.4.2	Hold/Resume Order .....	24
3.5	Radiopharmaceutical Preparation .....	28
3.6	Radiopharmaceutical Dispensing.....	38
3.7	Release Batch .....	41
3.8	Radiopharmaceutical Manufacturing .....	44
3.8.1	Generate Elution .....	44
3.8.2	KIT Preparation .....	54
3.8.3	Stock Verification.....	68
4.0	Acronyms .....	72
5.0	Links to Inventory Modules .....	72



## 1.0 Introduction

### 1.1 Overview of PhIS

Pharmacy Information System or better known as PhIS is a complete and comprehensive system that integrates pharmacy related services that geared toward pharmacy excellence. This implementation would transform most of the current manual process to electronic system would benefit facility end user in the health care sector.

There are 12 modules to assist service delivery by the health care sector which comprises of:

1. Order Management
2. Inpatient Pharmacy
3. Outpatient Pharmacy
4. Medication Counselling
5. Ward Pharmacy
6. Pharmacy Inventory
7. Manufacturing of Cytotoxic Drug Reconstitution (CDR), Parenteral Nutrition (PN), IV Admixture & Eye Drop, Radiopharmaceuticals and Extemporaneous
8. Adverse Drug Reaction & Drug Allergic Card (ADR & DAC)
9. Clinical Pharmacokinetics Services (TDM)
10. Drug Information & Consumer Education (DICE)
11. Medication Therapy Adherence Clinic (MTAC)
12. Data Mining (PhARM)

### 1.2 Purpose and Objectives

This user manual outlines the Radiopharmaceuticals module and its key features and functionalities. The primary objective is to guide users through the process of completing PhIS application process.

User will understand the following activities in details:

- Radiopharmaceuticals Appointment
- Radiopharmaceuticals Oder
- Radiopharmaceuticals Dispensing
- Radiopharmaceuticals Manufacturing









### 1.3 Organised Sections




These are the sections within this document:

- Section 1 : Introduction
- Section 2 : Application Standard Features
- Section 3 : Radiopharmaceuticals
- Section 4 : Acronyms
- Section 5 : Links to Clinical Modules

## 2.0 Application Standard Features

### 2.1 PhIS Legend

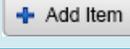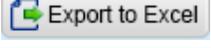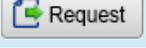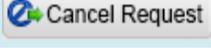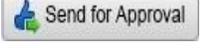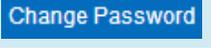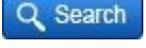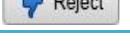Standard Legend			
	Login to PhIS		Logout from PhIS
	Reset Login Screen		Expand Menu
	Display Home Tab		Expand Module
	Close All Open Tabs		Refresh Screen
	Add/Create New Record		Show Help
*	Mandatory Field		Calendar Icon
	Close Window		Radio Button
	Edit Record		Cancel
	Save		Add Item to the list
	Export and Open Report in Excel Format		Automatically Display/Retrieve Box
	Request for Approval		Cancel the Request
	Send for Approval		Approve Transaction
	Change Login Password		Collapse Menu
	Collapse Module		Search Record
	Print		Search Icon
<input type="checkbox"/>	Checkbox		Delete Record
	Delete Item from the list		Empty Text Box
	Dropdown Box		Reject Transaction

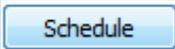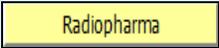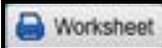Module Legend			
	Schedule button		Confirmation button
	Consent form		Verification
	Preparation		Dispensing
	Submission Approval		Printing/Editing Label
	Radiopharmaceuticals		Calculation Activity Radiopharmaceuticals
	Generate the worksheet		Release Batch

**Note**

To learn more about Login Information, kindly click [Login Information](#) Modules for descriptive step.

## 3.0 Radiopharmaceuticals

### Overview

This module will be used by the user at the

### User Group

This module is intended for Doctor, Specialist and Pharmacist (subject to user assigned by the facility).

### Functional Diagram

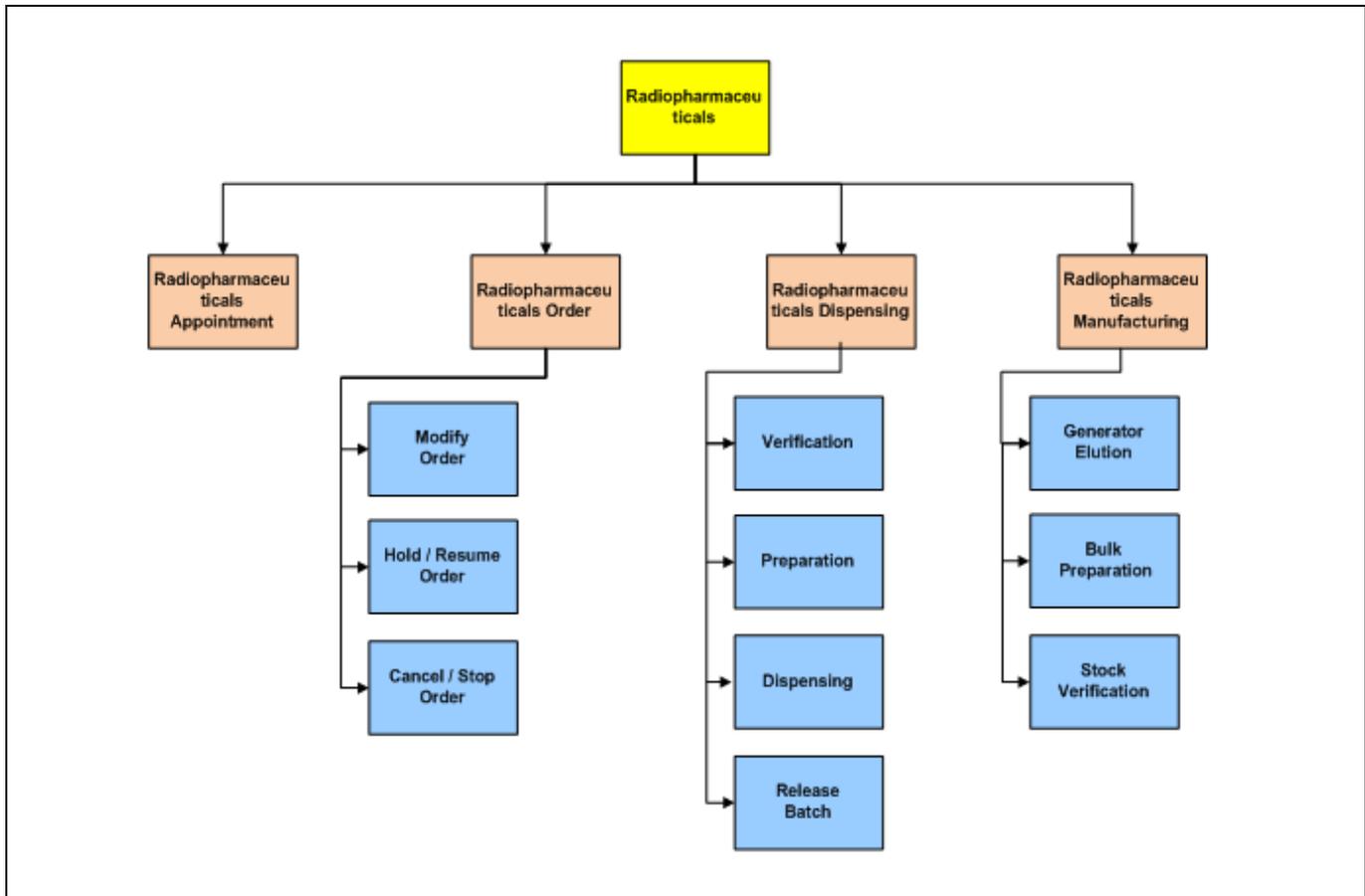


Figure 3.0-1

### Functional Description


Slow Moving comprises of three (4) main functions:

- **Radiopharmaceuticals Appointment**  
This function allows user to create an Appointment for the patient
- **Radiopharmaceuticals Order**  
This function is used for user to create an order for Radiopharmaceutical
- **Radiopharmaceuticals Dispensing**  
This function allows user to dispense Radiopharmaceutical order
- **Radiopharmaceuticals Manufacturing**

## 3.1 Radiopharmaceuticals

### 3.1.1 Create New Appointment

Radiopharmaceuticals Appointment screen will allow user to create appointment for the purpose of scanning and therapy for patient.

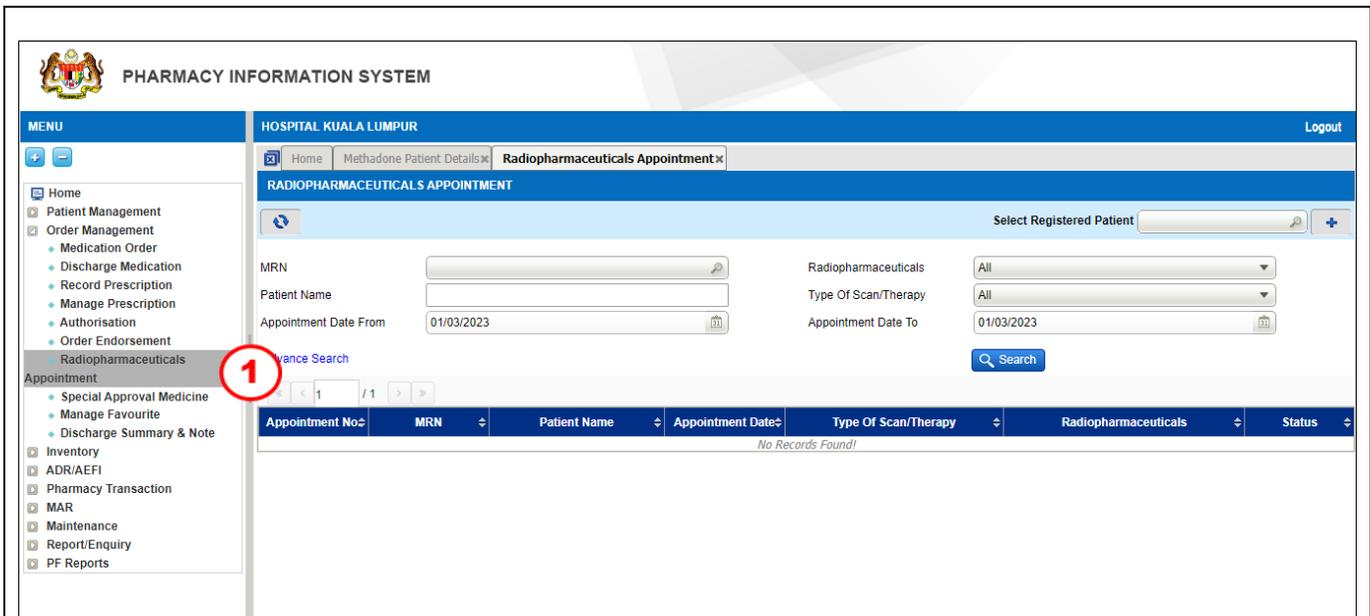

The screenshot displays the 'RADIOPHARMACEUTICALS APPOINTMENT' screen. The left sidebar menu has 'Radiopharmaceuticals' selected, with 'Appointment' highlighted. A red circle with the number '1' is placed over the 'Advanced Search' link. The main content area includes a search form with the following fields: MRN, Patient Name, Appointment Date From (01/03/2023), Appointment Date To (01/03/2023), Radiopharmaceuticals (All), and Type Of Scan/Therapy (All). A 'Search' button is located at the bottom right of the form. Below the form is a table with columns: Appointment No., MRN, Patient Name, Appointment Date, Type Of Scan/Therapy, Radiopharmaceuticals, and Status. The table is currently empty, showing 'No Records Found!'.


Figure 3.1-1 Radiopharmaceuticals Appointment Listing Page

#### STEP 1

Click on 'Order Management' and select 'Radiopharmaceuticals Appointment' sub-menu as per Figure 3.1-1

#### Note

- User have the options to view existing appointment records by clicking the check box :-
  - a) Show All
- Various search criteria is provided as below:
  - a) MRN
  - b) Radiopharmaceuticals
  - c) Patient Name
  - d) Appointment no
  - e) Status
  - f) Appointment Date From
  - g) Appointment Date To
- By default list of existing appointments will be displayed based on current system date
- Click on the [Advanced Search](#) hyperlink for advance search.

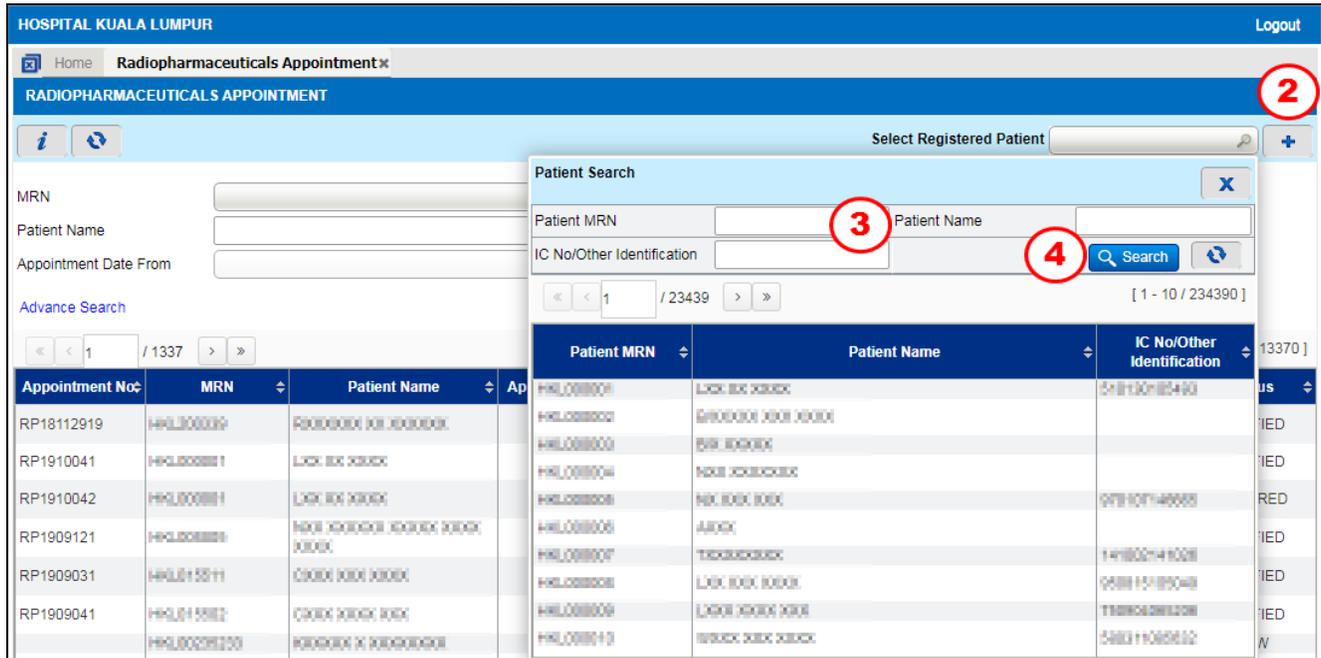

**HOSPITAL KUALA LUMPUR** Logout

Home **Radiopharmaceuticals Appointment**

**RADIOPHARMACEUTICALS APPOINTMENT**

Select Registered Patient

MRN

Patient Name

Appointment Date From

Advance Search

**Patient Search**

Patient MRN  Patient Name

IC No/Other Identification

< < 1 / 23439 > > [ 1 - 10 / 234390 ]

Appointment No	MRN	Patient Name	Ap	Patient MRN	Patient Name	IC No/Other Identification
RP18112919	HKL000039	RODODODI KE JODODOK		HKL000001	LAKI BAK BAKK	5-11-00-1854-00
RP1910041	HKL000081	LAKI BAK BAKK		HKL000002	PEREMBAHAN JODI JODOK	
RP1910042	HKL000081	LAKI BAK BAKK		HKL000003	PEREMBAHAN	
RP1909121	HKL000089	PEREMBAHAN JODOK JODOK JODOK		HKL000004	PEREMBAHAN	
RP1909031	HKL015511	PEREMBAHAN JODOK JODOK		HKL000005	PEREMBAHAN	978107146000
RP1909041	HKL015512	PEREMBAHAN JODOK JODOK		HKL000006	PEREMBAHAN	1-1-002141000
	HKL00208100	PEREMBAHAN JODOK JODOK		HKL000007	PEREMBAHAN	050815185040
				HKL000008	LAKI BAK BAKK	110000000000
				HKL000009	LAKI BAK BAKK	110000000000
				HKL000010	PEREMBAHAN JODOK JODOK	5-00-11-000010

**Figure 3.1.1-2 Radiopharmaceuticals Appointment**

**STEP 2**


Click on the  button to search for registered patient and system will be displayed as per Figure 3.1.1-2

**STEP 3**

Enter patient name at 'Patient name' field

**STEP 4**

Click on  button to search patient name and system will be displayed as per Figure 3.1.1-2

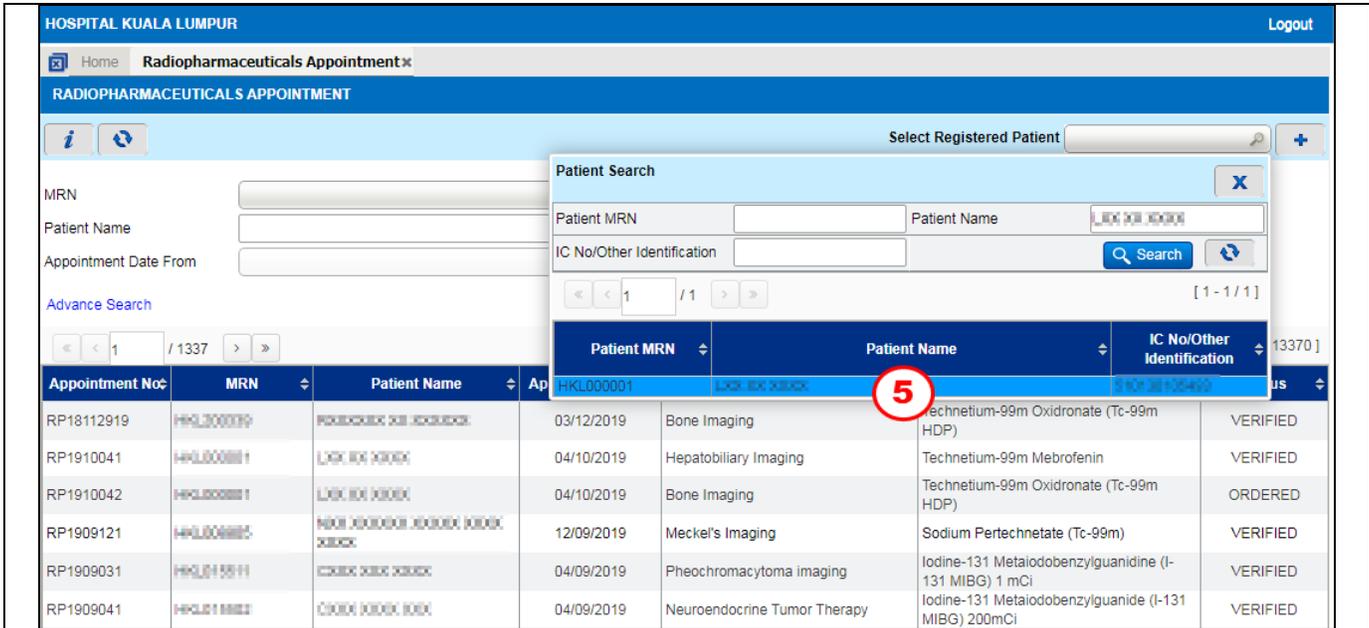

HOSPITAL KUALA LUMPUR Logout

Home Radiopharmaceuticals Appointment x

RADIOPHARMACEUTICALS APPOINTMENT

Select Registered Patient

MRN

Patient Name

Appointment Date From

Advance Search

1 / 1337

Patient Search

Patient MRN

Patient Name

IC No/Other Identification

Search


1 / 1

Appointment No	MRN	Patient Name	Ap	IC No/Other Identification	13370	
RP18112919	HKL000009	XXXXXXXXXX XXXXXXXX	03/12/2019	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	VERIFIED
RP1910041	HKL000011	LAK SA XXXX	04/10/2019	Hepatobiliary Imaging	Technetium-99m Mebrofenin	VERIFIED
RP1910042	HKL000011	LAK SA XXXX	04/10/2019	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	ORDERED
RP1909121	HKL000015	NAH XXXXXX XXXXXX XXXX XXXX	12/09/2019	Meckel's Imaging	Sodium Pertechnetate (Tc-99m)	VERIFIED
RP1909031	HKL010511	XXXXX XXXX XXXXX	04/09/2019	Pheochromacytoma imaging	Iodine-131 Metaiodobenzylguanidine (I-131 MIBG) 1 mCi	VERIFIED
RP1909041	HKL011002	XXXXX XXXXX XXXX	04/09/2019	Neuroendocrine Tumor Therapy	Iodine-131 Metaiodobenzylguanide (I-131 MIBG) 200mCi	VERIFIED

Figure 3.1.1-3 Radiopharmaceuticals Appointment

**STEP 5**

Double click on the patient name and system will be displayed as per Figure 3.1.1-3

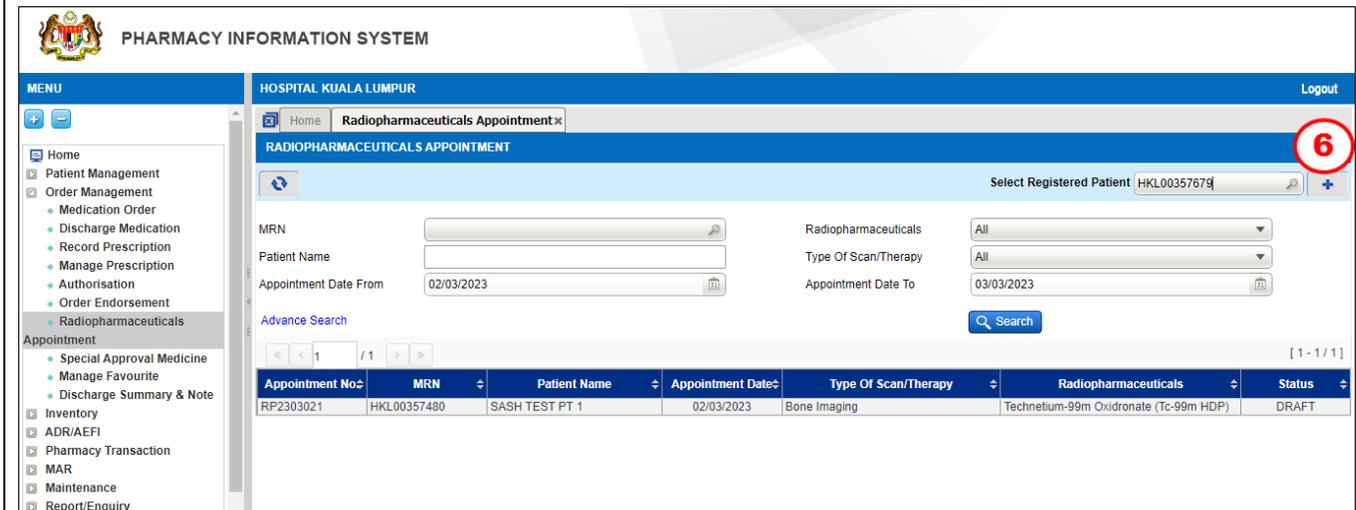

PHARMACY INFORMATION SYSTEM

HOSPITAL KUALA LUMPUR Logout

Home Radiopharmaceuticals Appointment x

RADIOPHARMACEUTICALS APPOINTMENT

Select Registered Patient

MRN

Patient Name

Appointment Date From

Radiopharmaceuticals

Type Of Scan/Therapy

Appointment Date To

Advance Search


1 / 1

Appointment No	MRN	Patient Name	Appointment Date	Type Of Scan/Therapy	Radiopharmaceuticals	Status
RP2303021	HKL00357480	SASH TEST PT 1	02/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	DRAFT

Figure 3.1.1-4 Radiopharmaceuticals Appointment

**STEP 6**

Click on  button to create new radiopharma appointment record and system will be displayed as per Figure 3.1.1-4

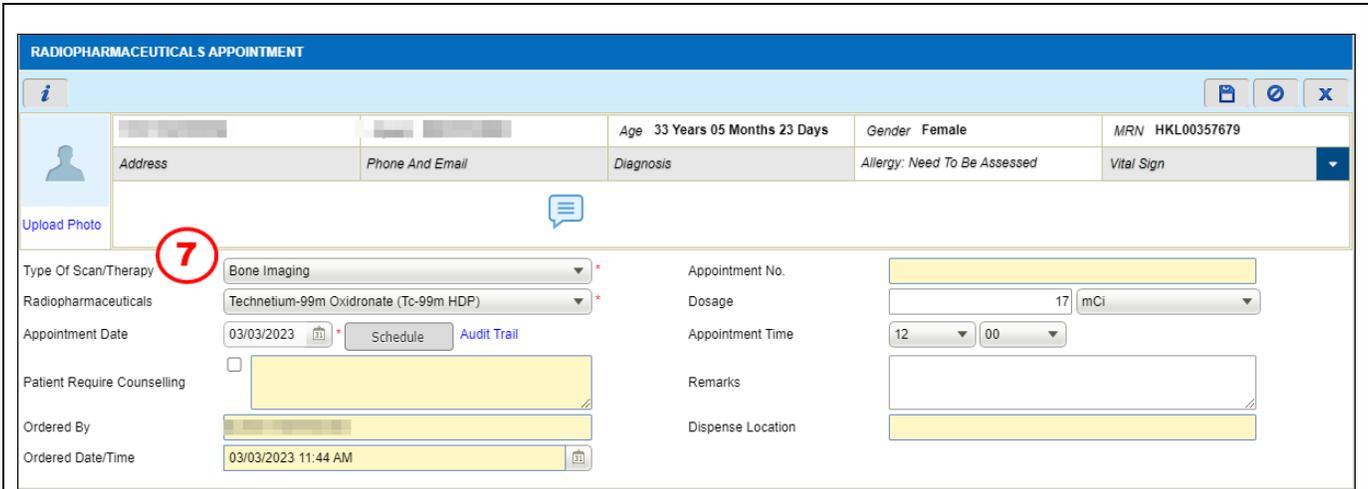

**RADIOPHARMACEUTICALS APPOINTMENT**

Age 33 Years 05 Months 23 Days | Gender Female | MRN HKL00357679

Address | Phone And Email | Diagnosis | Allergy: Need To Be Assessed | Vital Sign

Upload Photo

Type Of Scan/Therapy **7** Bone Imaging

Radiopharmaceuticals Technetium-99m Oxidronate (Tc-99m HDP)

Appointment Date 03/03/2023

Patient Require Counselling

Ordered By | Ordered Date/Time 03/03/2023 11:44 AM


Appointment No. | Dosage 17 mCi | Appointment Time 12:00

Remarks | Dispense Location

**Figure 3.1.1-5 Radiopharmaceuticals Appointment**

**STEP 7**

Select **Type Of Scan/Therapy** from the drop down box and system will be displayed as per Figure 3.1.1-5

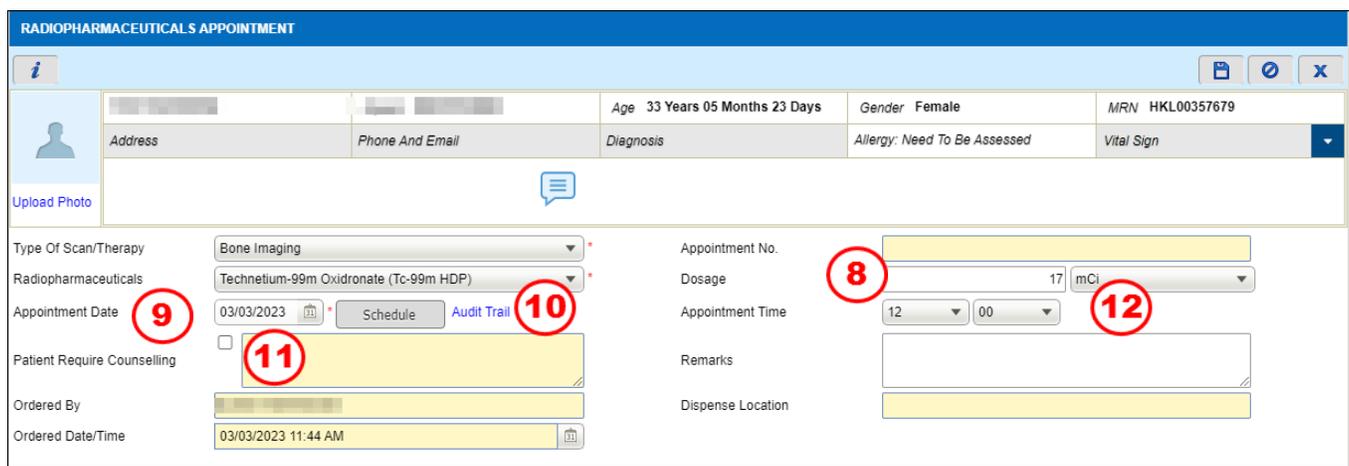

**RADIOPHARMACEUTICALS APPOINTMENT**

Age 33 Years 05 Months 23 Days | Gender Female | MRN HKL00357679

Address | Phone And Email | Diagnosis | Allergy: Need To Be Assessed | Vital Sign

Upload Photo

Type Of Scan/Therapy Bone Imaging

Radiopharmaceuticals Technetium-99m Oxidronate (Tc-99m HDP)

Appointment Date **9** 03/03/2023   **10**

Patient Require Counselling  **11**

Ordered By | Ordered Date/Time 03/03/2023 11:44 AM

Appointment No. | Dosage **8** 17 **12** mCi | Appointment Time 12:00

Remarks | Dispense Location

**Figure 3.1.1-6 Radiopharmaceuticals Appointment**

**Note**

After enter **Type of Scan/Therapy** field, **Radiopharmaceuticals** field will be automatically displayed as per Figure 3.1.1-6.

**STEP 8**

Enter **Dosage**. The UOM displayed will be based on the Radiopharmaceuticals SKU in Drug Master

**STEP 9**

Select **Appointment Date** field for appointment

**Note**

**Appointment Date** field is mandatory

**STEP 10**

Click on the  button to check for available date for appointment

**STEP 11**

Check on the  **Counseling Indicator** checkbox if applicable

### STEP 12

Select **Appointment Time** for appointment (hours and minutes) from the drop down box

**RADIOPHARMACEUTICALS APPOINTMENT**

Age 33 Years 05 Months 23 Days Gender Female MRN HKL001

Address Phone And Email Diagnosis Allergy: Need To Be Assessed Vital Sign

Type Of Scan/Therapy: Bone Imaging  
Radiopharmaceuticals: Technetium-99m Oxidronate (Tc-99m HDP)  
Appointment Date: 03/03/2023  
Appointment Time: 12:00  
Remarks:   
Dispense Location:   
Ordered By:   
Ordered Date/Time: 03/03/2023 11:44 AM

**Figure 3.1.1-7 Radiopharmaceuticals Appointment**

### STEP 13

Enter **Remarks** if applicable

#### Note

- **Ordered By** and **Ordered Date/Time** will automatically generate by system.
- **Dispense Location** will only be automatically generated in Medication Order Radiopharmaceuticals screen. Not in Radiopharmaceuticals screen.

### STEP 14

Click on the save  button to save record and alert message will be displayed as per Figure 3.1.1-7

**CONFIRMATION**

CONF0002: Are you sure you want to save?

Cancel Yes No

**Figure 3.1.1-8 Save Record Alert Message**

### STEP 15

Click on the  button to generate number for scan appointment and alert message will be displayed as per Figure 3.1.1-8

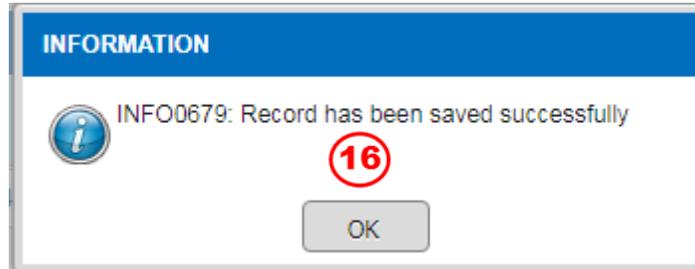
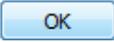


Figure 3.1.1-9 Save Record Alert Message

**STEP 16**

Click on the  button to save and update Radiopharmaceuticals Appointment order

**Note**

Record status is updated to Draft.

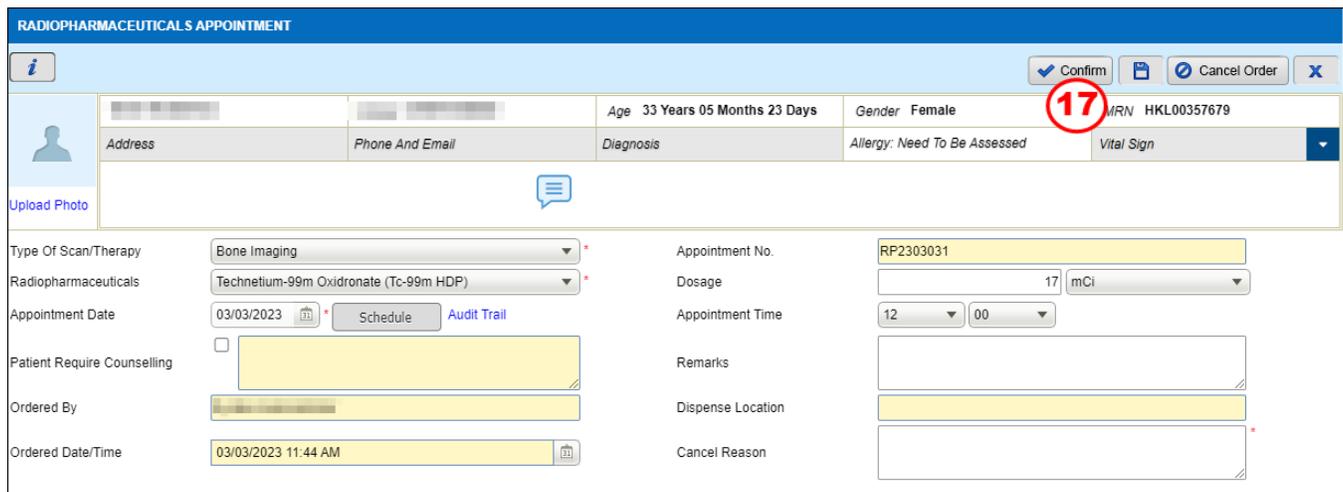


Figure 3.1.1-10 Radiopharmaceuticals Appointment

**STEP 17**

Click on the  button and system will display alert message as per Figure 3.1.1-11


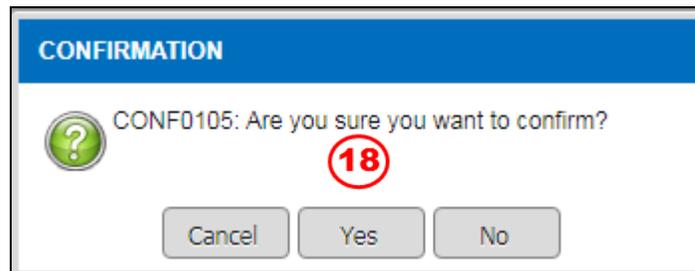


Figure 3.1.1-11 Confirm Order Alert Message

**STEP 18**

Click on the  button to confirm order

**RADIOPHARMACEUTICALS APPOINTMENT**

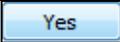

 Confirm


	<input type="text"/>	Age 33 Years 05 Months 23 Days	Gender Female	MRN HKL00357 <span style="border: 2px solid red; border-radius: 50%; padding: 2px;">19</span>
Address	Phone And Email	Diagnosis	Allergy: Need To Be Assessed	Vital Sign <span style="float: right;">▼</span>

Upload Photo 

Type Of Scan/Therapy	<input type="text" value="Bone Imaging"/>	Appointment No.	<input type="text" value="RP2303031"/>
Radiopharmaceuticals	<input type="text" value="Technetium-99m Oxidronate (Tc-99m HDP)"/>	Dosage	<input type="text" value="17"/> mCi
Appointment Date	<input type="text" value="03/03/2023"/> <input type="button" value="Schedule"/> <input type="button" value="Audit Trail"/>	Appointment Time	<input type="text" value="12"/> <input type="text" value="00"/>
Patient Require Counselling	<input type="checkbox"/>	Remarks	<input type="text"/>
Ordered By	<input type="text"/>	Dispense Location	<input type="text"/>
Ordered Date/Time	<input type="text" value="03/03/2023 11:44 AM"/>	Cancel Reason	<input type="text"/> <span style="border: 2px solid red; border-radius: 50%; padding: 2px;">20</span>

**Figure 3.1-12 Radiopharmaceuticals Appointment**

**STEP 19**

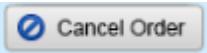Click on the  button to cancel Radiopharmaceuticals Appointment

**STEP 20**

Enter cancel reason and record status will be updated as per Figure 3.1-13



**PHARMACY INFORMATION SYSTEM**

Logout

**MENU**

- Home
- ⊖ Patient Management
- ⊖ Order Management
  - ⊖ Medication Order
  - ⊖ Discharge Medication
  - ⊖ Record Prescription
  - ⊖ Manage Prescription
  - ⊖ Authorisation
  - ⊖ Order Endorsement
  - ⊖ Radiopharmaceuticals
- ⊖ Appointment
  - ⊖ Special Approval Medicine
  - ⊖ Manage Favourite
  - ⊖ Discharge Summary & Note
- ⊖ Inventory
- ⊖ ADR/AEFI
- ⊖ Pharmacy Transaction
- ⊖ MAR
- ⊖ Maintenance
- ⊖ Report/Enquiry
- ⊖ PE Reports

HOSPITAL KUALA LUMPUR

Home Radiopharmaceuticals Appointment

**RADIOPHARMACEUTICALS APPOINTMENT**

MRN 
Radiopharmaceuticals

Patient Name 
Type Of Scan/Therapy

Appointment Date From 
Appointment Date To


Advance Search

Appointment No.	MRN	Patient Name	Appointment Date	Type Of Scan/Therapy	Radiopharmaceuticals	Status
RP2303031	HKL00357679	TEST PATIENT04	03/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	CANCELLED
RP2303021	HKL00357480	SASH TEST PT 1	02/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	DRAFT

**Figure 3.1.1-13 Radiopharmaceuticals Appointment**

### 3.2 Search Registered Patient

To view an existing patient record, please perform the steps below.

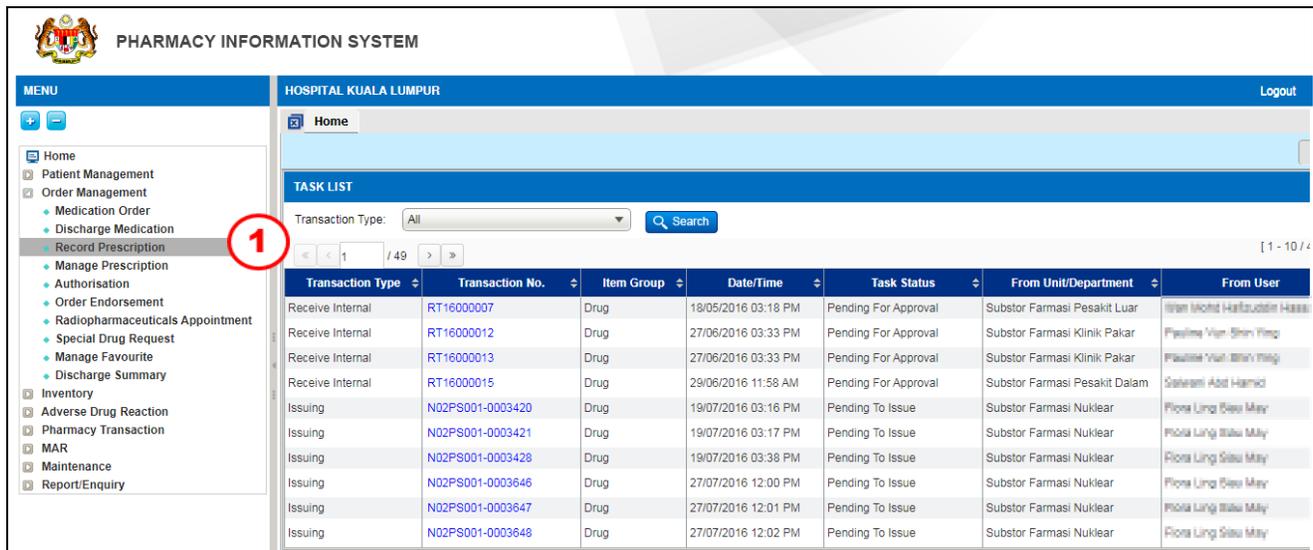

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR Logout

Home

**TASK LIST**

Transaction Type: All

[ 1 - 10 / 4 ]


Transaction Type	Transaction No.	Item Group	Date/Time	Task Status	From Unit/Department	From User
Receive Internal	RT16000007	Drug	18/05/2016 03:18 PM	Pending For Approval	Substor Farmasi Pesakit Luar	Wan Mohd Hafizudin Hassan
Receive Internal	RT16000012	Drug	27/06/2016 03:33 PM	Pending For Approval	Substor Farmasi Klinik Pakar	Pasime Wan Shu Ting
Receive Internal	RT16000013	Drug	27/06/2016 03:33 PM	Pending For Approval	Substor Farmasi Klinik Pakar	Pasime Wan Shu Ting
Receive Internal	RT16000015	Drug	29/06/2016 11:58 AM	Pending For Approval	Substor Farmasi Pesakit Dalam	Saleem Abd Hamid
Issuing	N02PS001-0003420	Drug	19/07/2016 03:16 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May
Issuing	N02PS001-0003421	Drug	19/07/2016 03:17 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May
Issuing	N02PS001-0003428	Drug	19/07/2016 03:38 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May
Issuing	N02PS001-0003646	Drug	27/07/2016 12:00 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May
Issuing	N02PS001-0003647	Drug	27/07/2016 12:01 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May
Issuing	N02PS001-0003648	Drug	27/07/2016 12:02 PM	Pending To Issue	Substor Farmasi Nuklear	Fiona Ling Siew May

Figure 3.2-1 PhIS Home page

This page is used to search for a registered patient

#### STEP 1

Click on 'Order Management' menu followed by 'Record Prescription' and screen will be displayed as per Figure 3.1-1

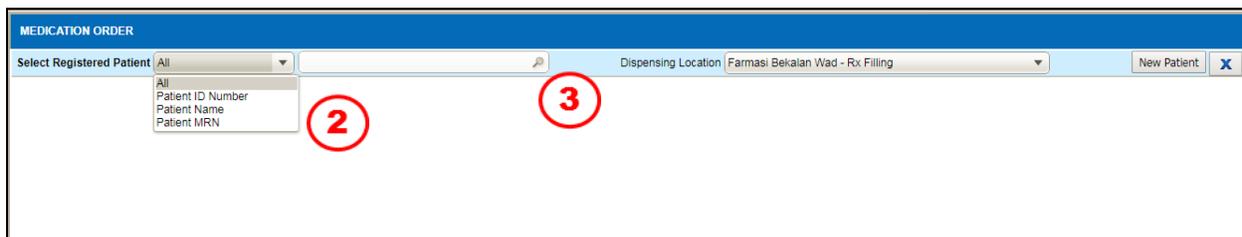

**MEDICATION ORDER**

Select Registered Patient: All  Dispensing Location: Farmasi Bekalan Wad - Rx Filling

All  
 Patient ID Number  
 Patient Name  
 Patient MRN

Figure 3.2-2 Medication Order

#### STEP 2

Select registration type from dropdown box

#### STEP 3

Click on the  button to search for a registered patient based on the type selected

#### Note

If the patient is not registered, click on the  button and the 'Patient Registration' screen will be displayed.

**MEDICATION ORDER**

Select Registered Patient: All | Dispensing Location: Farmasi Klinik Pakar | New Patient X

Patient Search [1 - 10 / 234394]

MRN	ID Number	Patient Name
HKL1796229	5546237283208	5546237283208
HKL228475	5546237283208	5546237283208
HKL2078237	5546237283208	5546237283208
HKL254362	5546237283208	5546237283208
HKL098421	5546237283208	5546237283208
HKL2078237	5546237283208	5546237283208
HKL107549	5546237283208	5546237283208
HKL154319	5546237283208	5546237283208
HKL188432	5546237283208	5546237283208
HKL213473	5546237283208	5546237283208

Figure 3.2-3 Search Patient

**STEP 4**

Double-click on the selected patient and the details will be displayed as per Figure 3.1-4

**Note**

- Click on the 'Renew Visit' hyperlink to renew a visit for the same location.
- Select the 'Visit Type' and 'Location' to create a new visit.
- Click on the Visit Management button to access to 'Visit Management' screen.

**MEDICATION ORDER**

Select Registered Patient: All | HKL00357679 | Dispensing Location: Farmasi Nuklear | Visit Management | New Patient X

Age: 33 Years 05 Months 23 Days | Gender: Female | MRN: HKL00357679

Address: | Phone And Email: | Diagnosis: | Allergy: Need To Be Assessed | Vital Sign: [v]

Height: [ ] cm | Weight: 55 kg | BMI/BSA: 0/0 m<sup>2</sup> | G6PD Deficiency: [ ] | Update (Last Updated: ) | Nationality: Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | Other Appointment | VAS History |

**VISIT/ADMISSION**

Encounter No: [ ] | Visit Type: Inpatient | Location: [ ] | Department: General Medicine | Clear | Save | Patient Label | Patient sticker

Full Paying Patient | Facility RN Number: [ ] | Bed No: [ ] |  Breast Feeding |  Pregnant | Week: [ ] | MFlex In: [ ]

[1 - 3 / 7]

Encounter No	Visit/Admission Date	Visit Type	Renew Visit	Visit Location	Department	Discharge Date
0001	26/09/2022	Outpatient	Renew Visit	Farmasi Klinik Pakar	Pharmacy	
0010	12/08/2022	Outpatient	Renew Visit	Farmasi Klinik Pakar	Pharmacy	
0001	09/08/2022	Outpatient	Renew Visit	Farmasi Klinik Pakar	Pharmacy	

**PARTIAL SUPPLY**

[1 - 1 / 1]

Order Date	Order Location	Visit Type	Department	Description No	Original Prescription	Next Collection Date	Next Collection Mode	Dispense Location
------------	----------------	------------	------------	----------------	-----------------------	----------------------	----------------------	-------------------

Figure 3.2-4 Medication Order Screen

**MEDICATION ORDER**

Select Registered Patient: All | HKL00357679 | Dispensing Location: Farmasi Nuklear | Visit Management | New Patient

---

Address | Phone And Email | Age: 33 Years 05 Months 23 Days | Gender: Female | MRN: HKL00357679  
Diagnosis | Allergy: Need To Be Assessed | Vital Sign  
 Height:  cm | Weight: 55 kg | BMI/BSA: 0/0 m<sup>2</sup> | G6PD Deficiency:  | Update (Last Updated: ) | Nationality: Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | Other Appointment | VAS History

---

**VISIT/ADMISSION**

Encounter No:  | Visit Type: Inpatient | Location:  | Department: General Medicine | Clear | Save | Patient Label | Patient sticker

Full Paying Patient | Facility RN Number:  | Bed No:  |  Breast Feeding |  Pregnant | Week:  |  MFlex In

Encounter No	Visit/Admission Date	Visit Type	Renew Visit	Visit Location	Department	Discharge Date
0001	26/09/2022	Outpatient	Renew Visit <span style="border: 2px solid red; border-radius: 50%; padding: 2px;">5</span>	Farmasi Klinik Pakar	Pharmacy	
0010	12/08/2022	Outpatient	Renew Visit	Farmasi Klinik Pakar	Pharmacy	
0001	09/08/2022	Outpatient	Renew Visit	Farmasi Klinik Pakar	Pharmacy	

---

**PARTIAL SUPPLY**

Order Date | Order Location | Visit Type | Department | Prescription No | Original Prescription | Next Collection Date | Next Collection Mode | Dispense Location

**Figure 3.2-5 Visit/Admission Detail**

### STEP 5

Double-click on the selected visit and the details will be displayed as per Figure 3.1-6

**MEDICATION ORDER**

Select Registered Patient: All | HKL00357679 | Dispensing Location: Farmasi Nuklear | Visit Management | New Patient

---

Address | Phone And Email | Age: 33 Years 05 Months 23 Days | Gender: Female | MRN: HKL00357679  
Diagnosis | Allergy: Need To Be Assessed | Vital Sign  
 Height:  cm | Weight: 55 kg | BMI/BSA: 0/0 m<sup>2</sup> | G6PD Deficiency:  | Update (Last Updated: ) | Nationality: Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | Other Appointment | VAS History

---

**VISIT/ADMISSION**

Encounter No: 0001 | Visit Type: Outpatient | Location: Farmasi Nuklear | Department: Nuclear Medicine | Clear | Patient Label | Patient sticker

Breast Feeding |  Pregnant | Week:  |  MFlex In

---

**VISIT RECORD**

Original Prescription Serial Number:

Patient Diagnosis

Allergy

Adverse Drug Reaction

Nutritional Assessment

Normal Order

**NORMAL ORDER - DRUG DETAILS**

Total Price: RM 0.00

Original Prescriber Name:  | Advanced Search


Drug:

Alert	Drug Name	Order Details	Order Qty	Drug Info

**Figure 3.2-6 Medication Order Detail**

### 3.3 Radiopharmaceuticals

Order Medication order for Radiopharmaceuticals screen will allow user to Confirm radiopharmaceuticals order

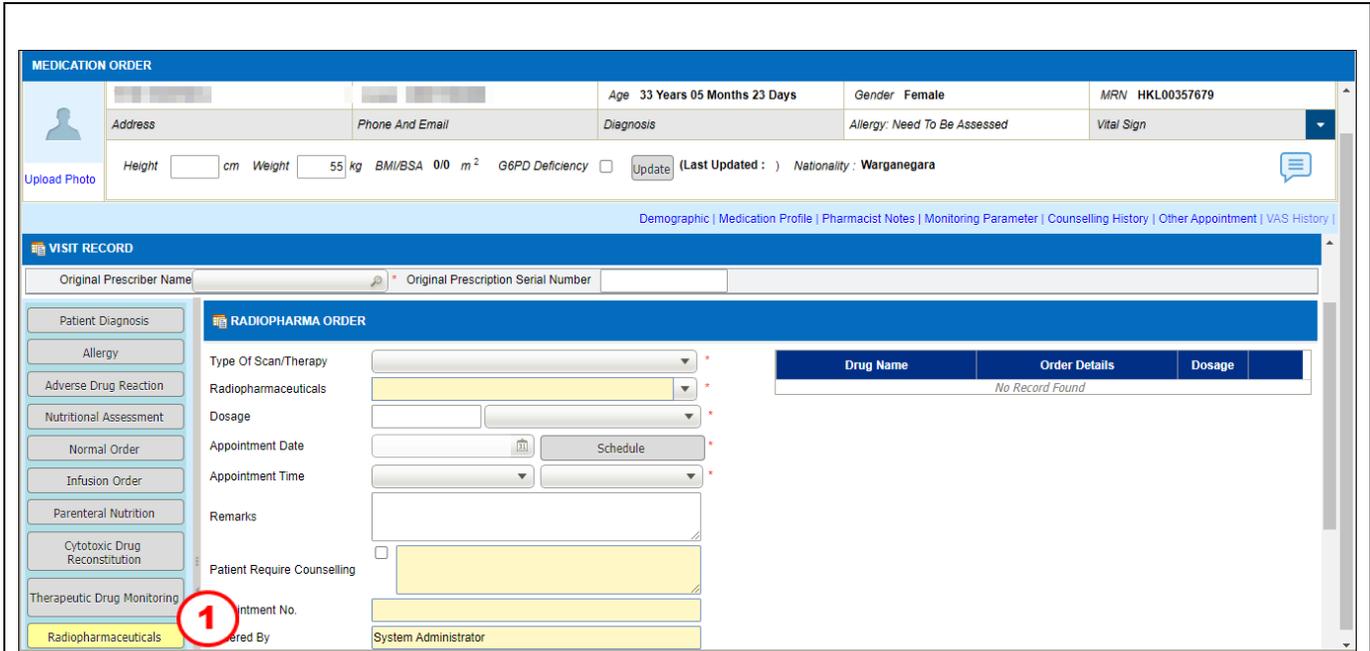

**MEDICATION ORDER**

Age 33 Years 05 Months 23 Days Gender Female MRN HKL00357679

Address Phone And Email Diagnosis Allergy: Need To Be Assessed Vital Sign

Height cm Weight 55 kg BMI/BSA 0/0 m<sup>2</sup> G6PD Deficiency  Update (Last Updated : ) Nationality : Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | Other Appointment | VAS History |

**VISIT RECORD**

Original Prescriber Name Original Prescription Serial Number

**RADIOPHARMA ORDER**

Type Of Scan/Therapy Radiopharmaceuticals Dosage Appointment Date Appointment Time Remarks

Patient Require Counselling

Appointment No. Entered By System Administrator

Drug Name	Order Details	Dosage
No Record Found		

Figure 3.3-1 Medication Order –Radiopharmaceuticals


**STEP 1**

Click on the **Radiopharmaceuticals** button to display Radiopharmaceuticals medication order screen as per Figure 3.3-1

**Note**

**Current Appointment section :**

This section will automatically display appointment records created in Appointment screen (if any). Only appointment records with 'New' status is displayed in this section.

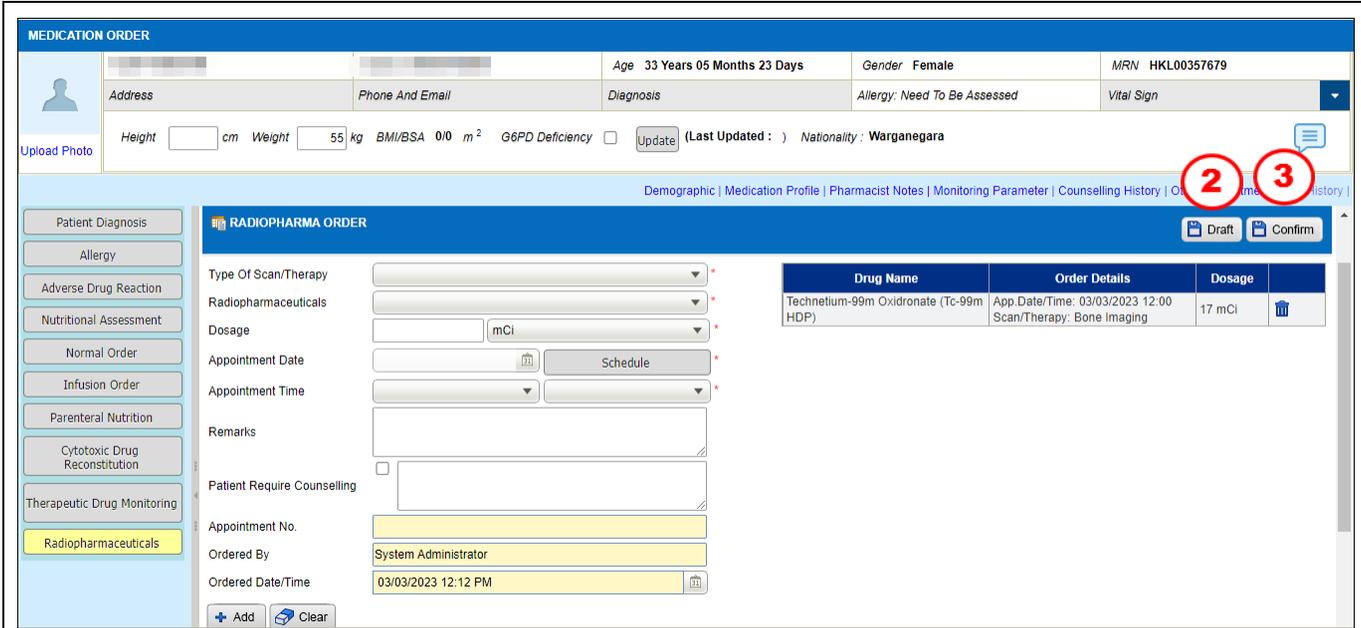

**MEDICATION ORDER**

Age 33 Years 05 Months 23 Days Gender Female MRN HKL00357679

Address Phone And Email Diagnosis Allergy: Need To Be Assessed Vital Sign

Height  cm Weight 55 kg BMI/BSA 0/0 m<sup>2</sup> G6PD Deficiency  Update (Last Updated : ) Nationality : Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | O...me...history |

**RADIOPHARMA ORDER** Draft Confirm

Drug Name	Order Details	Dosage
Technetium-99m Oxidronate (Tc-99m HDP)	App.Date/Time: 03/03/2023 12:00 Scan/Therapy: Bone Imaging	17 mCi

Type Of Scan/Therapy  Radiopharmaceuticals  Dosage  mCi

Appointment Date  Schedule  Appointment Time

Remarks

Patient Require Counselling

Appointment No.

Ordered By System Administrator

Ordered Date/Time 03/03/2023 12:12 PM


+ Add Clear

Figure 3.3-2 Medication Order -Radiopharmaceuticals

**Note**

**Consent Form:**

- Right click in the  button and select Consent Form button.
- Click on the Consent Form link and screen will display as per Figure 3.3-3

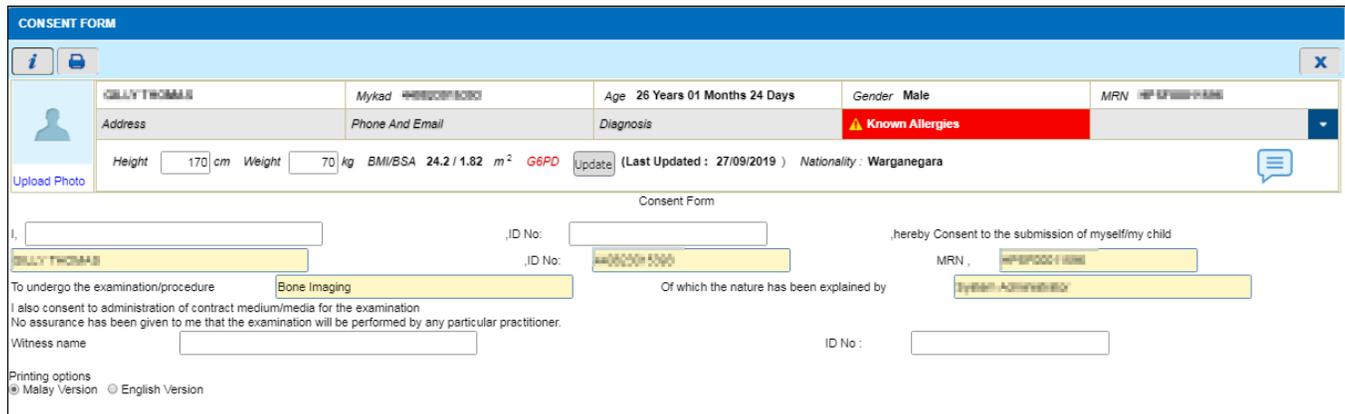

**CONSENT FORM**

Age 26 Years 01 Months 24 Days Gender Male MRN

Address Phone And Email Diagnosis **Known Allergies**

Height 170 cm Weight 70 kg BMI/BSA 24.2 / 1.82 m<sup>2</sup> G6PD Update (Last Updated : 27/09/2019 ) Nationality : Warganegara

Consent Form

I,  ID No:  hereby Consent to the submission of myself/my child

ID No:  MRN,

To undergo the examination/procedure  Of which the nature has been explained by

I also consent to administration of contract medium/media for the examination  
No assurance has been given to me that the examination will be performed by any particular practitioner.

Witness name  ID No:


Printing options  
 Malay Version  English Version

Figure 3.3-3 Consent Form

**Note**

**Consent Form:**

- Consent Form is an optional function.
- Enter Patient Name, ID. No., Witness Name and ID. No. of witness in the consent form.
- Consent form has two printing options:
  - Malay Version
  - English Version

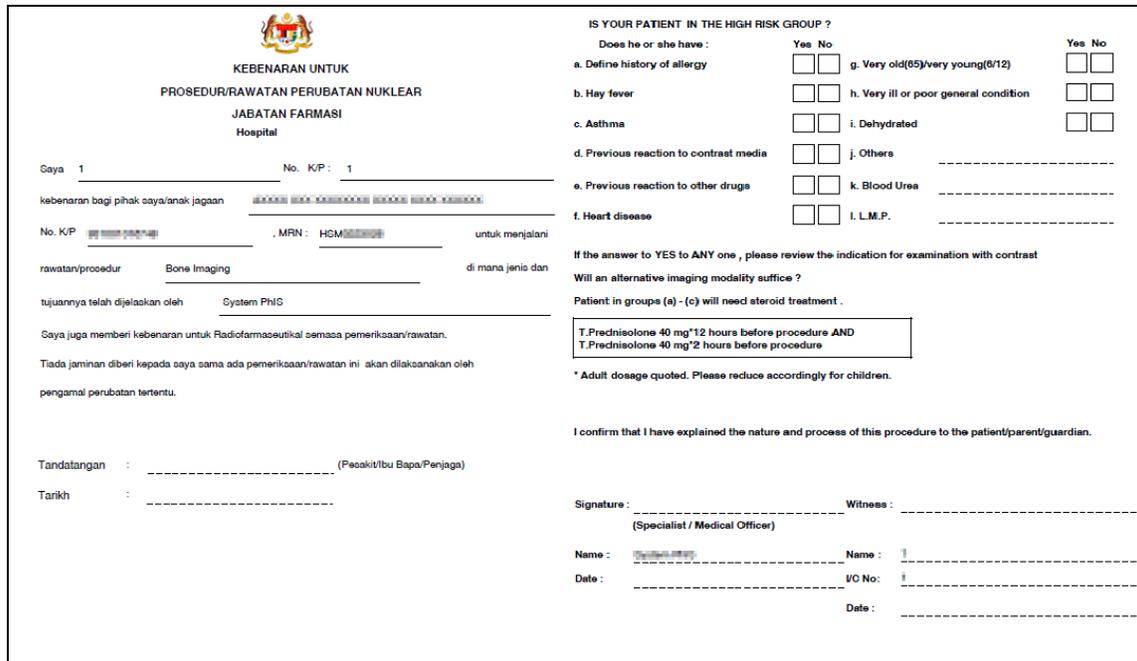

The form is titled 'KEBENARAN UNTUK PROSEDUR/RAWATAN PERUBATAN NUKLEAR JABATAN FARMASI Hospital'. It contains fields for patient identification (Name, No. K.P., MRN, HSM) and procedure details (Bone Imaging, System PhIS). A section titled 'IS YOUR PATIENT IN THE HIGH RISK GROUP?' lists various conditions (a-j) with 'Yes/No' checkboxes. Below this, there are instructions for steroid treatment (Prednisolone) and a confirmation statement. Signature lines for the patient/parent and the medical officer are provided at the bottom.

Figure 3.3-4 Sample Consent Form

**Note**

Sample of Consent Form as per Figure 3.3-4.

**STEP 2**

Click on the  button to save record and alert message will be displayed as per Figure 3.3-5

**STEP 3**

Click on the  button to confirm the record and alert message will be displayed as per Figure 3.3-6

**Note**

- Alert message will be displayed *be displayed after save the record.*

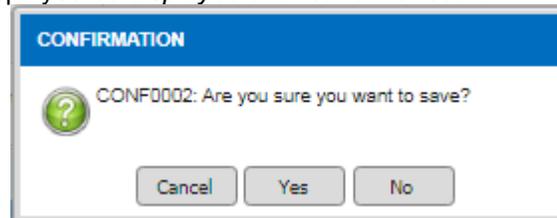


Figure 3.3-5 Alert Message

- Alert message will be displayed *be displayed after confirm the record.*


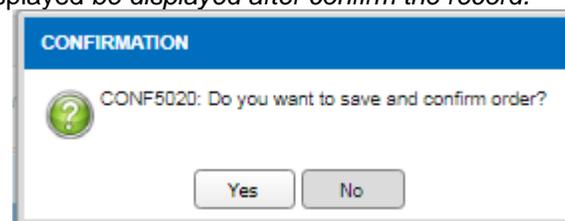


Figure 3.3-6 Alert Message

**Note**

**Cancel function:**

- *Cancel is an optional function.*
- *Click on the  button to delete radiopharmaceuticals order*
- *Need to enter cancel reason once user click cancel function:*

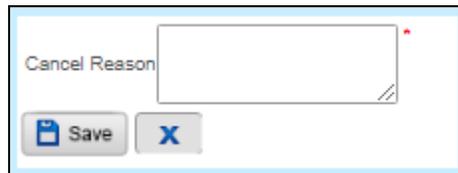


The screenshot shows a dialog box titled "Cancel Reason". It contains a text input field for entering the reason. Below the input field are two buttons: "Save" (with a floppy disk icon) and "X" (with a blue 'X' icon).

**Figure 3.3-7 Cancel Order**

- *Once confirm the order, drug details will auto appear in screening stage*

### 3.4 Radiopharmaceuticals Verification

Radiopharmaceuticals Verify screen allow user to screen prescription details and verify order.

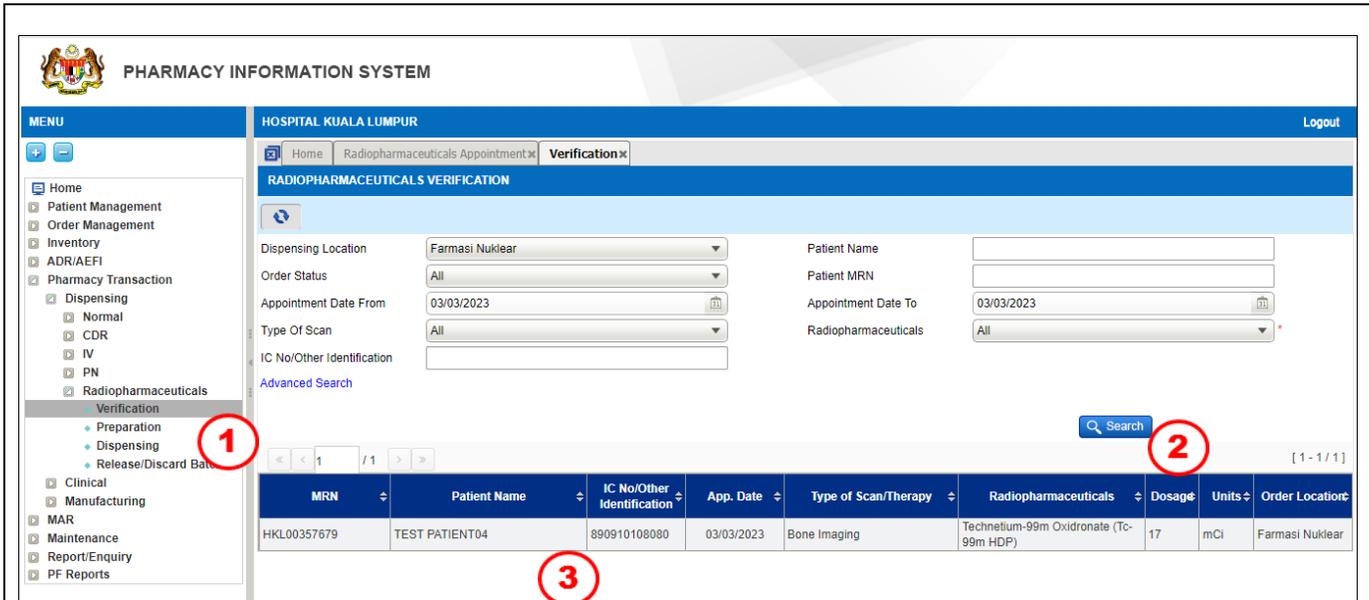

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Home | Radiopharmaceuticals Appointment | **Verification** | Logout

**RADIOPHARMACEUTICALS VERIFICATION**

Dispensing Location: Farmasi Nuklear  
 Order Status: All  
 Appointment Date From: 03/03/2023  
 Type Of Scan: All  
 IC No/Other Identification:

Patient Name:   
 Patient MRN:   
 Appointment Date To: 03/03/2023  
 Radiopharmaceuticals: All

Advanced Search

Search [ 1 - 1 / 1 ]

MRN	Patient Name	IC No/Other Identification	App. Date	Type of Scan/Therapy	Radiopharmaceuticals	Dosage	Units	Order Location
HKL00357679	TEST PATIENT04	890910108080	03/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	17	mCi	Farmasi Nuklear

Figure 3.4-1 Radiopharmaceuticals Verification Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu followed by 'Dispensing' and click on 'Radiopharmaceuticals Verification' sub-menu

#### Note

- Various search criteria is provided as below:
  - Dispense Location
  - Patient Name
  - Order Status
  - Patient MRN
  - Appointment Date From
  - Appointment Date To
  - Type of Scan
  - Radiopharmaceuticals
  - IC No/Other Identification
- Click on the [Advanced Search](#) hyperlink for advance search.

#### STEP 2

Click on the [Search](#) button and system will display related order

#### STEP 3

Double click on the selected patient record and system will display Radiopharmaceuticals Verification screen as per Figure 3.4-2

**VERIFICATION**

TEST PATIENT04
Mykad 890910108080
Age 33 Years 05 Months 23 Days
Gender Female
MRN HKL00357679

Address
Phone And Email
Diagnosis
Allergy: Need To Be Assessed
Vital Sign

Height  cm
Weight  kg
BMI/BSA  m<sup>2</sup>
 (Last Updated : )
Nationality: Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counseling History | Other Appointment

Encounter No: 0001
Actual Visit Date: 03/03/2023
Attending Practitioner:
Visit Location: Farmasi Nuklear
Rx No: FNUK0001371631
Order By:
Prescriber Dept.:

**PRESCRIPTION**

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Technetium-99m Oxidronate (Tc-99m HDP)	App Date/Time: 03/03/2023 12:00 Scan/Therapy: Bone Imaging	<input checked="" type="checkbox"/> Details	17 mCi	0 mCi	Remarks	Ordered


**Figure 3.4-2 Radiopharmaceuticals Verification**

**STEP 4**

Click on  button to verify order

**Note**


- **Drug Counselling** column is an optional function.
- Choose **Yes** if the patient need counselling or **No** if the patient do not need counselling.

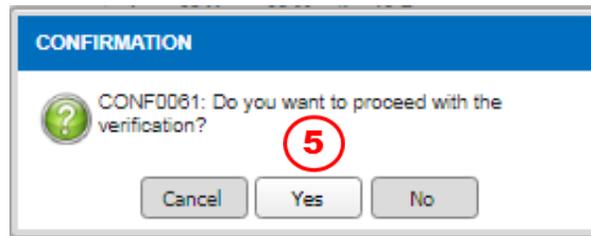

**Figure 3.4-3 Information Alert Message**

**STEP 5**

Click on the  button to proceed with the verification

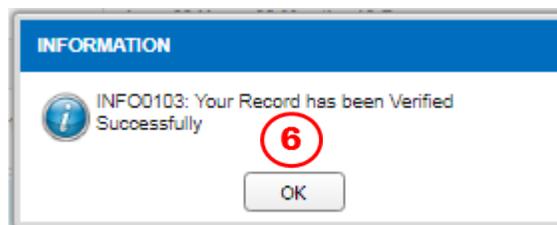


**Figure 3.4-4 Verify Order Alert Message**

**STEP 6**

Click on the  button to verify successfully

### 3.4.1 Intervention

Radiopharmaceuticals verification screen will allow user to verify prescription, before preparation and dispensing stage

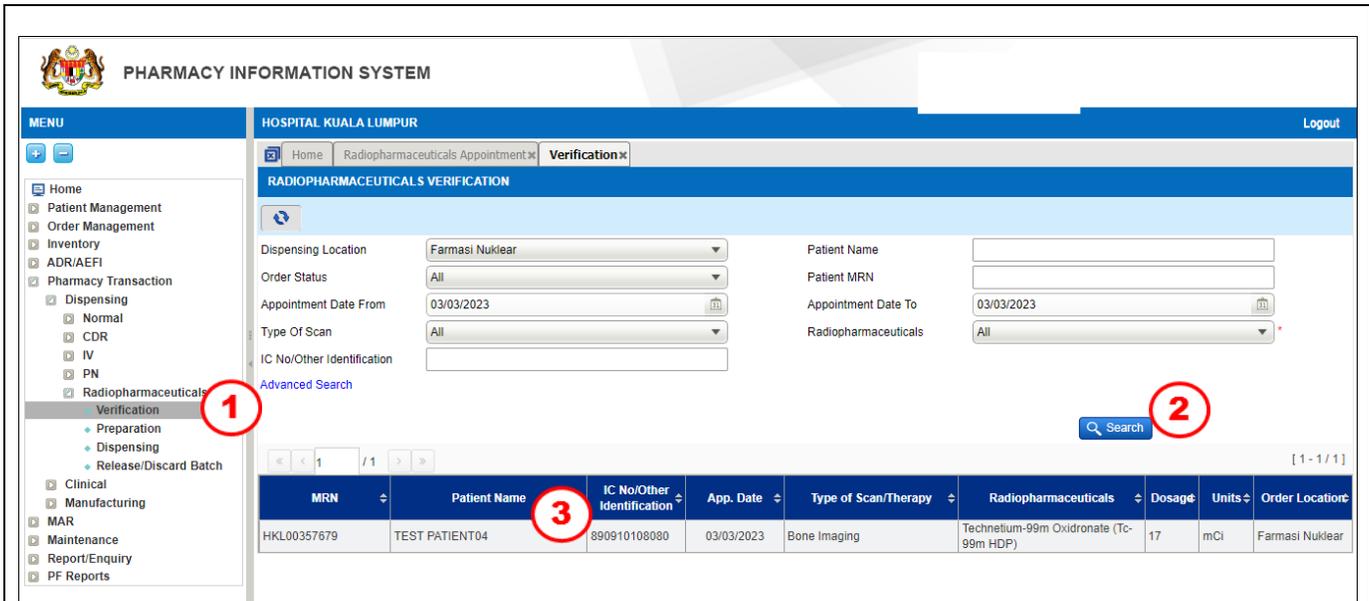

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Home Radiopharmaceuticals Appointment Verification

**RADIOPHARMACEUTICALS VERIFICATION**

Dispensing Location: Farmasi Nuklear  
 Order Status: All  
 Appointment Date From: 03/03/2023  
 Type Of Scan: All  
 IC No/Other Identification:

Patient Name:   
 Patient MRN:   
 Appointment Date To: 03/03/2023  
 Radiopharmaceuticals: All

Advanced Search

Search

MRN	Patient Name	IC No/Other Identification	App. Date	Type of Scan/Therapy	Radiopharmaceuticals	Dosage	Units	Order Location
HKL00357679	TEST PATIENT04	890910108080	03/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	17	mCi	Farmasi Nuklear

Figure 3.4.1-1 Radiopharmaceuticals Verification Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu followed by 'Dispensing' and click on 'Radiopharmaceuticals Verification' sub-menu

#### Note

- Various search criteria is provided as below:
  - Dispense Location
  - Patient Name
  - Order Status
  - Patient MRN
  - Appointment Date From
  - Appointment Date To
  - Type of Scan
  - Radiopharmaceuticals
  - IC No/Other Identification
- Click on the [Advanced Search](#) hyperlink for advance search.

#### STEP 2

Click on the  button and system will display related order

#### STEP 3

Double click on the selected patient record and system will be displayed Radiopharmaceuticals Verification screen as per Figure 3.3.1-2

**VERIFICATION**

TEST PATIENT04    Mykad 890910108080    Age 33 Years 05 Months 23 Days    Gender Female    MRN HKL00357679

Address    Phone And Email    Diagnosis    Allergy: Need To Be Assessed    Vital Sign

Height  cm    Weight  55 kg    BMI/BSA 0/0 m<sup>2</sup>    Update (Last Updated : )    Nationality : Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counselling History | Other Appointment |

Encounter No: 0001    Actual Visit Date: 03/03/2023    Attending Practitioner:    Visit Location: Farmasi Nuklear    Rx No: FNUK0001371631    Order By:    Prescriber Dept.:

**PRESCRIPTION**    Verify

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Technetium-99m HDP (Tc-99m HDP)	App Date/Time: 03/03/2023 12:00 Scan/Therapy: Bone Imaging	<input checked="" type="checkbox"/> Details	17 mCi	0 mCi	Remarks	Ordered

Figure 3.4.1-2 Radiopharmaceuticals Verification Listing Page

**STEP 4**

Right click in the  button to select either **Intervention** or **Hold Order** as per Figure 343.1-2

**STEP 5**

Click on the  button to modify drug details

**INTERVENTION**    10

PHARMACY REMARKS    [ 1 - 10 / 12 ]

Intervention Reason	Remarks	Intervention Accepted	Intervention Not Accepted
<input type="checkbox"/> Authenticity		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Contraindication		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Drug Interaction		<input type="radio"/>	<input type="radio"/>
<input checked="" type="checkbox"/> Inappropriate Dose    6	7	<input checked="" type="radio"/> 8	<input type="radio"/>
<input type="checkbox"/> Inappropriate Drug		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Inappropriate Duration		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Inappropriate Frequency		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Incompatibility		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Others		<input type="radio"/>	<input type="radio"/>
<input type="checkbox"/> Polypharmacy		<input type="radio"/>	<input type="radio"/>

**MODIFY ORDER**    9

Appointment No. RP160524634    Type Of Scan/Therapy Lung V/Q Perfusion

Radiopharmaceuticals Technetium-99m Macroaggregated Albumin (Tc-99m MAA)    Dosage 115 MBq

Date 24/05/2016    View Schedule    Confirmed By MAIKAL FARDHA

Time 14 : 15    Confirmed Date 24/05/2016

Figure 3.4.1-3 Radiopharmaceuticals Intervention

**STEP 6**

Select **Intervention Reason** from list.


**STEP 7**

Enter **Remarks** if applicable.

**STEP 8**


Select **Intervention Status** radio button either **Accepted** or **Not Accepted** and edit **Amount** field

**STEP 9**

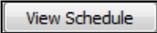Click on the  button check appointment if fully booked for a particular date.


**STEP 10**


Click on the  button to save intervention record and alert message will display as per Figure 3.4.1-4.

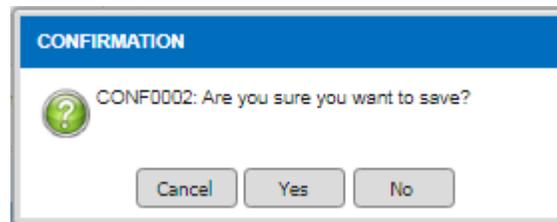

**Figure 3.4.1-4 Save record Alert Message**

**Note**

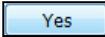- Click on the  button to save intervention record.

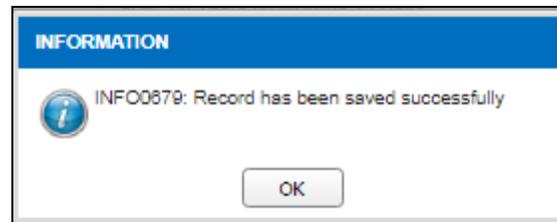


**Figure 3.4.1-5 Modify Order Alert Message**

**Note**

- Click on the  button to save modify order record.

### 3.4.2 Hold/Resume Order

Radiopharmaceuticals verification screen allow user to hold order and add reason for hold. User able to resume the order.

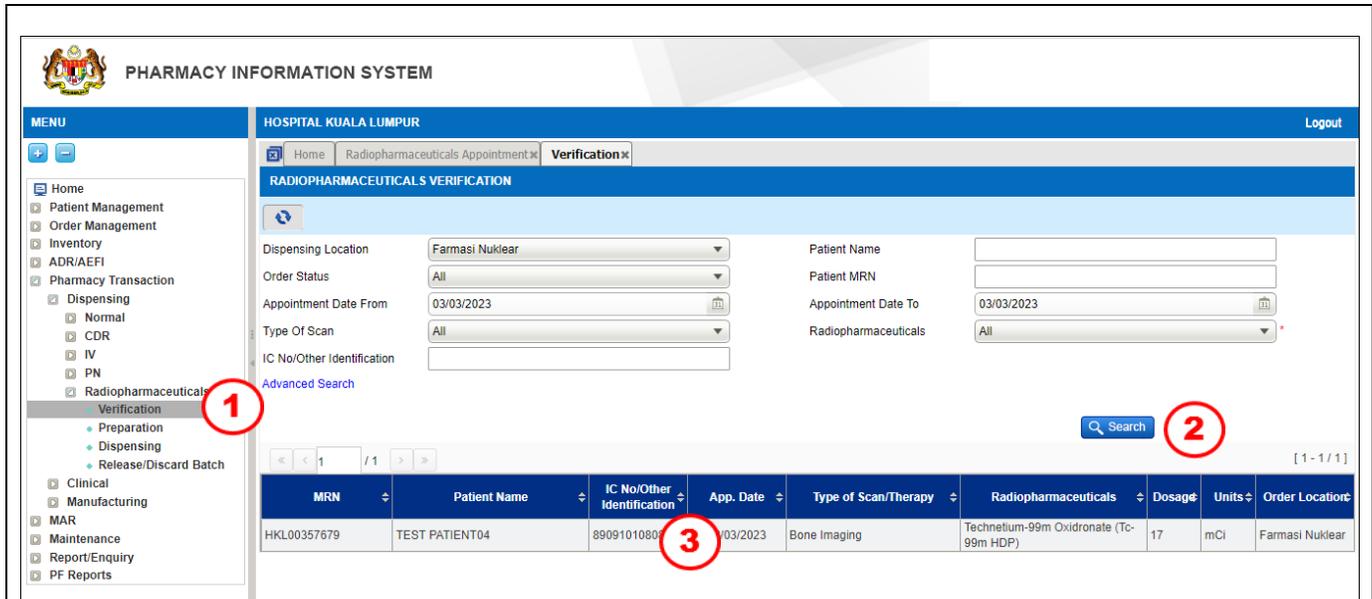

The screenshot displays the 'RADIOPHARMACEUTICALS VERIFICATION' page. The left sidebar contains a 'MENU' with 'Verification' highlighted (1). The main area has search filters for Dispensing Location, Order Status, Appointment Date, Type Of Scan, and IC No/Other Identification. A 'Search' button (2) is located below the filters. Below the filters is a table with the following data:

MRN	Patient Name	IC No/Other Identification	App. Date	Type of Scan/Therapy	Radiopharmaceuticals	Dosage	Units	Order Location
HKL00357679	TEST PATIENT04	8909101080	03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	17	mCi	Farmasi Nuklear

The patient record in the table is highlighted with a red circle (3).

Figure 3.4.2-1 Radiopharmaceuticals Verification Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu followed by 'Dispensing' and click on 'Radiopharmaceuticals Verification' sub-menu

#### Note

Various search criteria is provided as below:

- Dispense Location
- Patient Name
- Order Status
- Patient MRN
- Appointment Date From
- Appointment Date To
- Type of Scan
- Radiopharmaceuticals
- IC No/Other Identification

#### STEP 2

Click on the  button and system will display related order

#### STEP 3

Double click on the selected patient record and system will be displayed Radiopharmaceuticals Verification screen as per Figure 3.4.2-2

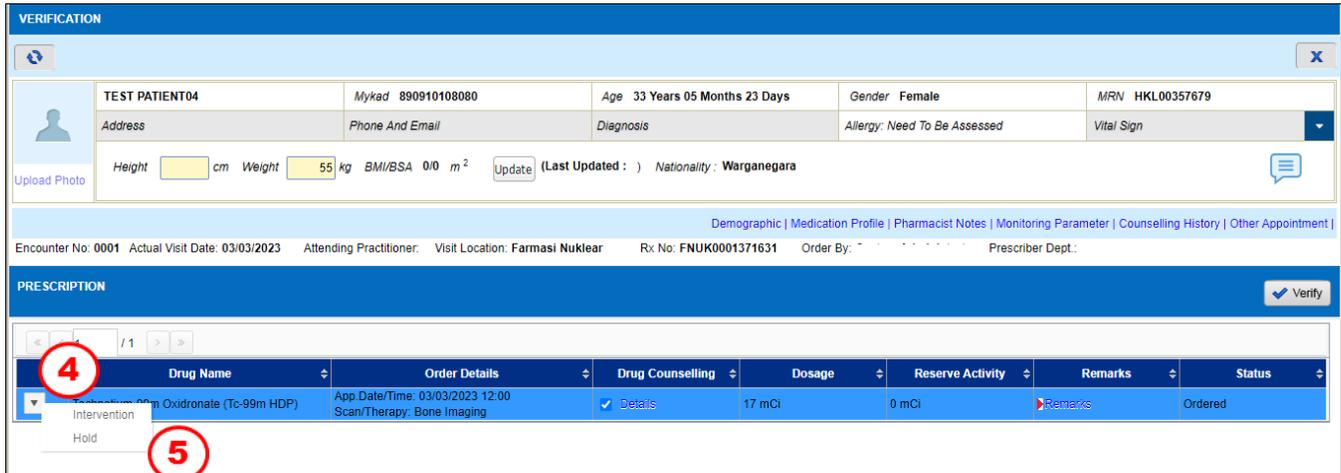


Figure 3.4.2-2 Radiopharmaceuticals Verification Listing Page

**STEP 4**

Right click in the  button to select either **Intervention** or **Hold Order**

**STEP 5**

Click on the  button to hold order

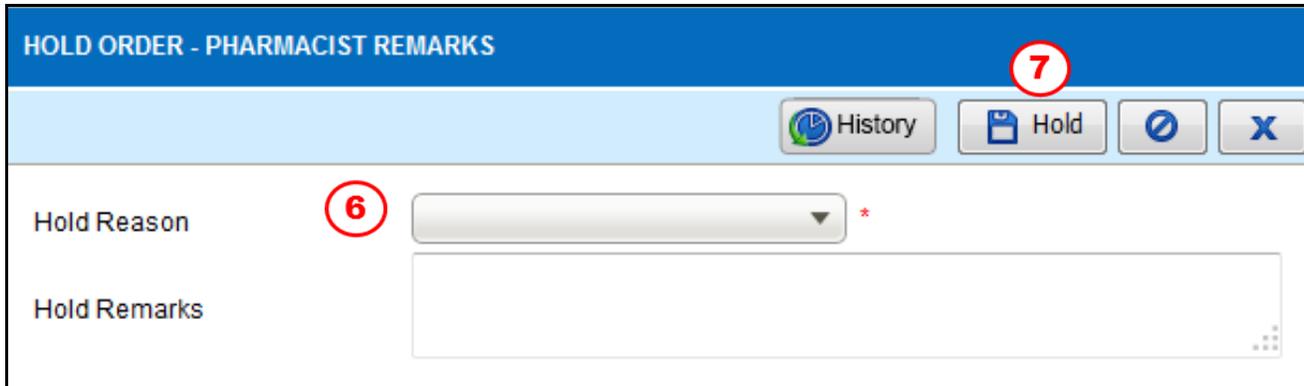


Figure 3.4.2-3 Hold Order – Pharmacist Remarks

**STEP 6**


Select the **Hold Reason** from dropdown box:

- Authenticity
- Contraindication
- Drug Interaction
- Inappropriate Dose
- Inappropriate Drug
- Inappropriate Duration
- Inappropriate Frequency
- Incompatibility
- Intervention
- Others
- Polypharmacy

- Suggest to Monitor Vital Signs/Lab Result
- Wrong Patient

**STEP 7**

Enter **Hold Remarks** field. Click on the  button to hold order and alert message will display as per Figure 3.3.2-3

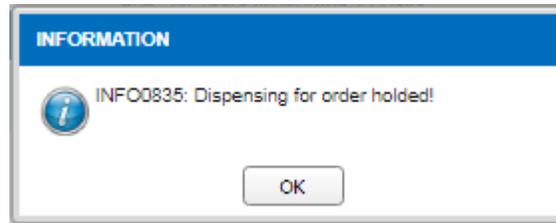


**Figure 3.4.2-4 Information Alert Message**



**Figure 3.4.2-5 Radiopharmaceuticals Verification**

**Note**

- Click on the  button to hold order.
- The  button is to cancel the entered data.
- Click on the  button to close Hold Order – Pharmacist Remarks screen.
- Status will change to Ordered Hold.
- Click on the  link to resume order as per Figure 3.3.2-6.

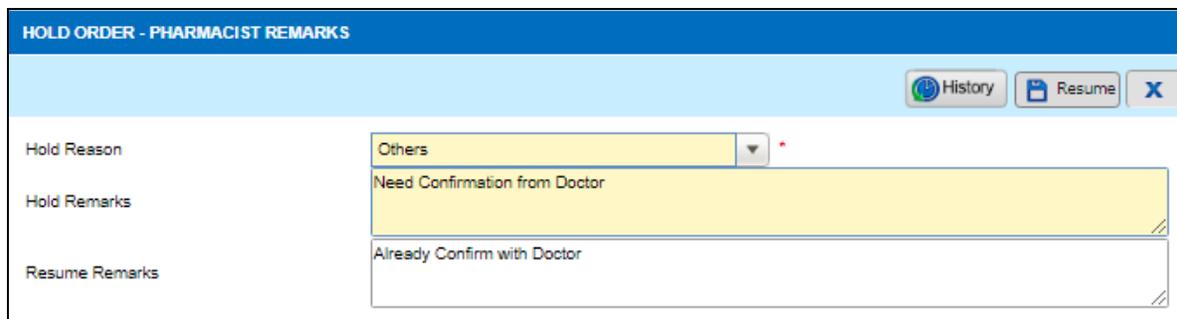


**Figure 3.4.2-6 Hold Order – Pharmacist Remarks**

**Note**


- Enter reason at **Resume Remarks** field


- Click on the  button to resume order.

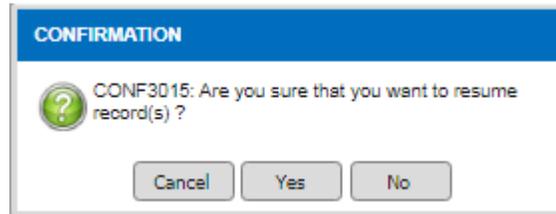

**Figure 3.4.2-7 Hold/Resume Alert Message**

**Note**

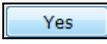
- Click on the  button to resume order.

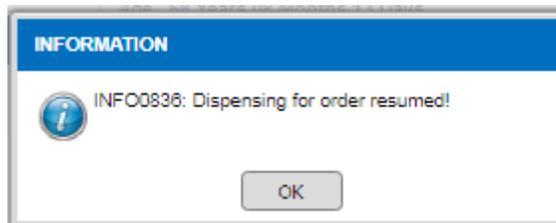

**Figure 3.4.2-8 Information Alert Message**

### 3.5 Radiopharmaceutical Preparation

Radiopharmaceuticals preparation screen allow user to allocate batch quantity before drug dispensed / administered to patient.

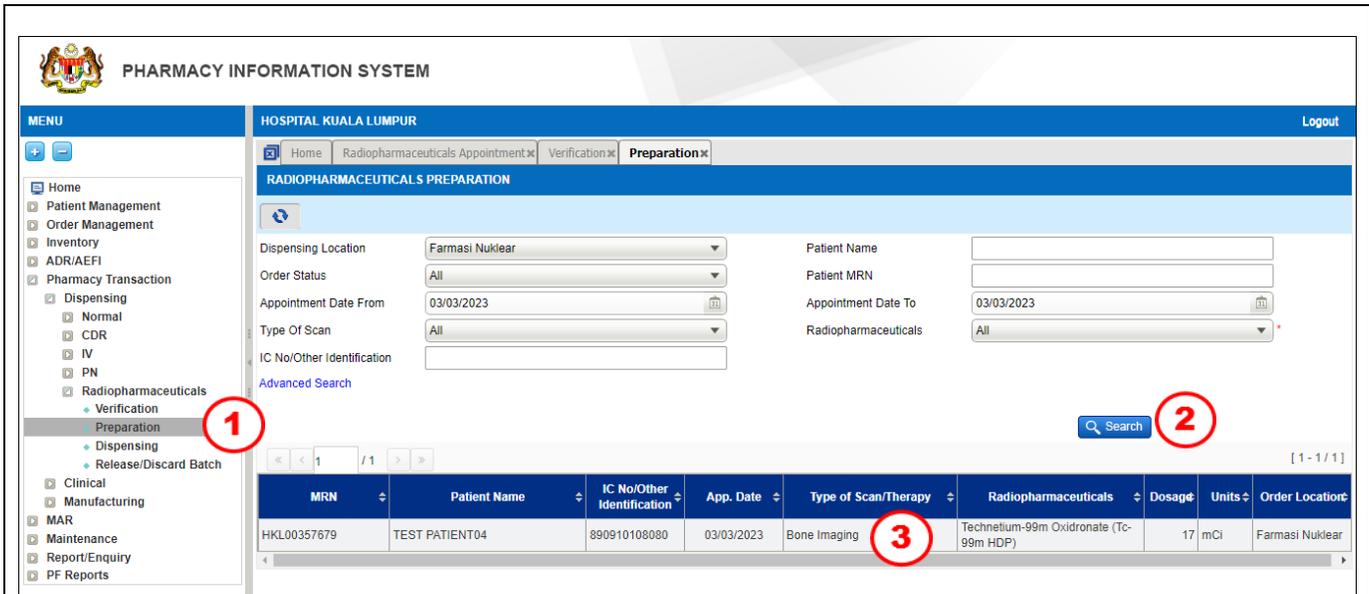

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Logout

Home Radiopharmaceuticals Appointment\* Verification\* Preparation\*

**RADIOPHARMACEUTICALS PREPARATION**

Dispensing Location: Farmasi Nuklear  
 Order Status: All  
 Appointment Date From: 03/03/2023  
 Type Of Scan: All  
 IC No/Other Identification:

Patient Name:   
 Patient MRN:   
 Appointment Date To: 03/03/2023  
 Radiopharmaceuticals: All

Advanced Search

Search

[ 1 - 1 / 1 ]

MRN	Patient Name	IC No/Other Identification	App. Date	Type of Scan/Therapy	Radiopharmaceuticals	Dosage	Units	Order Location
HKL00357679	TEST PATIENT04	890910108080	03/03/2023	Bone Imaging	Technetium-99m Oxidronate (Tc-99m HDP)	17	mCi	Farmasi Nuklear

Figure 3.5-1 Radiopharmaceuticals Preparation Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu follow by 'Dispensing' and click on 'Radiopharmaceuticals Preparation' sub-menu

#### Note

- Various search criteria is provided as below:
  - Dispense Location
  - Patient Name
  - Order Status
  - Patient MRN
  - Appointment Date From
  - Appointment Date To
  - Type of Scan
  - Radiopharmaceuticals
  - IC No/Other Identification
- Click on the [Advanced Search](#) hyperlink for advance search.

#### STEP 2

Click on the  button and system will be displayed related order

#### STEP 3

Double click on the selected patient record and system will be displayed Radiopharmaceuticals Preparation screen as per Figure 3.5-2

**PREPARATION**

Calculate Activity X

TEST PATIENT04 Mykad 890910108080 Age 33 Years 05 Months 23 Days Gender Female MRN HKL00357679

Address Phone And Email Diagnosis Allergy: Need To Be Assessed Vital Sign

Height cm Weight 55 kg BMI/BSA 0/0 m<sup>2</sup> Update (Last Updated : ) Nationality : Warganegara

Demographic | Medication Profile | Pharmacist Notes | Monitoring Parameter | Counseling History | Other Appointment

Encounter No: 0001 Actual Visit Date: 03/03/2023 Attending Practitioner: Visit Location: Farmasi Nuklear Rx No: FNUK0001371631 Order By: Prescriber Dept.:

**PRESCRIPTION** Prepare Label Online Outsource

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Techneium-99m Oxidronate (Tc-99m HDP)	App Date/Time: 03/03/2023 12:00 Scan/Therapy: Bone Imaging	<input checked="" type="checkbox"/> Details	17 mCi	0 mCi	Remarks	Verified

Figure 3.5-2 Radiopharmaceuticals Preparation

**STEP 4**

Click on the hyperlink at **Reserve Activity** to allocate quantity and Item Batch Details screen will be displayed as per Figure 3.5-3

**ITEM BATCH DETAILS**

Drug Code RPH000003 Order Detail 15.0 mCi App.Date/Time:05/05/2016 13 :0 Scan/Therapy:Diacyoscintigraphy Default SKU vial

Drug Name Sodium Perchnetate (Tc-99m) Dispense UOM mCi

Element Techneium-99m Half Life 6.01 H Allocation Date/Time 05/05/2016 11:29:23 AM

**BATCH DETAIL**

Batch No	Available Qty	Allocated Qty	Expiry Date/Time	Available Activity	Available Volume	Batch Details	Allocated Activity	Allocated Volume
45RP99	194 vial	0 vial	20/12/2016 06:03 PM	0 mCi	555 ml	Batch Details	0 mCi	ml
Total	194 vial	0 vial					0 mCi	

Figure 3.5-3 Item Batch Details

**Note**

- Item Batch Details screen will be displayed as per Figure 3.5-4.

**BATCH DETAILS**

1 / 1

1 - 1 / 1

Calibration Activity	Calibration Date/Time	Dispensed Activity	Measured Activity	Measurement Date Time	Dispensed Volume	Volume	Allocated Volume
27 mCi	28/01/2016 11:00 AM	0 mCi	26.99 mCi	28/01/2016 11:00 AM	0 ml	555 ml	0 ml

Figure 3.5-4 Edit Batch Details

**STEP 5**

Click on **Edit Details** quantity to allocate quantity

**RADIOPHARMACEUTICALS VERIFICATION**

Date: 28/01/2016      Receipt No: R116018636

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	200	val	45RP99	02/01/2021	Accepted

**STOCK DETAILS**

Calibration Activity: 999 MBq \* 27 mCi \*      Calibration Date/Time: 28/01/2016 11:00 AM \*  
 Expected Activity: 999 MBq \* 27 mCi \*      Expiry Date/Time: 20/12/2016 6:03 PM \*  
 Measured Activity: 998.56 MBq \* 26.99 mCi \*      Measured Date/Time: 28/01/2016 11:00 AM \*  
 Activity Variation: -0.04 %      Volume: 555 ml \*

**MEASUREMENT DETAILS**

Label Packaging Type: White I [Guideline](#) **6**  
 Survey on Surface: 0  
 Survey at 1 Meter Distance: 0  
 Wipe Test: 0  
 Remarks:

**VERIFICATION DETAILS**

Accept       Reject  
 Reject Reason:  
 Done By: Pharmacist Radio      Date/Time: 28/01/2016 11:00 AM  
 Updated By: MAIKAL FARDHA      Updated Date/Time: 24/05/2016 9:59 PM

Figure 3.5-5 Sample of Item Batch Details

**STEP 6**

Click on guideline link to get details about measurement details


**Note**

- Label Packaging Type Guidelines as per Figure 6.0-6
- Verification can be Accept and Reject based verification

**LABEL PACKAGING TYPE GUIDELINES**

Unit	Label Packaging Type	Description
mSv/hr	White I	Survey on Surface = Not more than (or equal) 0.005 mSv/hr
	Yellow II	Survey on Surface = More than 0.005 mSv/hr, not more than (or equal) 0.5 mSv/hr
	Yellow III	Survey on Surface = More than 0.5 mSv/hr, not more than (or equal) 2 mSv/hr
	Yellow III(a)	Survey on Surface = More than 2 mSv/hr, not more than (or equal) 10 mSv/hr
mR/hr	White I	White I : Survey on Surface = Not more than (or equal) 0.5 mR/hr
	Yellow II	Yellow II : Survey on Surface = More than 0.5 mR/hr, not more than (or equal) 50 mR/hr
	Yellow III	Yellow III : Survey on Surface = More than 50 mR/hr, not more than (or equal) 200 mR/hr
	Yellow III(a)	Yellow III(a) : Survey on Surface = More than 200 mR/hr, not more than (or equal) 1000 mR/hr

Figure 3.5-6 Label Packaging Type Guidelines

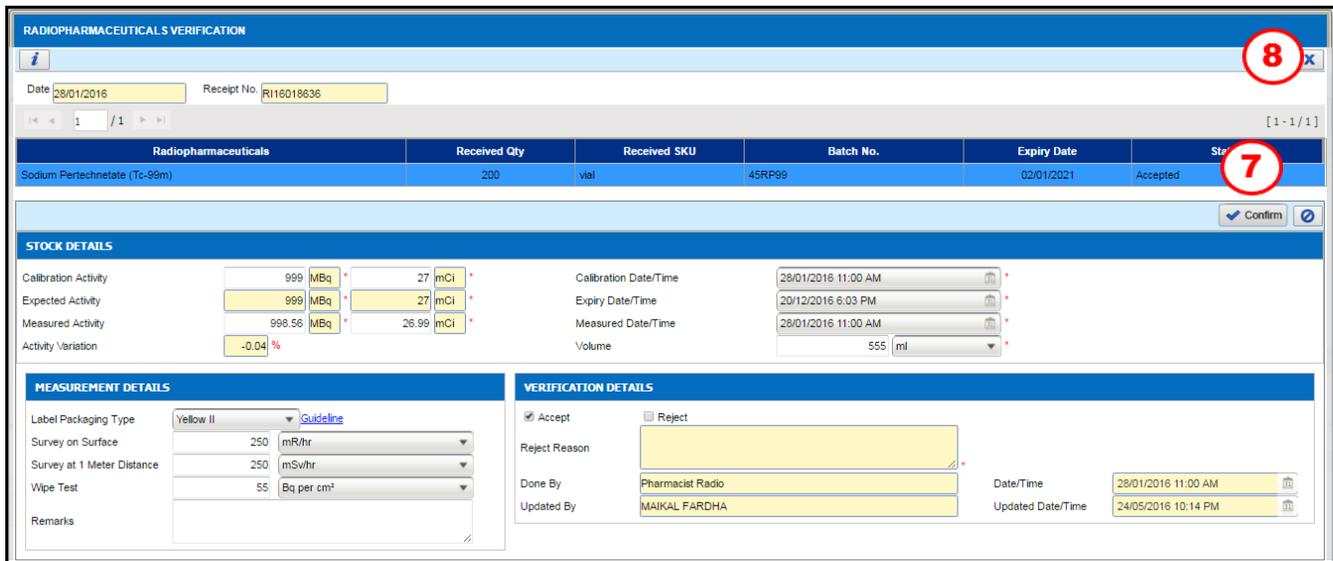

**RADIOPHARMACEUTICALS VERIFICATION**

Date: 28/01/2016 Receipt No: R116018636

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	200	vial	45RP99	02/01/2021	Accepted

**STOCK DETAILS**

Calibration Activity: 999 MBq, 27 mCi  
 Expected Activity: 999 MBq, 27 mCi  
 Measured Activity: 998.56 MBq, 26.99 mCi  
 Activity Variation: -0.04 %

Calibration Date/Time: 28/01/2016 11:00 AM  
 Expiry Date/Time: 20/12/2016 6:03 PM  
 Measured Date/Time: 28/01/2016 11:00 AM  
 Volume: 555 ml

**MEASUREMENT DETAILS**

Label Packaging Type: Yellow II  
 Survey on Surface: 250 mR/hr  
 Survey at 1 Meter Distance: 250 mSv/hr  
 Wipe Test: 55 Bq per cm<sup>2</sup>

**VERIFICATION DETAILS**

Accept  Reject  
 Done By: Pharmacist Radio  
 Updated By: MAIKAL FARDHA  
 Date/Time: 28/01/2016 11:00 AM  
 Updated Date/Time: 24/05/2016 10:14 PM

Figure 3.5-7 Batch Details

**STEP 7**

Click on  button to finalize batch details

**STEP 8**

Click on  button to close Item Batch Details screen



Encounter No: 0632 Admission Date: 19/01/2016 Attending Practitioner: Visit Location: RADIOLOGI Rx No: RADIO0000180122 Order By: Pharmacist Radio

**PRESCRIPTION**

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Sodium Pertechnetate (Tc-99m)	App.Date/Time: 05/05/2016 13:00PM Scan/Therapy: Dacryoscintigraphy	Details	15 mCi	25 mCi (05/05/2016 11:29 AM)	Remarks	Verified

Figure 3.5-8 Radiopharmaceuticals Preparation

**Note**

- User has the option to measure radiopharmaceuticals radioactivity in preparation screen, user need to click the calculate activity button.
- Click on the  button and new window will be displayed as figure 3.4-6.

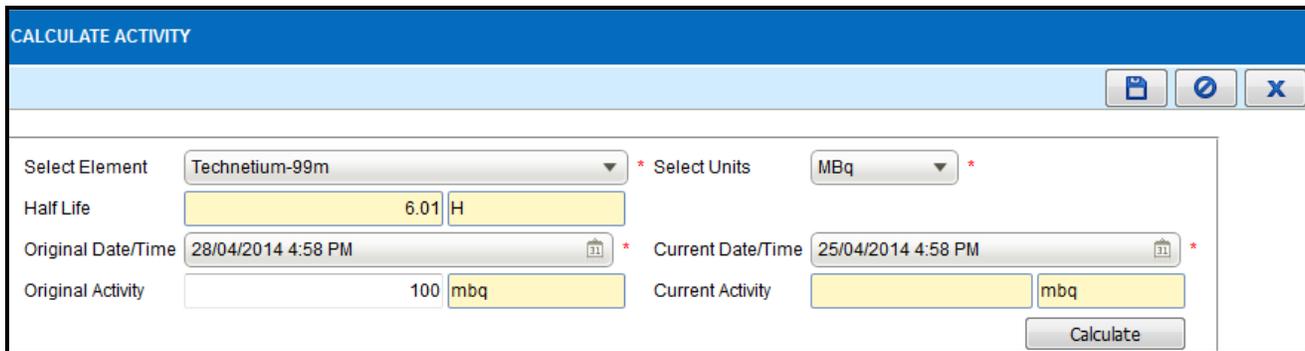


Figure 3.5-9 Radiopharmaceuticals Preparation

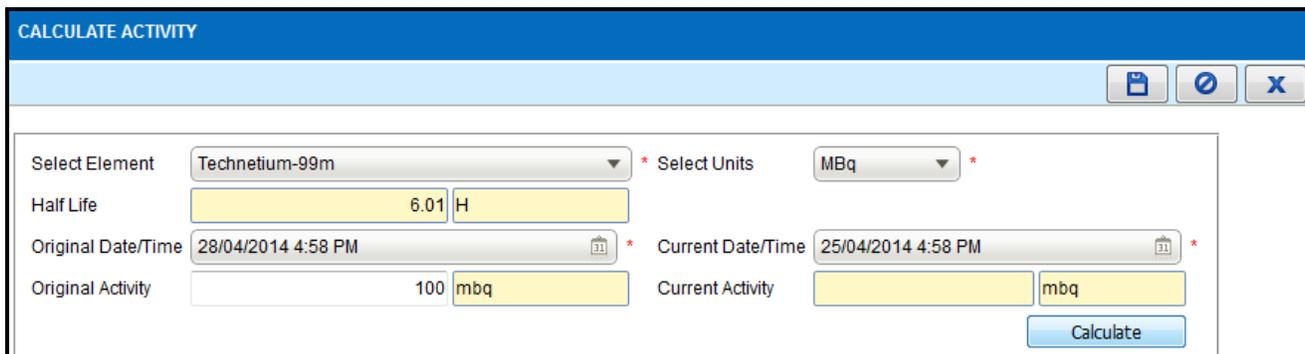
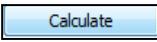


Figure 3.5-10 Calculate Activity

**Note**

- Select value for **Select Element** drop down box.
- **Half Life** field will be automatically displayed by system after user selects value in **Select Element** field.
- **Original Date/Time** is mandatory field.
- **Select Units** from drop down box.
- **Select Units** is mandatory field.
- Select **Current Date/Time** field.
- **Current Date/Time** is mandatory field.
- **Original Activity** field will automatically display value based on the ordered activity but user is allowed to edit the value.
- Click on the  button and **Current Activity** field will be automatically calculated by system as per Figure 3.5-11.

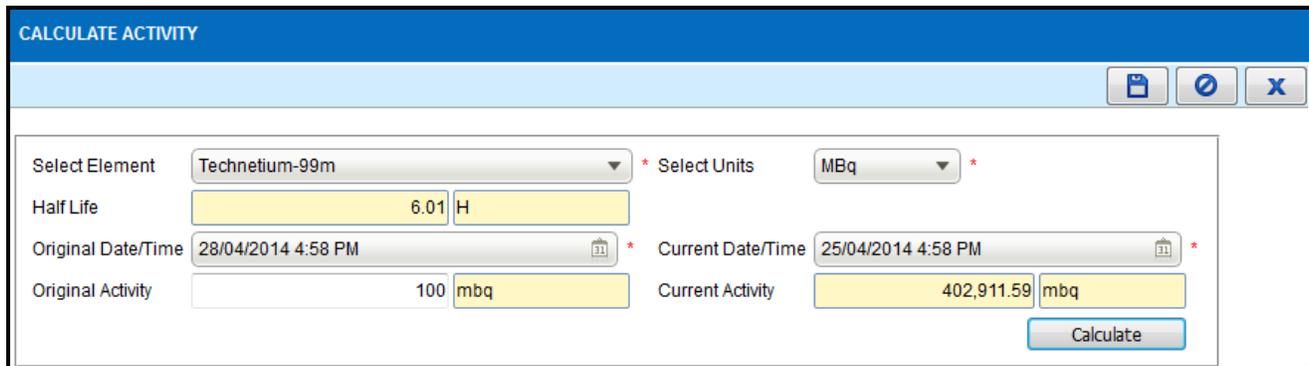


Figure 3.5-11 Calculate Activity

**Note**

- Click on the  button to save record as per Figure 3.5-12.

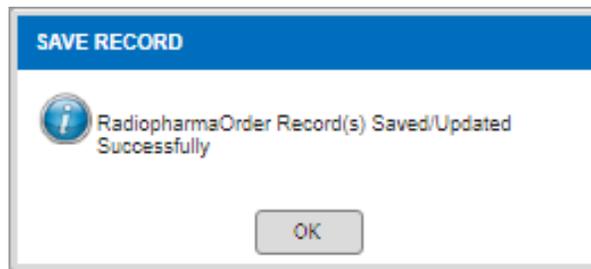
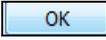


Figure 3.5-13 Save Record Alert Message

**Note**

- Click on the  button to save record as per Figure 3.5-12.

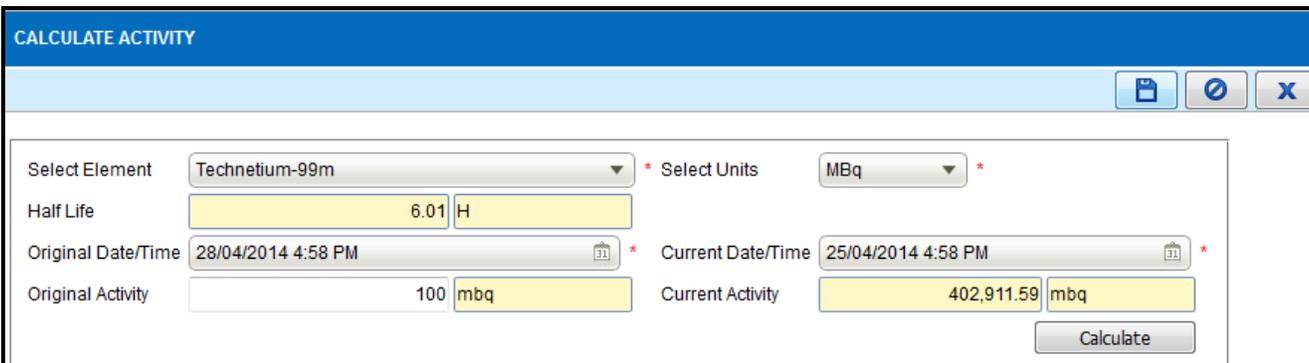


Figure 3.5-14 Calculate Activity

**Note**

Click on the  button to close Calculate Activity screen as per Figure 3.5-14.



Encounter No: 0632 Admission Date: 19/01/2016 Attending Practitioner: Visit Location: RADIOLOGI Rx No: RADIO0000180122 Order By: Pharmacist Radio

Lab Parameter | Counseling History | Demographic | Pharmacist Notes | Medication Profile | Other Appointment


PRESCRIPTION Prepare Label Online Outsource

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Sodium Pertechnetate (Tc-99m)	App.Date/Time: 05/05/2016 13:00PM Scan/Therapy: Dacryoscintigraphy	Details	15 mCi	25 mCi (05/05/2016 11:29 AM)	Remarks	Verified

Figure 3.5-15 Radiopharmaceuticals Preparation

**Note**

- Click on the  button to edit label as per Figure 3.5-16.
- Drug Label screen will be displayed as per Figure 3.5-17.

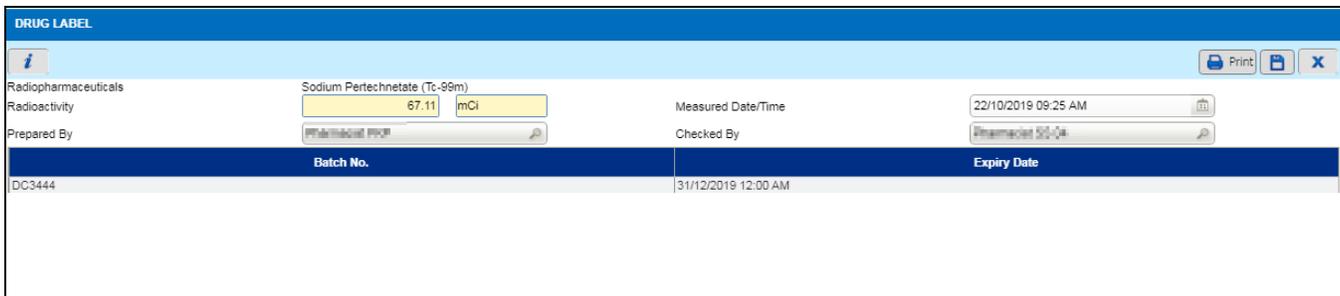

DRUG LABEL Print Save Close

Radiopharmaceuticals Sodium Pertechnetate (Tc-99m)

Radioactivity  mCi Measured Date/Time


Prepared By  Checked By

Batch No.	Expiry Date
DC3444	31/12/2019 12:00 AM

Figure 3.5-16 Drug Label


**Note**


- Radiopharmaceuticals, Batch No., Expiry Date and Administered field will be automatically displayed by system based on the selected drug batch.
- Enter **Radioactivity** amount. **Radioactivity UOM** will be automatically displayed based on the drug SKU.
- Select **Measured Date/Time** field.
- Select **Prepared By** and **Checked By** field.
- Click on the  button to save record as per Figure 3.5-16.

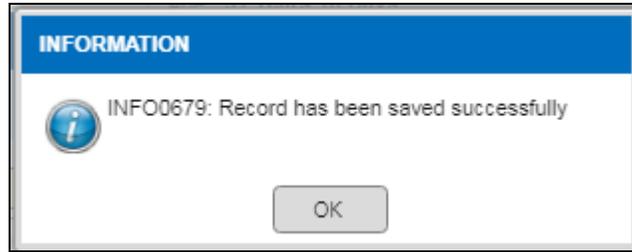

**Figure 3.5-17 Alert Message**

**Note**

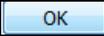
- Click on the  button to save updated Radiopharmaceuticals Order Record as per Figure 3.5-18.

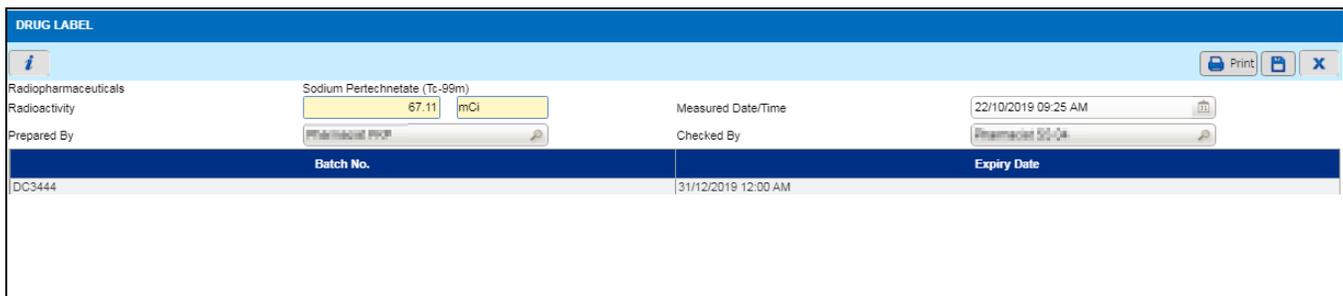

**Figure 3.5-18 Drug Label**

**Note**


- Click on the  button to print label as per Figure 3.4-14 and label will be displayed as per Figure 3.4-18.

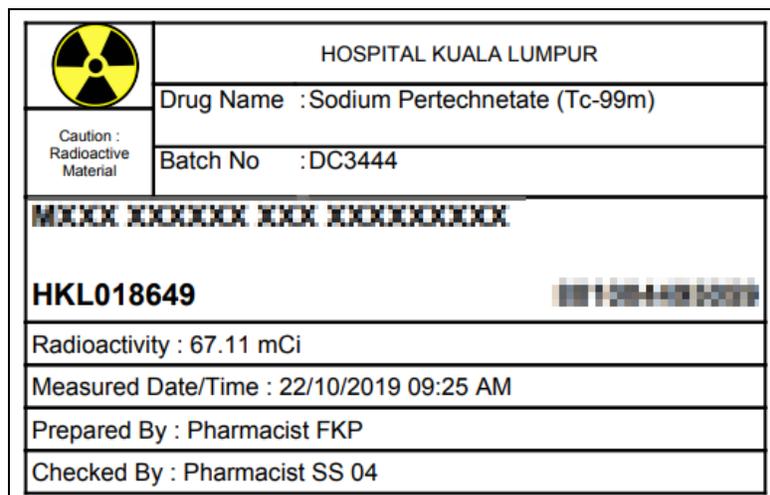

**Figure 3.5-19 Sample Radiopharmaceuticals Label**

**Note**

Sample of Radiopharmaceuticals Label as per Figure 3.5-19.

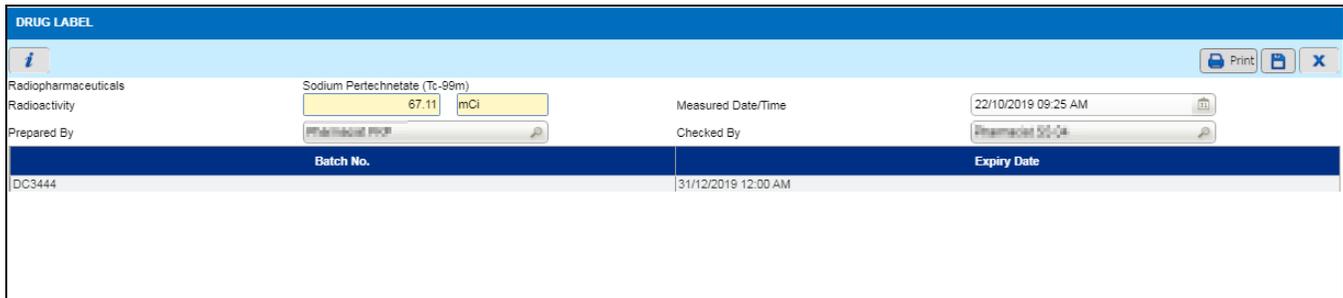

Batch No.	Expiry Date
DC3444	31/12/2019 12:00 AM

Figure 3.5-20 Drug Label

**Note**


Click on the  button to close Drug Label screen as per Figure 3.5-20.



Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Sodium Pertechnetate (Tc-99m)	App.Date/Time: 05/05/2016 13:00PM Scan/Therapy; Dacryoscintigraphy	Details	15 mCi	25 mCi (05/05/2016 11:29 AM)	Remarks	Verified

Figure 3.5-21 Radiopharmaceuticals Preparation

**STEP 8** Click on the  button to prepare Radiopharmaceuticals order and alert message will be displayed as per Figure 3.5-21

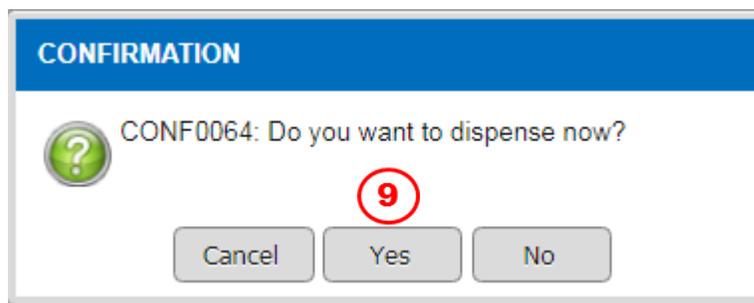

**CONFIRMATION**


CONF0064: Do you want to dispense now?

Buttons: Cancel, Yes, No

Figure 3.5-22 Ordered Prepared Alert Message

**STEP 9**

Click on the  button

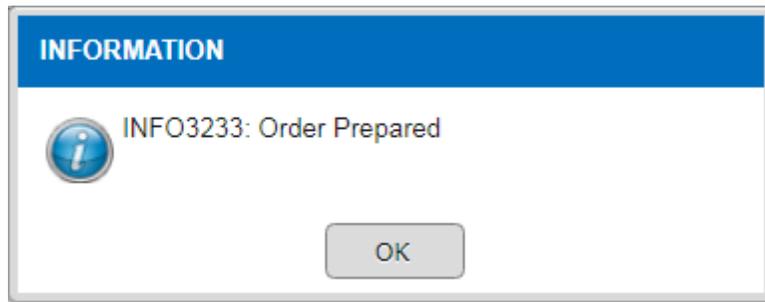


**Figure 3.5-23 Ordered Prepared Alert Message**

**Note**

- Click on the  button and order status will be updated to Prepared
- *User still allowed doing intervention and holding order in preparation stage.*

### 3.6 Radiopharmaceutical Dispensing

Radiopharmaceuticals dispensing screen allow user to dispense medication to patient, reprint drug label and record collection details.

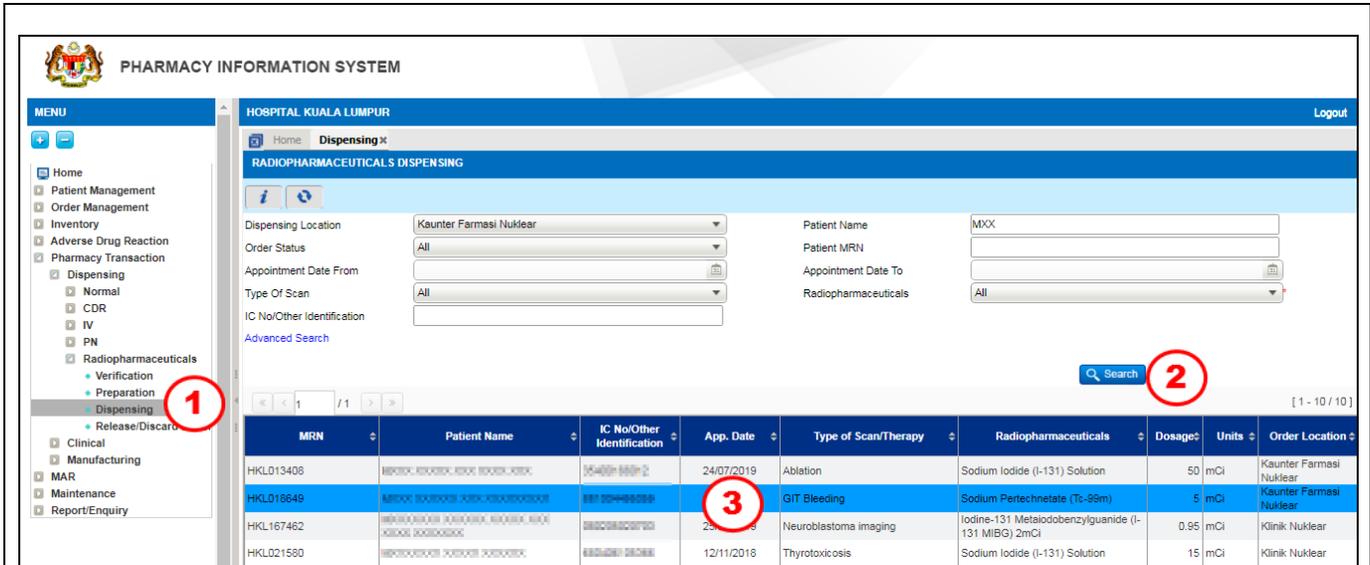

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Logout

Home Dispensing

**RADIOPHARMACEUTICALS DISPENSING**

Dispensing Location: Kaunter Farmasi Nuklear  
 Order Status: All  
 Appointment Date From:   
 Appointment Date To:   
 Type Of Scan: All  
 IC No/Other Identification:   
 Patient Name: MXX  
 Patient MRN:   
 Radiopharmaceuticals: All

Advanced Search

Search

MRN	Patient Name	IC No/Other Identification	App. Date	Type of Scan/Therapy	Radiopharmaceuticals	Dosage	Units	Order Location
HKLD13408	XXXXXXXXXX XXXX XXXX	XXXXXXXXXX	24/07/2019	Ablation	Sodium Iodide (I-131) Solution	50	mCi	Kaunter Farmasi Nuklear
HKLD18649	XXXXXXXXXX XXXX XXXX	XXXXXXXXXX	24/07/2019	GIT Bleeding	Sodium Pertechnetate (Tc-99m)	5	mCi	Kaunter Farmasi Nuklear
HKLD167462	XXXXXXXXXX XXXX XXXX	XXXXXXXXXX	24/07/2019	Neuroblastoma imaging	Iodine-131 Metaiodobenzylguanide (I-131 MIBG) 2mCi	0.95	mCi	Klinik Nuklear
HKLD021580	XXXXXXXXXX XXXX XXXX	XXXXXXXXXX	12/11/2018	Thyrotoxicosis	Sodium Iodide (I-131) Solution	15	mCi	Klinik Nuklear

Figure 3.6-1 Radiopharmaceuticals Dispensing Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu follow by 'Dispensing' and click on 'Radiopharmaceuticals Dispensing' sub-menu

#### Note


- Various search criteria is provided as below: Dispense Location
  - a) Patient Name
  - b) Order Status
  - c) Patient MRN
  - d) Appointment Date From
  - e) Appointment Date To
  - f) Type of Scan
  - g) Radiopharmaceuticals
- Click on the [Advanced Search](#) hyperlink for advance search.

#### STEP 2

Click on the [Search](#) button and system will be displayed related order

#### STEP 3

Double click on the selected patient record and system will be displayed Radiopharmaceuticals Dispensing screen as per Figure 3.6-1

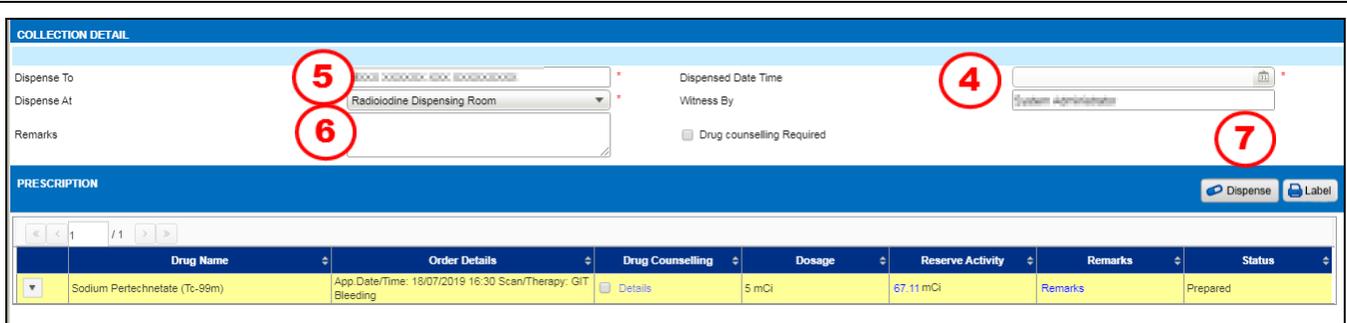

**COLLECTION DETAIL**

Dispense To: [5] [XXXXXXXXXX XXXXXXXXXXXXX] \*

Dispense At: [6] Radiiodine Dispensing Room \*

Remarks: [7]

Dispensed Date Time: [4] [ ] \*

Witness By: [ ] \*

Drug counselling Required

**PRESCRIPTION** [Dispense] [Label]

Drug Name	Order Details	Drug Counselling	Dosage	Reserve Activity	Remarks	Status
Sodium Pertechnetate (Tc-99m)	App Date/Time: 15/07/2019 16:30 Scan/Therapy: GIT Bleeding	<input type="checkbox"/> Details	5 mCi	67.11 mCi	Remarks	Prepared

Figure 3.6-2 Radiopharmaceuticals

**STEP 4**

Select **Dispensed Date Time** field.

**STEP 5**

Select **Dispense At** from drop down box:

- Administration Room
- Radioiodine Dispensing Room

**STEP 6**


Enter **Remarks** field

**Note**

Check on the  **Drug Counseling Required** checkbox if counseling required.

**STEP 7**

Click on the  button to dispense Radiopharmaceuticals and Collection Detail screen will be displayed as per Figure 3.6-2

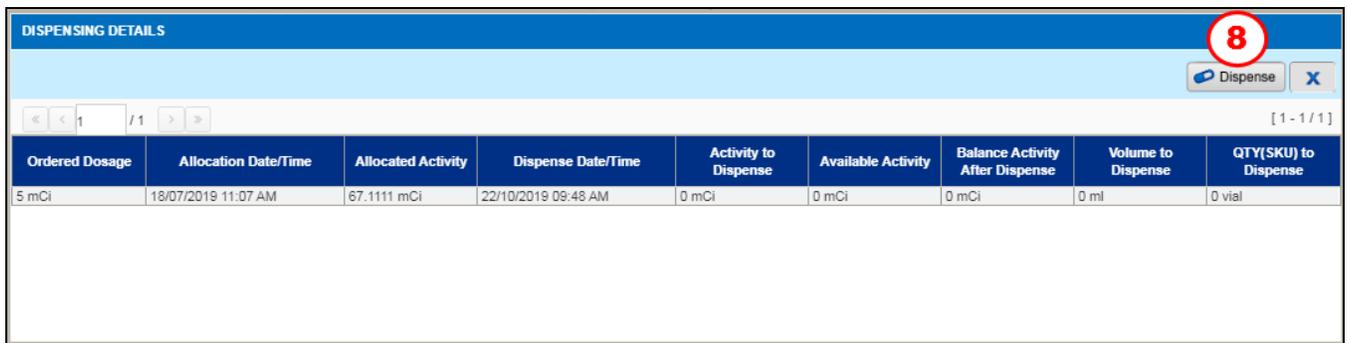

**DISPENSING DETAILS** [8] [Dispense] [X]


[1 - 1 / 1]

Ordered Dosage	Allocation Date/Time	Allocated Activity	Dispense Date/Time	Activity to Dispense	Available Activity	Balance Activity After Dispense	Volume to Dispense	QTY(SKU) to Dispense
5 mCi	18/07/2019 11:07 AM	67.1111 mCi	22/10/2019 09:48 AM	0 mCi	0 mCi	0 mCi	0 ml	0 vial

Figure 3.6-3 Dispensing Details Screen

**STEP 8**

Click on the  button to dispense Radiopharmaceuticals record as per Figure 3.6-3.

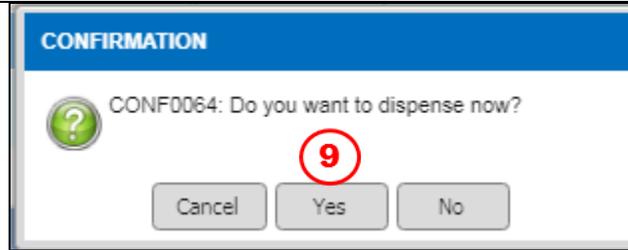


**Figure 3.6-4 Information Alert Message**

**STEP 9**

Click on the  button to proceed dispensing

### 3.7 Release Batch

Radiopharmaceuticals Release Batch Screen will allow user to release any reserved drug for the patient back into inventory.

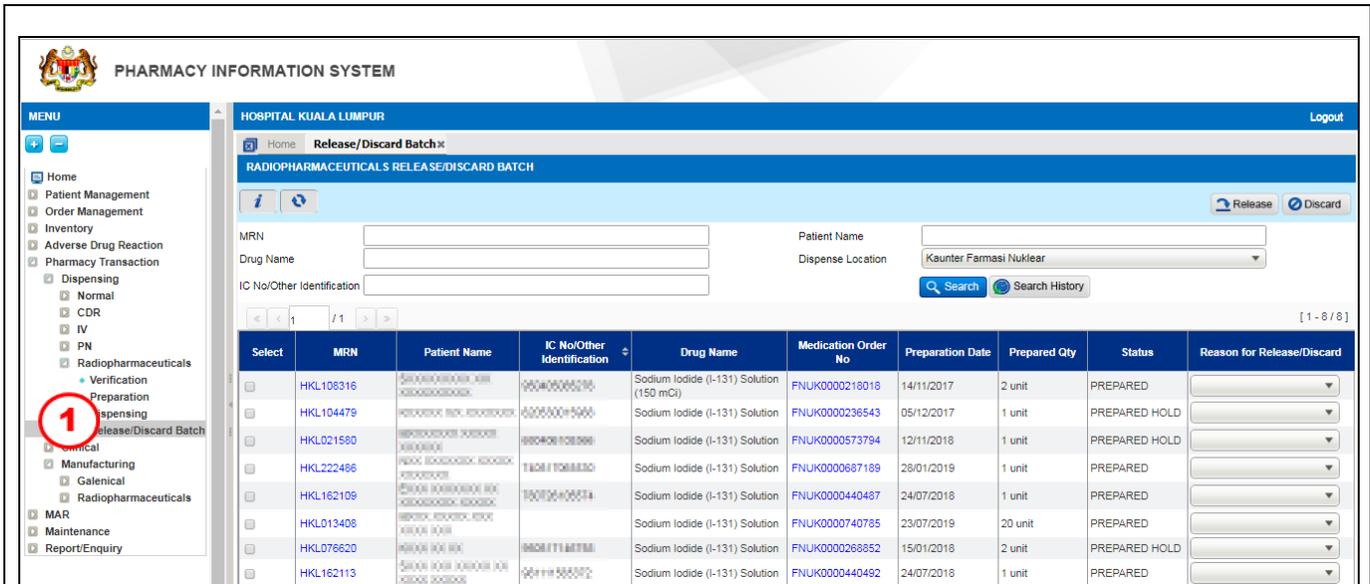

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Home Release/Discard Batch

RADIOPHARMACEUTICALS RELEASE/DISCARD BATCH

MRN:  Patient Name:   
 Drug Name:  Dispense Location: Kaunter Farmasi Nuklear  
 IC No/Other Identification:  Search Search History

Select	MRN	Patient Name	IC No/Other Identification	Drug Name	Medication Order No	Preparation Date	Prepared Qty	Status	Reason for Release/Discard
<input type="checkbox"/>	HKL108316	SRIKANDIAH, ANN	99040000290	Sodium Iodide (I-131) Solution (150 mCi)	FNUK0000218018	14/11/2017	2 unit	PREPARED	
<input type="checkbox"/>	HKL104479	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000236543	05/12/2017	1 unit	PREPARED HOLD	
<input type="checkbox"/>	HKL021580	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000573794	12/11/2018	1 unit	PREPARED HOLD	
<input type="checkbox"/>	HKL222486	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000687189	28/01/2019	1 unit	PREPARED	
<input type="checkbox"/>	HKL162109	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000440487	24/07/2018	1 unit	PREPARED	
<input type="checkbox"/>	HKL013408	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000740785	23/07/2019	20 unit	PREPARED	
<input type="checkbox"/>	HKL078620	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000268852	15/01/2018	2 unit	PREPARED HOLD	
<input type="checkbox"/>	HKL162113	ABDULLAH, MOHAMMAD	00000000000	Sodium Iodide (I-131) Solution	FNUK0000440492	24/07/2018	1 unit	PREPARED	

Figure 3.7-1 Radiopharmaceuticals Release Batch Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu follow by 'Dispensing' and click on 'Radiopharmaceuticals Release Batch' sub-menu

#### Note

Various search criteria is provided as below:

- MRN
- Patient Name
- Drug Name
- Dispense Location
- IC No/Other Identification

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Home Release/Discard Batch

RADIOPHARMACEUTICALS RELEASE/DISCARD BATCH

MRN: [ ] Patient Name: [ ]  
 Drug Name: [ ] Dispense Location: Kaunter Farmasi Nuklear  
 IC No/Other Identification: [ ]

[ Search ] [ Search History ] [ 1 - 8 / 8 ]

Select	MRN	Patient Name	IC No/Other Identification	Drug Name	Medication Order No	Preparation Date	Prepared Qty	Status	Reason for Release/Discard
<input checked="" type="checkbox"/>	HKL108316	[ ]	[ ]	Sodium Iodide (I-131) Solution (150 mCl)	FNUK0000218018	14/11/2017	2 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL104479	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000236543	05/12/2017	1 unit	PREPARED HOLD	Change of Regimen Faulty Product Intervention Patient Refused Uncollected Wrong Allocation
<input type="checkbox"/>	HKL021580	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000573794	12/11/2018	1 unit	PREPARED HOLD	[ ]
<input type="checkbox"/>	HKL222486	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000687189	28/01/2019	1 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL162109	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000440487	24/07/2018	1 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL013408	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000740785	23/07/2019	20 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL078620	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000268852	15/01/2018	2 unit	PREPARED HOLD	[ ]
<input type="checkbox"/>	HKL162113	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000440492	24/07/2018	1 unit	PREPARED	[ ]

**Figure 3.7-2 Radiopharmaceuticals Release Batch**

**STEP 2**

Select **Reason for Release** from drop down box:

- Uncollected
- Wrong Allocation
- Change of Regimen
- Patient Refused
- Intervention
- Faulty Product

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Home Release/Discard Batch

RADIOPHARMACEUTICALS RELEASE/DISCARD BATCH

MRN: [ ] Patient Name: [ ]  
 Drug Name: [ ] Dispense Location: Kaunter Farmasi Nuklear  
 IC No/Other Identification: [ ]

[ Search ] [ Search History ] [ 1 - 8 / 8 ]


Select	MRN	Patient Name	IC No/Other Identification	Drug Name	Medication Order No	Preparation Date	Prepared Qty	Status	Reason for Release/Discard
<input checked="" type="checkbox"/>	HKL108316	[ ]	[ ]	Sodium Iodide (I-131) Solution (150 mCl)	FNUK0000218018	14/11/2017	2 unit	PREPARED	Change of Regimen
<input type="checkbox"/>	HKL104479	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000236543	05/12/2017	1 unit	PREPARED HOLD	[ ]
<input type="checkbox"/>	HKL021580	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000573794	12/11/2018	1 unit	PREPARED HOLD	[ ]
<input type="checkbox"/>	HKL222486	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000687189	28/01/2019	1 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL162109	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000440487	24/07/2018	1 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL013408	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000740785	23/07/2019	20 unit	PREPARED	[ ]
<input type="checkbox"/>	HKL078620	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000268852	15/01/2018	2 unit	PREPARED HOLD	[ ]
<input type="checkbox"/>	HKL162113	[ ]	[ ]	Sodium Iodide (I-131) Solution	FNUK0000440492	24/07/2018	1 unit	PREPARED	[ ]

**Figure 3.7-3 Radiopharmaceuticals Release Batch**

**STEP 3**


Check the  **Select** checkbox

#### STEP 4

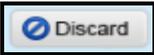Click on the  button to release batch and  button to discard

#### **Note**

- Release function is to release and update back stock into inventory.
- Discard function is to discard stock from inventory even not doing dispensing.

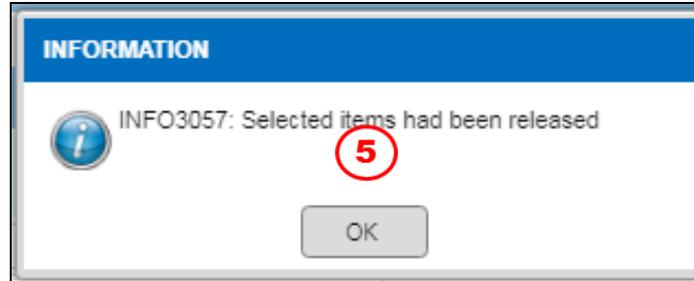

**Figure 3.7-4 Information Alert Message**

#### STEP 5

Click on the  button to release selected items


#### **Note**

*Released order will flow back to Preparation screen for user to allocate new batch quantity.*

## 3.8 Radiopharmaceutical Manufacturing

### 3.8.1 Generate Elution

Generator Elution screen allow user to create Technetium (99mTc) Sodium Pertechnetate which can be used to manufacture Kit based Radiopharmaceuticals in kit Preparation screen.

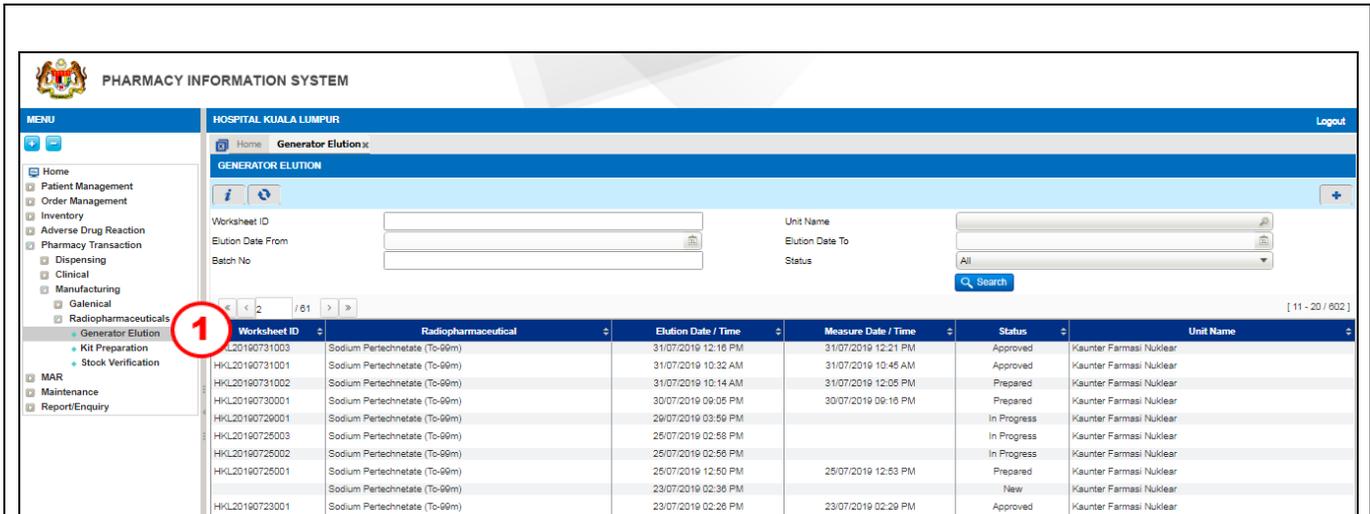

Worksheet ID	Radiopharmaceutical	Elution Date / Time	Measure Date / Time	Status	Unit Name
HLK20190731003	Sodium Pertechnetate (Tl-99m)	31/07/2019 12:16 PM	31/07/2019 12:21 PM	Approved	Kaunter Farmasi Nuklear
HLK20190731001	Sodium Pertechnetate (Tl-99m)	31/07/2019 10:32 AM	31/07/2019 10:45 AM	Approved	Kaunter Farmasi Nuklear
HLK20190731002	Sodium Pertechnetate (Tl-99m)	31/07/2019 10:14 AM	31/07/2019 12:05 PM	Prepared	Kaunter Farmasi Nuklear
HLK20190730001	Sodium Pertechnetate (Tl-99m)	30/07/2019 09:05 PM	30/07/2019 09:16 PM	Prepared	Kaunter Farmasi Nuklear
HLK20190729001	Sodium Pertechnetate (Tl-99m)	29/07/2019 03:59 PM		In Progress	Kaunter Farmasi Nuklear
HLK20190725003	Sodium Pertechnetate (Tl-99m)	25/07/2019 02:58 PM		In Progress	Kaunter Farmasi Nuklear
HLK20190725002	Sodium Pertechnetate (Tl-99m)	25/07/2019 02:58 PM		In Progress	Kaunter Farmasi Nuklear
HLK20190725001	Sodium Pertechnetate (Tl-99m)	25/07/2019 12:50 PM	25/07/2019 12:53 PM	Prepared	Kaunter Farmasi Nuklear
HLK20190723001	Sodium Pertechnetate (Tl-99m)	23/07/2019 02:38 PM		New	Kaunter Farmasi Nuklear
HLK20190723001	Sodium Pertechnetate (Tl-99m)	23/07/2019 02:28 PM	23/07/2019 02:29 PM	Approved	Kaunter Farmasi Nuklear

Figure 3.8.1-1 Generator Elution Listing Page

#### STEP 1

Click on 'Pharmacy Transaction' menu followed by 'Manufacturing' and click on 'Radiopharmaceuticals Generator Elution' sub-menu

#### Note

- Various search criteria is provided as below:
  - Worksheet ID
  - Date From
  - Date To
  - Batch No
  - Status

Click on the  button and system will display related record.

PHARMACY INFORMATION SYSTEM

HOSPITAL KUALA LUMPUR

Logout

Home Generator Elution

GENERATOR ELUTION

Worksheet ID: [ ] Elution Date From: [ ] Elution Date To: [ ] Unit Name: [ ] Status: [ All ]

Batch No: [ ] [ Search ]

[ 11 - 20 / 602 ]

Worksheet ID	Radiopharmaceutical	Elution Date / Time	Measure Date / Time	Status	Unit Name
HKL20190731003	Sodium Pertechnetate (Tc-99m)	31/07/2019 12:18 PM	31/07/2019 12:21 PM	Approved	Kaunter Farmasi Nuklear
HKL20190731001	Sodium Pertechnetate (Tc-99m)	31/07/2019 10:32 AM	31/07/2019 10:45 AM	Approved	Kaunter Farmasi Nuklear
HKL20190731002	Sodium Pertechnetate (Tc-99m)	31/07/2019 10:14 AM	31/07/2019 12:05 PM	Prepared	Kaunter Farmasi Nuklear
HKL20190730001	Sodium Pertechnetate (Tc-99m)	30/07/2019 09:05 PM	30/07/2019 09:18 PM	Prepared	Kaunter Farmasi Nuklear
HKL20190729001	Sodium Pertechnetate (Tc-99m)	29/07/2019 03:59 PM		In Progress	Kaunter Farmasi Nuklear
HKL20190725003	Sodium Pertechnetate (Tc-99m)	25/07/2019 02:58 PM		In Progress	Kaunter Farmasi Nuklear
HKL20190725002	Sodium Pertechnetate (Tc-99m)	25/07/2019 02:56 PM		In Progress	Kaunter Farmasi Nuklear
HKL20190725001	Sodium Pertechnetate (Tc-99m)	25/07/2019 12:50 PM	25/07/2019 12:53 PM	Prepared	Kaunter Farmasi Nuklear
HKL20190723001	Sodium Pertechnetate (Tc-99m)	23/07/2019 02:38 PM		New	Kaunter Farmasi Nuklear
HKL20190723001	Sodium Pertechnetate (Tc-99m)	23/07/2019 02:26 PM	23/07/2019 02:29 PM	Approved	Kaunter Farmasi Nuklear

Figure 3.8.1-2 Generator Elution Listing Page

STEP 2

Click on button to create new generator elution record and system will be automatically displayed screen as per Figure 3.8.1-3

Note

- System will be automatically displayed Preparation screen as per Figure 3.8.1-3.
- The **Product** field will be automatically displayed by system.

GENERATOR ELUTION

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 29/07/2019 03:56:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: In Progress

Preparation

ITEM BATCH DETAILS

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand	Volume	Allocated Qty
Generator	Tc-99m Sterile Generator 1000 mCi (HKL)	222	31/01/2020 00:00 AM	DRYTEC/GE Healthcare		
Elution Vial	Elution Vial DRYTEC@@	4	31/12/2019	DRYTEC		1 EA
Saline Vial	Saline Vial DRYTEC (10ml)@@		31/12/2019	DRYTEC	10 ml*	2 EA

PREPARATION DETAILS

Worksheet Prepared By: [ ]  
Worksheet Checked By: [ ]  
Remarks: [ ]  
Storage: Keep below 25 °C

Measured Radioactivity: 1 mCi  
Radioactivity Concentration: 0.1 mCi/ml  
Measured Date/Time: 12/10/2019 3:31 PM  
Product Batch No.: [ ]  
Expiry Date: 25/10/2019 3:31 PM

Figure 3.8.1-3 Generator Elution

STEP 3


Select **Unit Name**, **Generator** from drop down box and **Batch Number** will be automatically displayed based on current available stock

STEP 4

Select **Elution Vial** from drop down box and **Batch Number** will be automatically displayed based on current available stock

**STEP 5**

Select **Saline Vial** from drop down box and **Batch Number** will be automatically displayed base on current available stock

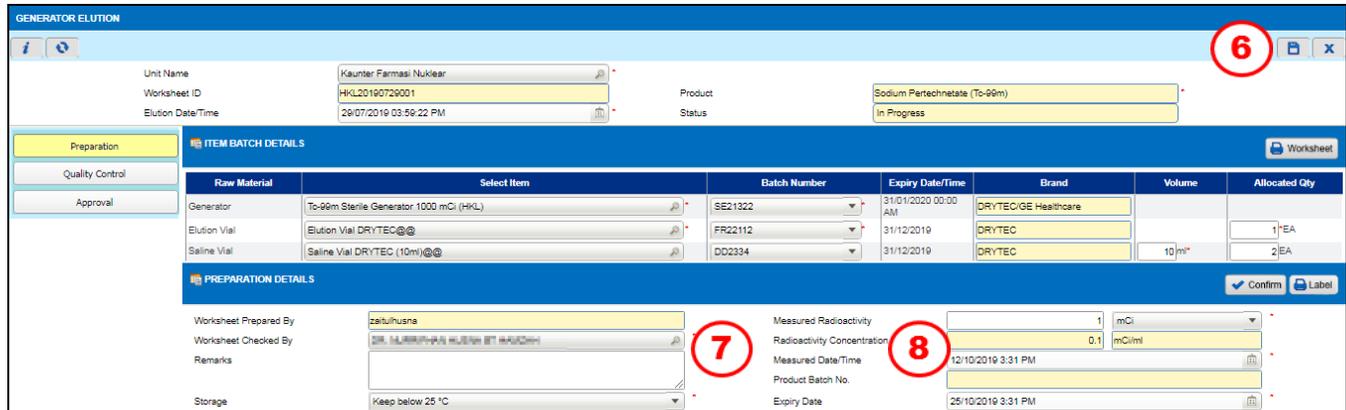

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 29/07/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: In Progress

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand	Volume	Allocated Qty
Generator	To-99m Sterile Generator 1000 mCi (HKL)	SE21322	31/01/2020 00:00 AM	DRYTEC-GE Healthcare		
Elution Vial	Elution Vial DRYTEC	FR22112	31/12/2019	DRYTEC		1 EA
Saline Vial	Saline Vial DRYTEC (10ml)	DD2334	31/12/2019	DRYTEC	10 ml*	2 EA

**PREPARATION DETAILS**

Worksheet Prepared By: zafuhusna  
Worksheet Checked By: DR. MURRAY HUSAIN BT HAJIDIN  
Remarks:  
Storage: Keep below 25 °C

Measured Radioactivity: 1 mCi  
Radioactivity Concentration: 0.1 mCi/ml  
Measured Date/Time: 12/10/2019 3:31 PM  
Product Batch No.:  
Expiry Date: 25/10/2019 3:31 PM

**Figure 3.8.1-4 Generator Elution**

**STEP 6**

Click on the  button to save record

**STEP 7**

Enter **Measured Radioactivity** field and **Radioactivity Concentration** will be automatically calculated and displayed as per Figure 3.8.1-4

**Note**

- Enter **Volume** and **Allocated Qty** for Elution Vial and Saline Vial.
- **Status** will changed to **New** after saving.
- The  button will automatically activate after saving.

**STEP 8**

Enter **Measured Date/Time** field

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 29/07/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: In Progress

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand	Volume	Allocated Qty
Generator	To-99m Sterile Generator 1000 mCi (HKL)	SE21322	31/01/2020 00:00 AM	DRYTEC/GE Healthcare		1 EA
Elution Vial	Elution Vial DRYTEC	FR22112	31/12/2019	DRYTEC		2 EA
Saline Vial	Saline Vial DRYTEC (10ml)	DD2334	31/12/2019	DRYTEC	10 ml*	

**PREPARATION DETAILS**

Worksheet Prepared By: Zatulhuana  
Worksheet Checked By: DR. NURRIZKIYAH HUSAIN BT HUSAINI  
Remarks:  
Storage: Keep below 25 °C

Measured Radioactivity: 1 mCi  
Radioactivity Concentration: 0.1 mCi/ml  
Measured Date/Time: 12/10/2019 3:31 PM  
Product Batch No.:  
Expiry Date: 25/10/2019 3:31 PM

**Figure 3.8.1-5 Generator Elution**


**STEP 9**



Click on the button to generate the worksheet and alert message will be displayed as per Figure 3.8.1-5

**Note**

Enter **Remarks** if applicable.

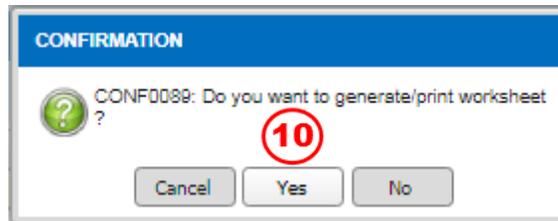

**Figure 3.8.1-6 Generate Worksheet Alert Message**

**STEP 10**

Click on the button to generate/print worksheet and worksheet will be displayed as per Figure 3.8.1-6

RADIOPHARMACEUTICALS WORKSHEET					
MINISTRY OF HEALTH HOSPITAL KUALA LUMPUR					
Generator Elution records for Sodium Pertechnetate (Tc-99m)					
Eluate Batch No:					
Worksheet ID: HKL20190729001					
Preparation of Sodium Pertechnetate (Tc-99m) Eluate					
Elution Date/Time	29/07/2019 03:59 PM		Generator Brand/Manufacturer	Worksheet Prepared By : zali/husna	
Generator Batch No.	SE21322	Exp. Date	31/01/2020 00:00	Signature Name Date/Time	
Elution Vial Batch No.	FR22112	Exp. Date	31/12/2019	Signature Name Date/Time	
Saline Vial Batch No.	DD2334	Exp. Date	31/12/2019	Remarks :	Worksheet Checked By : DR. NURRIFHAN HUSNA BT HAMZAH
Volume	10 ml		Signature Name Date/Time		
Elution BN	Signature Name Date/Time				
Measured Activity	Preparation Done By :				
Measured Date/Time	Signature Name Date/Time				
Radioactivity Concentration	0.1 mCi/ml				
Quality Control of Eluate					
Measure Mo-99 Activity		Limit		Pass/Fail	Quality Control Done By
Quality Control Test	Mo-99 Breakthrough	<0.15 µCi of Mo-99 per mCi of Tc-99m			Signature Name Date/Time
	Eluate pH	4.0-8.0			
	Alumina Breakthrough	Intensity less than Standard			
Approval of Eluate				Label	
Checklist for approval :			Approved By :		
Preparation done accordingly			Signature Name Date/Time		
QC done accordingly and pass					
Label done accordingly					
Remarks	Signature Name Date/Time				

Figure 3.8.1-7 Sample Of Worksheet

**Note**

Sample of worksheet as per Figure 3.7.1-7.

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 29/07/2019 03:59:22 PM

Product: Sodium Pertechnetate (Tc-99m)  
Status: In Progress

Preparation	ITEM BATCH DETAILS						
Quality Control	Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand	Volume	Allocated Qty
Approval	Generator	Tc-99m Sterile Generator 1000 mCi (HKL)	SE21322	31/01/2020 00:00 AM	DRYTEC/GE Healthcare		1 EA
	Elution Vial	Elution Vial DRYTEC	FR22112	31/12/2019	DRYTEC		1 EA
	Saline Vial	Saline Vial DRYTEC (10ml)	DD2334	31/12/2019	DRYTEC	10 ml	2 EA

**PREPARATION DETAILS**

Worksheet Prepared By: zali/husna  
Worksheet Checked By: DR. NURRIFHAN HUSNA BT HAMZAH  
Remarks:   
Storage: Keep below 25 °C


Measured Radioactivity: 1 mCi  
Radioactivity Concentration: 0.1 mCi/ml  
Measured Date/Time: 12/10/2019 3:31 PM  
Product Batch No.:   
Expiry Date: 25/10/2019 3:31 PM

Confirm  Label

Figure 3.8.1-8 Generator Elution

**Note**

- **Worksheet ID** will automatically display after generate worksheet as per Figure 3.8.1-8.
- **Status** will change to *In Progress* after generate worksheet.
- The  Confirm button will activate after generate worksheet.

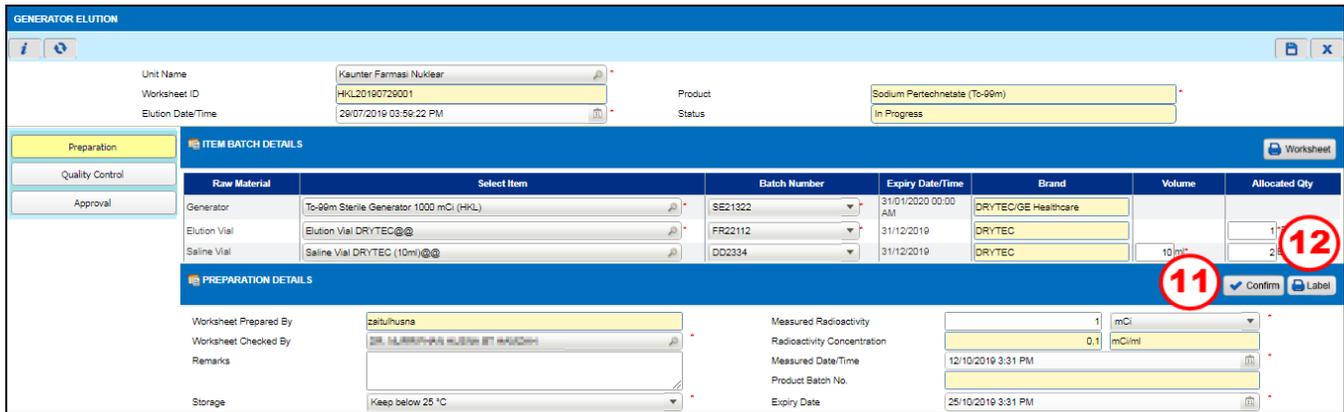

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 29/07/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: In Progress

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand	Volume	Allocated Qty
Generator	To-99m Sterile Generator 1000 mCi (HKL)	SE21322	31/01/2020 00:00 AM	DRYTEC/GE Healthcare		1
Elution Vial	Elution Vial DRYTEC	FR22112	31/12/2019	DRYTEC		1
Saline Vial	Saline Vial DRYTEC (10ml)	DO2234	31/12/2019	DRYTEC	10 ml	2

**PREPARATION DETAILS**

Worksheet Prepared By: zafrihusna  
Worksheet Checked By: DR. HUSRIYATI HUSRIYATI HUSRIYATI  
Remarks:  
Storage: Keep below 25 °C


Measured Radioactivity: 1 mCi  
Radioactivity Concentration: 0.1 mCi/ml  
Measured Date/Time: 12/10/2019 3:31 PM  
Product Batch No.:  
Expiry Date: 25/10/2019 3:31 PM

**11** Confirm Label

Figure 3.8.1-9 Generator Elution

**STEP 11**

Click on the  button to confirm information preparation details

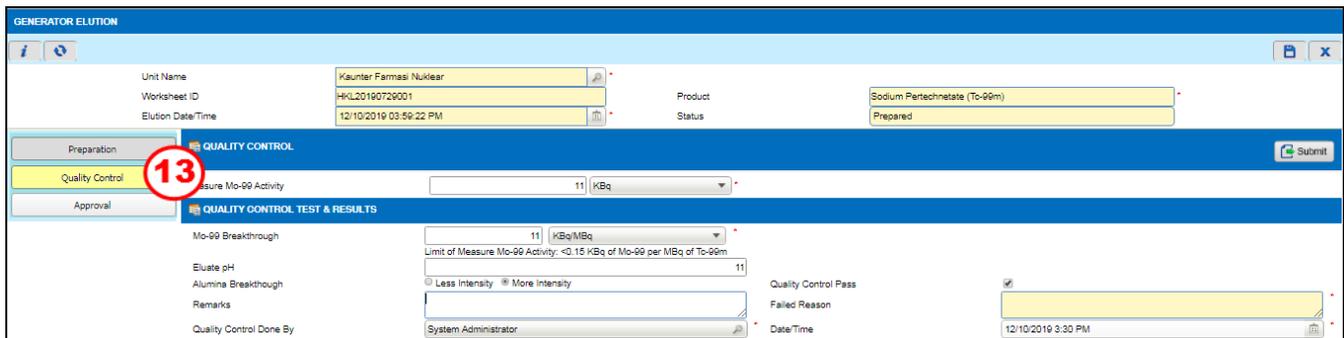

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 12/10/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: Prepared

**QUALITY CONTROL**

Measure Mo-99 Activity: 11 KBq

**QUALITY CONTROL TEST & RESULTS**

Mo-99 Breakthrough: 11 KBq/MBq  
Limit of Measure Mo-99 Activity: <0.15 KBq of Mo-99 per MBq of Tc-99m  
Eluate pH: 11  
Alumina Breakthrough: Less Intensity More Intensity  
Remarks:  
Quality Control Done By: System Administrator  
Quality Control Pass:   
Failed Reason:  
Date/Time: 12/10/2019 3:30 PM

**13** Submit

Figure 3.8.1-10 Generator Elution

**Note**

- The  button will activate after user confirmed the preparation details.
- Status** will change to Prepared and the  button will be activated.

**STEP 12**

Click on the  button to edit label

**STEP 13**

Click on the  button and system will be displayed screen as per Figure 3.8.1-10

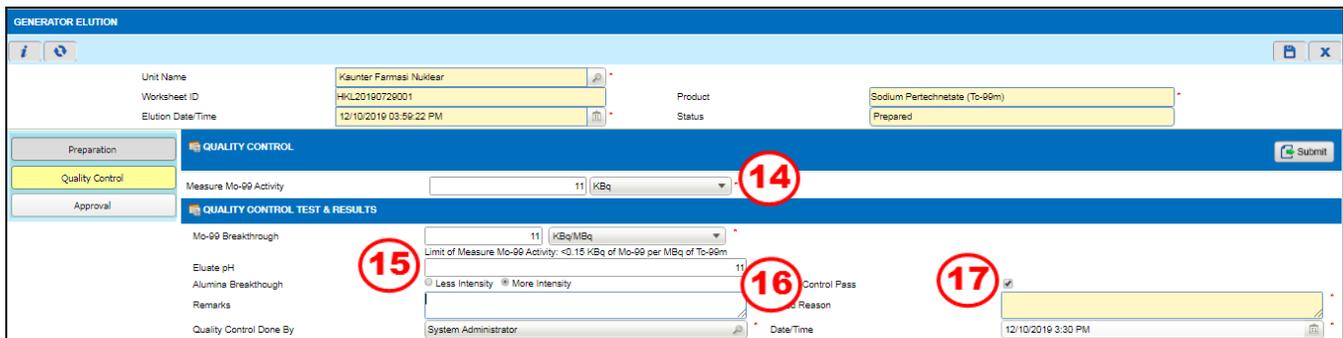


Figure 3.8.1-11 Generator Elution

#### STEP 14

Enter **Measure Mo-99 Activity kBq** field

#### Note

- The **Mo-99 Breakthrough kBq per MBq** field will be automatically displayed by system after enter **Measure Mo-99 Activity kBq**.
- **Date/Time** field will be automatically displayed by system.

#### STEP 15

Enter **Eluate pH** field

#### STEP 16


Select **Alumina Breakthrough** either **Less Intensity** or **More Intensity** radio button

#### STEP 17

Select **Quality Control Pass** either **Yes** or **No**

#### Note

If the selection of **Quality Control Pass, No; Failed Reason** field will be activated and user have to enter the reason of failure as per Figure 3.8.1-11.

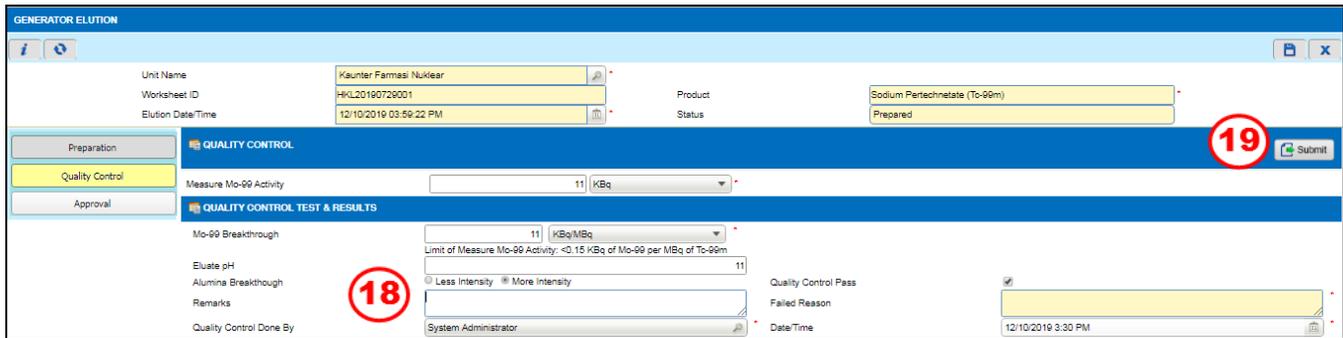

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 12/10/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: Prepared

**19** Submit

**18**

Mo-99 Breakthrough: 11 KBq/MBq  
Limit of Measure Mo-99 Activity: <1.15 KBq of Mo-99 per MBq of Tc-99m  
Eluate pH: 11  
Alumina Breakthrough:  Less Intensity  More Intensity  
Quality Control Pass:   
Remarks:   
Quality Control Done By: System Administrator  
Failed Reason:   
Date/Time: 12/10/2019 3:30 PM

**Figure 3.8.1-12 Generator Elution**

**Note**

- If the selection of **Quality Control Pass, Yes; Failed Reason** field will not activate as per Figure 3.8.1-12.


**STEP 18**

Enter **Remarks** if applicable

**STEP 19**



Click on the  button to submit quality control results

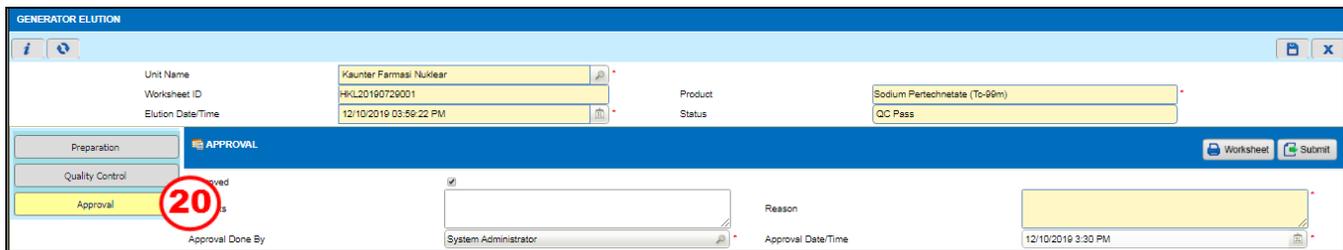

**GENERATOR ELUTION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20190729001  
Elution Date/Time: 12/10/2019 03:59:22 PM  
Product: Sodium Pertechnetate (Tc-99m)  
Status: QC Pass

**20** Submit


Approval Done By: System Administrator  
Reason:   
Approval Date/Time: 12/10/2019 3:30 PM

**Figure 3.8.1-13 Generator Elution**

**Note**

- The **Quality Control Done By** field will be automatically displayed after click on the  button.
- The  button will disable after done submitting the Quality Control record.
- The  button will be activated after the submission.

**STEP 20**

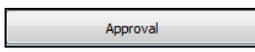Click on the  button and system will be displayed screen as per Figure 3.8.1-13

**Note**

Status will be changed to QC Pass at Approval screen.

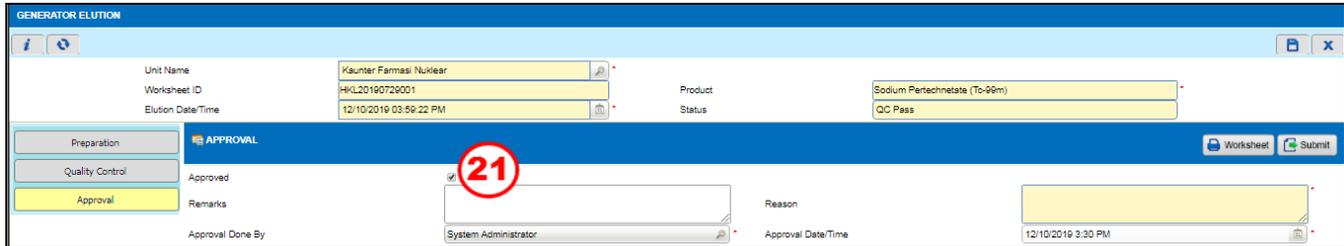
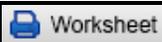


Figure 3.8.1-14 Generator Elution

**STEP 21**

Select **Approved** either **Yes** or **No**

**Note**

- If the selection of **Approved, No; Reason** field will be activated and user have to enter the reason of not approve as per Figure 3.8.1-14.
- If the selection of **Approved, Yes; Reason** field will not activate as per Figure 3.8.1-15.
- Enter **Remarks** if applicable.
- The **Date/Time** field will automatically default by system.
- The  **Worksheet** button is for user to reprint the worksheet.

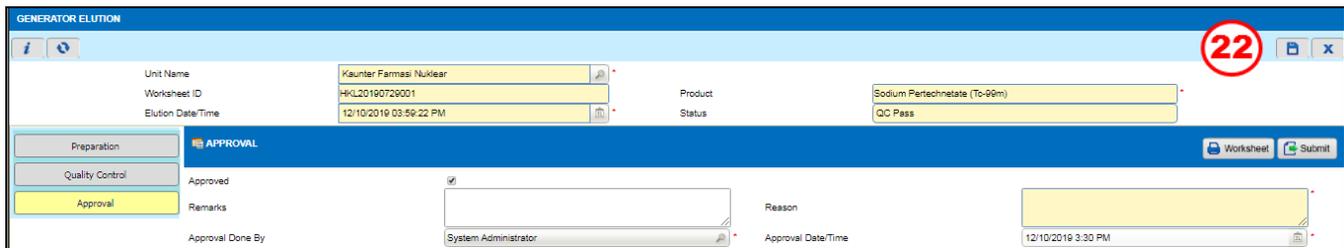


Figure 3.8.1-15 Generator Elution

**STEP 22**

Click on the  **Submit** button to submit the approval details

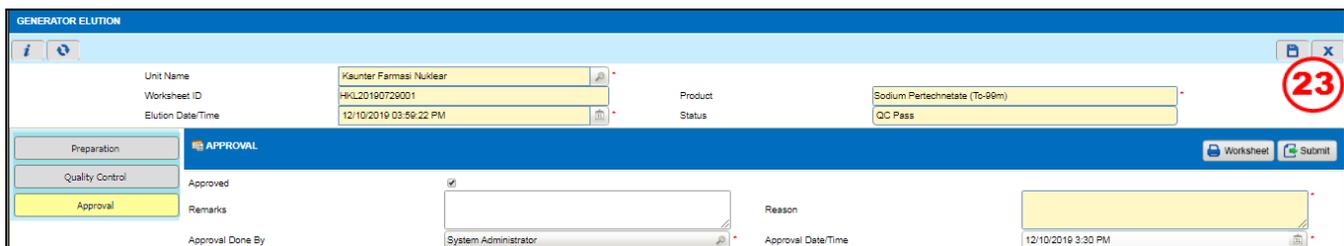


Figure 3.8.1-16 Generator Elution




**Note**

- **Status** will change to *Approved* after the submission had done.
- The **Approval Done By** field will be automatically displayed after the submission had done.

**STEP 23**

Click on the  button to close generator elution screen

### 3.8.2 KIT Preparation

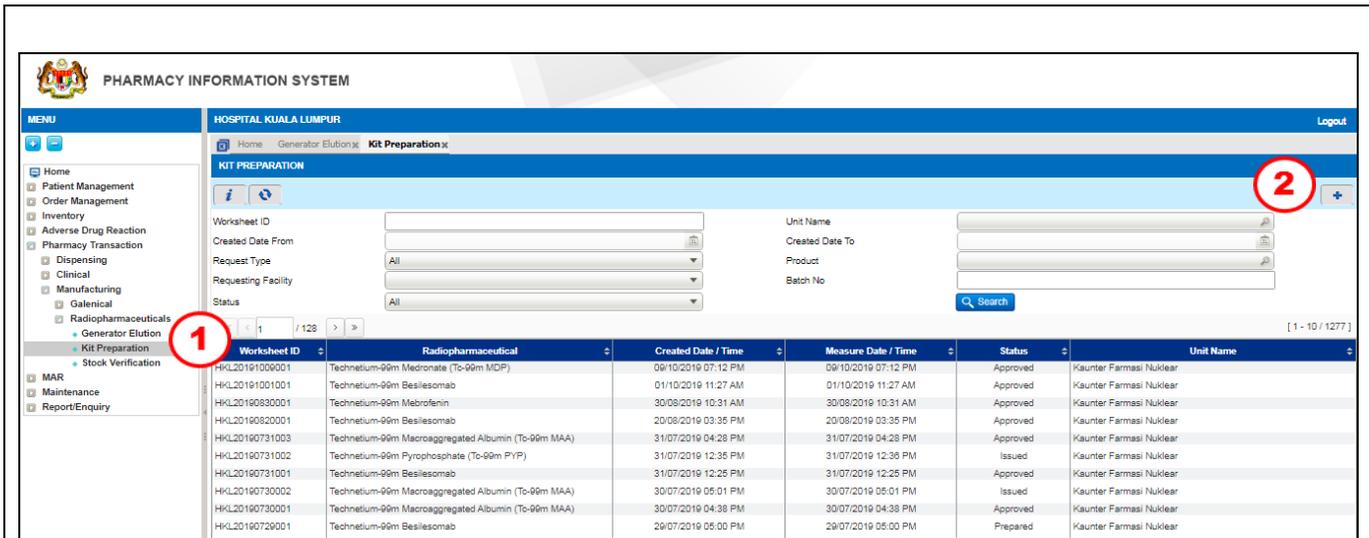

**PHARMACY INFORMATION SYSTEM**

HOSPITAL KUALA LUMPUR

Logout

Home Generator Elution Kit Preparation x

**KIT PREPARATION**

Worksheet ID:  Unit Name:

Created Date From:  Created Date To:

Request Type:  Product:

Requesting Facility:  Batch No:

Status:

[ 1 - 10 / 1277 ]

Worksheet ID	Radiopharmaceutical	Created Date / Time	Measure Date / Time	Status	Unit Name
HKL20191000001	Techneium-99m Medronate (Tc-99m MDP)	09/10/2019 07:12 PM	09/10/2019 07:12 PM	Approved	Kaunter Farmasi Nuklear
HKL20191001001	Techneium-99m Besilesomab	01/10/2019 11:27 AM	01/10/2019 11:27 AM	Approved	Kaunter Farmasi Nuklear
HKL20190830001	Techneium-99m Mebrofenin	30/08/2019 10:31 AM	30/08/2019 10:31 AM	Approved	Kaunter Farmasi Nuklear
HKL20190830001	Techneium-99m Besilesomab	20/08/2019 03:35 PM	20/08/2019 03:35 PM	Approved	Kaunter Farmasi Nuklear
HKL20190731003	Techneium-99m Macroaggregated Albumin (Tc-99m MAA)	31/07/2019 04:28 PM	31/07/2019 04:28 PM	Approved	Kaunter Farmasi Nuklear
HKL20190731002	Techneium-99m Pyrophosphate (Tc-99m PYP)	31/07/2019 12:35 PM	31/07/2019 12:35 PM	Issued	Kaunter Farmasi Nuklear
HKL20190731001	Techneium-99m Besilesomab	31/07/2019 12:25 PM	31/07/2019 12:25 PM	Approved	Kaunter Farmasi Nuklear
HKL20190730002	Techneium-99m Macroaggregated Albumin (Tc-99m MAA)	30/07/2019 05:01 PM	30/07/2019 05:01 PM	Issued	Kaunter Farmasi Nuklear
HKL20190730001	Techneium-99m Macroaggregated Albumin (Tc-99m MAA)	30/07/2019 04:38 PM	30/07/2019 04:38 PM	Approved	Kaunter Farmasi Nuklear
HKL20190720001	Techneium-99m Besilesomab	29/07/2019 05:00 PM	29/07/2019 05:00 PM	Prepared	Kaunter Farmasi Nuklear

Figure 3.8.2-1 Kit Preparation Listing Page

#### STEP 1


Click on 'Pharmacy Transaction' menu follow by 'Manufacturing' and click on 'Radiopharmaceuticals Kit Preparation' sub-menu

#### Note

- Various search criteria is provided as below:
  - Worksheet ID
  - Product
  - Ordered Date From
  - Ordered Date To
  - Batch No
  - Status
  - Request Type/Requesting Facility
- Click on the  button and system will be displayed related record.

#### STEP 2

Click on  button to create new bulk preparation record and system will be automatically displayed screen as per Figure 3.8.2-1

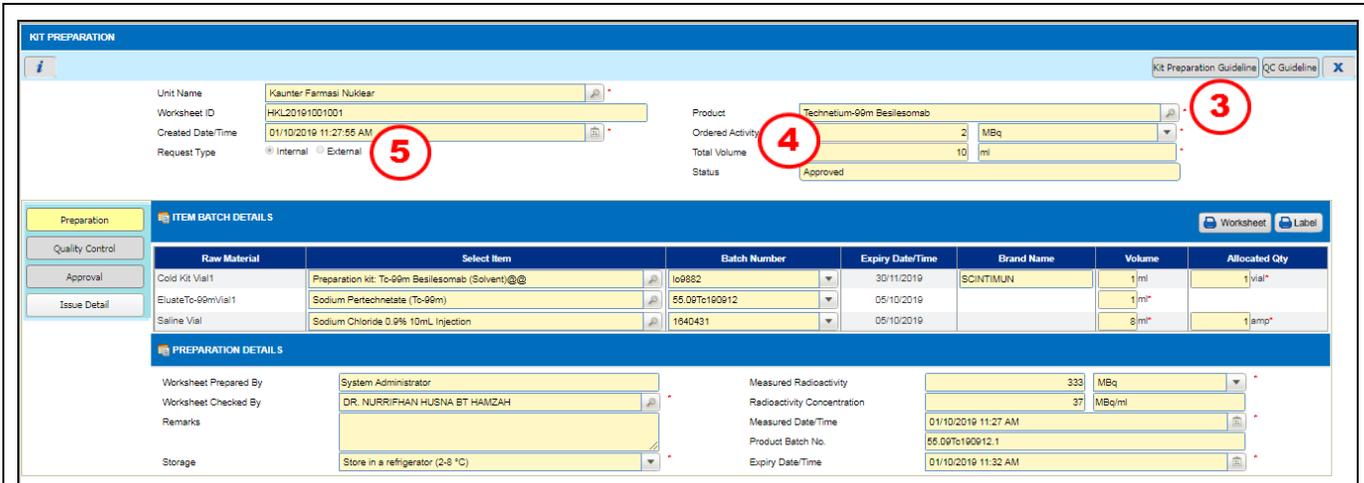

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20191001001  
Created Date/Time: 01/10/2019 11:27:55 AM  
Request Type:  Internal  External

Product: Technetium-99m Besilesomab  
Ordered Activity: 2 MBq  
Total Volume: 10 ml  
Status: Approved

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1	Preparation kit, Tc-99m Besilesomab (Solvent)@@	1c9882	30/11/2019	SCINTIMUN	1 ml	1 vial*
ElvateTc-99mVial1	Sodium Pertechnetate (Tc-99m)	55.097c190912	05/10/2019		1 ml*	
Saline Vial	Sodium Chloride 0.9% 10mL Injection	1640431	05/10/2019		8 ml*	1 amp*

**PREPARATION DETAILS**

Worksheet Prepared By: System Administrator  
Worksheet Checked By: DR. NURRIFHAN HUSNA BT HAMZAH  
Remarks:  
Storage: Store in a refrigerator (2-8 °C)

Measured Radioactivity: 333 MBq  
Radioactivity Concentration: 37 MBq/ml  
Measured Date/Time: 01/10/2019 11:27 AM  
Product Batch No.: 55.097c190912.1  
Expiry Date/Time: 01/10/2019 11:32 AM

Figure 3.8.2-2 Kit Preparation

### STEP 3

Click on the  button and select **Product** from the drop down box

### Note

- The **Unit Name** field will automatically default to current date by system.
- The **Created Date** field will automatically default to current date by system.
- The **Request Type** will automatically default to **Internal** radiobutton.

### STEP 4

Enter **Ordered Activity** and **Total Volume** field

### STEP 5

Select **Request Type** either **Internal** or **External** radiobutton

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: [Blank]  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: New

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1 Add	Preparation kit: Tc-99m Mebrofenin@@	HJU8754654587898776	31/10/2019	BRIDATEC	1 ml	1 vial*
Eluate Tc-99m Vial1 Add	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	
Saline Vial Add	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1 EA*

**PREPARATION DETAILS**

Worksheet Prepared By: [Blank]  
Worksheet Checked By: DR. NURRIHAN HUSNA BT HAMZAH  
Remarks: [Blank]  
Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq  
Radioactivity Concentration: 5.5 MBq/ml  
Measured Date/Time: 12/10/2019 4:02 PM  
Product Batch No.: [Blank]  
Expiry Date/Time: [Blank]

Figure 3.8.2-4 Kit Preparation

**STEP 6**

Select the **Item** from drop down box for every **Raw Material** details

**Note**


- **Batch Number, Expiry Date and Brand Name** will be automatically displayed after done the selection for **Item** field.
- **Enter Volume and Allocated Qty** field.
- **Cold Kit Vial1, Eluate Tc-99m Vial1 and Saline Vial** are mandatory field.

**STEP 7**

Enter **Measured Radioactivity** field and **Radioactivity Concentration** will be automatically calculated and displayed by system

**STEP 8**

Enter **Remarks** if applicable

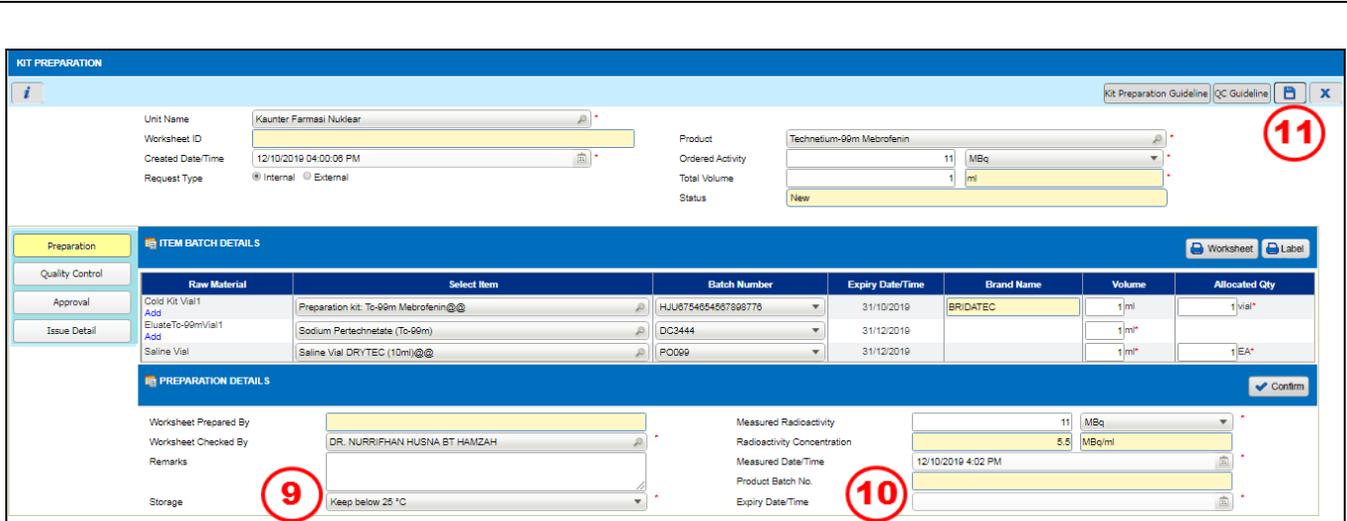

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: [Blank]  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: New

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1 Add	Preparation kit: Tc-99m Mebrofenin@@	HJU8754654587898776	31/10/2019	BRIDATEC	1 ml	1 vial*
BuadsTo-99mVial1 Add	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	
Saline Vial	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1 EA*

Worksheet Prepared By: [Blank]  
Worksheet Checked By: DR. NURRIHAN HUSNA BT HAMZAH  
Remarks: [Blank]  
Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq  
Radioactivity Concentration: 5.5 MBq/ml  
Measured Date/Time: 12/10/2019 4:02 PM  
Product Batch No.: [Blank]  
Expiry Date/Time: [Blank]

Figure 3.8.2-5 Kit Preparation

**Note**

- **Measured Date/Time** and **Worksheet Prepared By** will be automatically displayed by system.
- **Measured Date/Time** is mandatory field.

**STEP 9**

Select **Storage** from drop down box:

- Store in refrigerator (2-8°C )
- Keep below 25°C
- Keep in room temperature

**Note**


- **Storage** is mandatory field.

**STEP 10**

Select **Expiry Date** and this field is mandatory

**STEP 11**

Click on the  button to save record

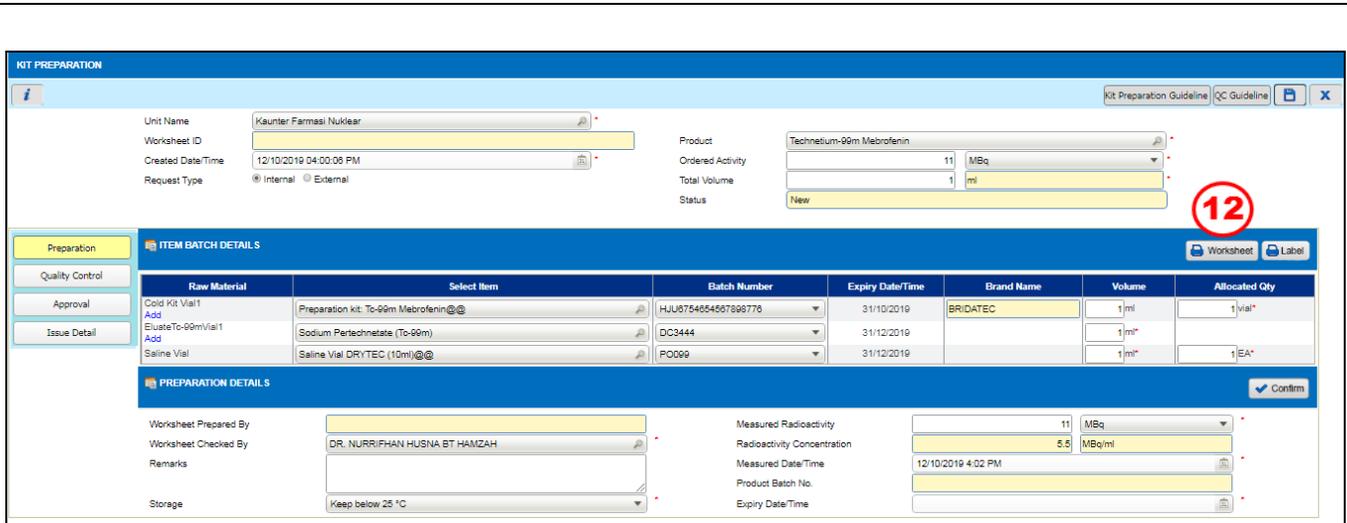

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: [Field]  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: New

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1 Add	Preparation kit: Tc-99m Mebrofenin@@	HJU8754654587898776	31/10/2019	BRIDATEC	1 ml	1 vial*
BudesTo-99mVial1 Add	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	
Saline Vial	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1 EA*

**PREPARATION DETAILS**

Worksheet Prepared By: [Field]  
Worksheet Checked By: DR. NURRIHAN HUSNA BT HAMZAH  
Remarks: [Field]  
Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq  
Radioactivity Concentration: 5.5 MBq/ml  
Measured Date/Time: 12/10/2019 4:02 PM  
Product Batch No.: [Field]  
Expiry Date/Time: [Field]

Buttons: Worksheet, Label


Figure 3.8.2-6 Kit Preparation

**Note**

**Status** will change to **New** after saving record.

- The  button will be activated after saving.
- The  button will be activated after generate worksheet.

**STEP 12**

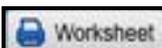Click on the  button to generate worksheet and alert message will be displayed as per Figure 3.7.2-6

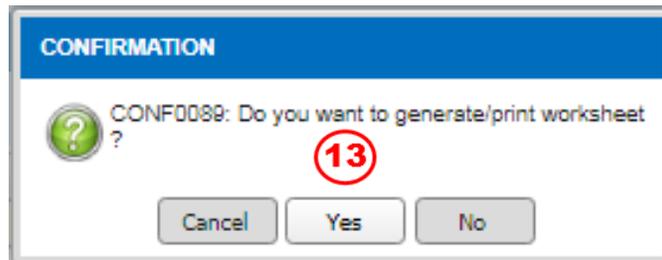


Figure 3.8.2-7 Bulk Preparation

**STEP 13**

Click on the  button to generate/print worksheet and worksheet will be displayed as per Figure 3.8.2-7

RADIOPHARMACEUTICALS WORKSHEET					
MINISTRY OF HEALTH HOSPITAL KUALA LUMPUR					
Kit Preparation records for Technetium-99m Mebrofenin					
Batch No:					
Worksheet ID: HKL20191012001					
Preparation of Technetium-99m Mebrofenin			Cold Kit Brand/Manufacturer Name:		Worksheet Prepared By : System Administrator  Signature Name Date/Time  Worksheet Checked By : DR. NURRIFHAN HUSNA BT HAMZAH  Signature Name Date/Time  Preparation Done By :  Signature Name Date/Time
Date	12/10/2019	Expected Activity	11 MBq		
Raw Material	Batch Number/Expiry		Volume		
Cold Kit Vial	HJU6754654567898776 -		Total : 1 ml 1 ml		
Eluate Tc-99m Vial	DC3444 - 31/12/2019		1 ml		
Saline Vial	PO099 - 31/12/2019		1 ml		
Measured Activity					
Measured Date/Time					
Radioactivity Concentration					
Quality Control of Technetium-99m Mebrofenin			Label		Approval of Technetium-99m Mebrofenin  Checklist for approval : Preparation done accordingly <input type="checkbox"/> QC done accordingly and pass <input type="checkbox"/> Label done accordingly <input type="checkbox"/> Remarks :  Approved By :  Signature Name Date/Time
QC Sampling Time	Measure		Actual	%	
First Strip	BG Reading	N/A			
	1st Portion	N/A	N/A	N/A	
	2nd Portion	N/A	N/A	N/A	
Second Strip	BG Reading	N/A			
	3rd Portion	N/A	N/A	N/A	
	4th Portion	N/A	N/A	N/A	
Radiochemical Purity					
Remarks :			Quality Control Done By		

Figure 3.8.2-8 Sample Worksheet

**Note**

Sample of worksheet as per Figure 3.8.2-8.

KIT PREPARATION

Unit Name: Kaunter Farmasi Nuklear

Worksheet ID: [Blank]

Created Date/Time: 12/10/2019 04:00:09 PM

Request Type:  Internal  External

Product: Technetium-99m Mebrofenin

Ordered Activity: 11 MBq

Total Volume: 1 ml

Status: New

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1	Preparation kit: Tc-99m Mebrofenin@@	HJU6754654567898776	31/10/2019	BRIDATEC	1 ml	1  vial*
EluateTc-99m/Vial1	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	
Saline Vial	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1  EA*

**PREPARATION DETAILS**

Worksheet Prepared By: [Blank]

Worksheet Checked By: DR. NURRIFHAN HUSNA BT HAMZAH

Remarks: [Blank]

Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq

Radioactivity Concentration: 5.5 MBq/ml

Measured Date/Time: 12/10/2019 4:02 PM

Product Batch No: [Blank]

Expiry Date/Time: [Blank]

14


Figure 3.8.2-9 Bulk Preparation

**STEP 14**

Click on the  button to edit label and **Drug Label** screen will be displayed as per Figure 3.8.2-9

**Note**

**Status** will be changed to **In Progress** after generated worksheet as per Figure 3.8.2-9.

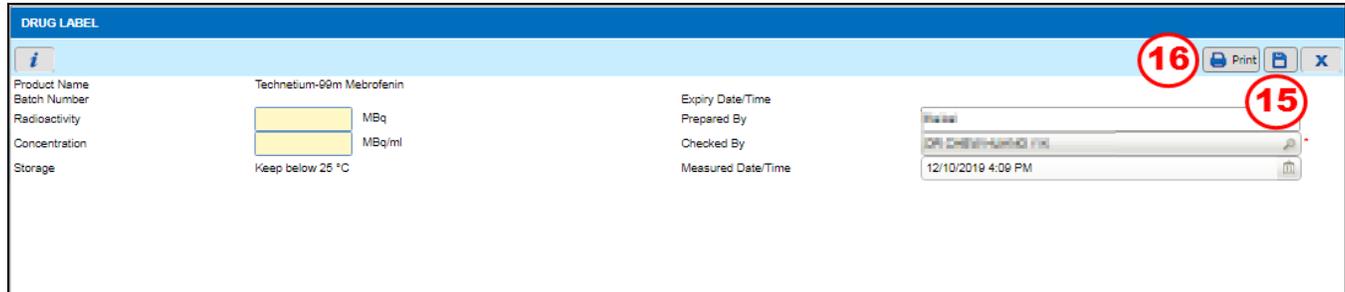

DRUG LABEL			
Product Name	Technetium-99m Mebrofenin		
Batch Number		Expiry Date/Time	
Radioactivity	<input type="text" value=""/>	Prepared By	
Concentration	<input type="text" value=""/> MBq	Checked By	<input type="text" value="DR CHEP-HUMDI PR"/>
Storage	<input type="text" value=""/> MBq/ml	Measured Date/Time	<input type="text" value="12/10/2019 4:09 PM"/>
	Keep below 25 °C		

**Figure 3.8.2-10 Drug Label**

**STEP 15**

Click on the  button to save edited label

**STEP 16**

Click on the  button to print label as per Figure 3.8.2-10

 Caution : Radioactive Material	<b>HOSPITAL KUALA LUMPUR</b>	
	<b>Product Name</b> :	Technetium-99m Mebrofenin
	<b>Batch No.</b> :	
<b>Radioactivity</b> :	0 MBq	<b>Volume</b> : 1 ml
<b>Concentration</b> :	0 MBq/ml	
<b>Measured Date</b> :	12/10/2019	<b>Time</b> : 04:09 PM
<b>Prepared By</b> :	Maikal	
<b>Checked By</b> :	DR CHEW HUANG YIK	
<b>Expiry Date</b> :		<b>Time</b> :
<b>Storage</b> :	Keep below 25 °C	

 Caution : Radioactive Material	<b>HOSPITAL KUALA LUMPUR</b>	
	<b>Product Name</b> :	Technetium-99m Mebrofenin
	<b>Batch No.</b> :	
<b>Radioactivity</b> :	0 MBq	<b>Volume</b> : 1 ml
<b>Concentration</b> :	0 MBq/ml	
<b>Measured Date</b> :	12/10/2019	<b>Time</b> : 04:09 PM
<b>Prepared By</b> :	Maikal	
<b>Checked By</b> :	DR CHEW HUANG YIK	
<b>Expiry Date</b> :		<b>Time</b> :
<b>Storage</b> :	Keep below 25 °C	

Figure 3.8.2-11 Drug Label

**Note**

Sample of Drug Label as per Figure 3.8.2-11.

**DRUG LABEL**

Product Name: Technetium-99m Mebrofenin

Batch Number:

Radioactivity:  MBq

Concentration:  MBq/ml

Storage: Keep below 25 °C

Expiry Date/Time:

Prepared By:

Checked By:

Measured Date/Time:

17

Figure 3.8.2-12 Drug Label

**STEP 17**

Click on the  button to close **Drug Label** screen

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: [Blank]  
Created Date/Time: 12/10/2019 04:00:08 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: New

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1	Preparation kit: Tc-99m Mebrofenin@@	HJU8754654567898776	31/10/2019	BRIDATEC	1 ml	1 vial*
Eluate Tc-99m Vial1	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	1
Saline Vial	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1

**PREPARATION DETAILS**

Worksheet Prepared By: [Blank]  
Worksheet Checked By: DR. NURRIFAHAN HUSNA BT HAMZAH  
Remarks: [Blank]  
Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq  
Radioactivity Concentration: 5.5 MBq/ml  
Measured Date/Time: 12/10/2019 4:02 PM  
Product Batch No.: [Blank]  
Expiry Date/Time: [Blank]

**18**

Confirm

Figure 3.8.2-13 Drug Label

**STEP 18**

Click on the  button to confirm information preparation details

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20191012001  
Created Date/Time: 12/10/2019 04:00:08 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: Prepared

**ITEM BATCH DETAILS**

Raw Material	Select Item	Batch Number	Expiry Date/Time	Brand Name	Volume	Allocated Qty
Cold Kit Vial1	Preparation kit: Tc-99m Mebrofenin@@	HJU8754654567898776	31/10/2019	BRIDATEC	1 ml	1 vial*
Eluate Tc-99m Vial1	Sodium Pertechnetate (Tc-99m)	DC3444	31/12/2019		1 ml*	1
Saline Vial	Saline Vial DRYTEC (10ml)@@	PO099	31/12/2019		1 ml*	1 EA*

**PREPARATION DETAILS**


Worksheet Prepared By: System Administrator  
Worksheet Checked By: DR. NURRIFAHAN HUSNA BT HAMZAH  
Remarks: [Blank]  
Storage: Keep below 25 °C

Measured Radioactivity: 11 MBq  
Radioactivity Concentration: 5.5 MBq/ml  
Measured Date/Time: 12/10/2019 4:02 PM  
Product Batch No.: DC3444.1  
Expiry Date/Time: 12/10/2019 4:16 PM

Figure 3.8.2-14 Kit Preparation

**Note**

- **Status** will change to **Prepared** after preparation details confirmation.
- The  button will be activated after preparation details confirmation
- The  button will disable after preparation details confirmation

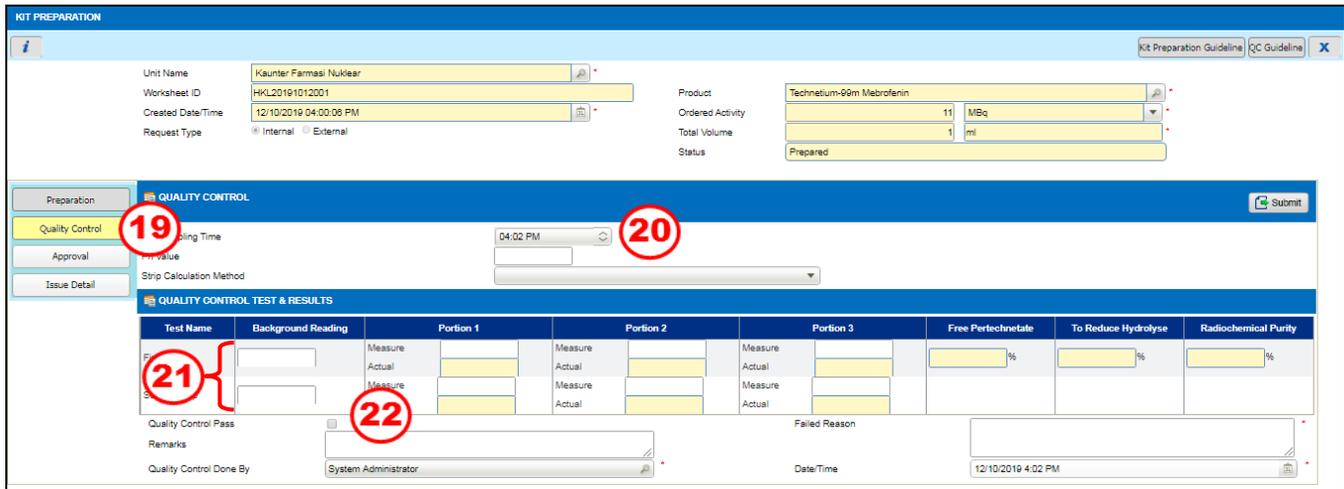

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20191012001  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: Prepared

**QUALITY CONTROL**

Quality Control (19) Sampling Time: 04:02 PM (20)

**QUALITY CONTROL TEST & RESULTS**

Test Name	Background Reading	Portion 1	Portion 2	Portion 3	Free Per technetate	To Reduce Hydrolyse	Radiochemical Purity
F	Measure	Measure	Measure	Measure			
S	Actual	Actual	Actual	Actual	%	%	%
S	Measure	Measure	Measure	Measure			
S	Actual	Actual	Actual	Actual			

Quality Control Pass (22)  
Remarks  
Quality Control Done By: System Administrator  
Failed Reason  
Date/Time: 12/10/2019 4:02 PM

Figure 3.8.2-15 Bulk Preparation

**STEP 19**

Click on the  button and system will be displayed screen as per Figure 3.8.2-15

**STEP 20**

Enter **Qc-Sampling Time** field. Default value is based on the current system time

**STEP 21**

Enter **First Strip, Second Strip, Radiochemical Purify** and **Ph Value** field as per Figure 3.8.2-15

**STEP 22**

Select **Quality Control Pass** either **approved** or **rejected**

**Note**

- If the selection of **Quality Control Pass, No; Failed Reason** field will be activated and user have to enter the reason of failure as per Figure 3.7.2-15.
- If the selection of **Quality Control Pass, Yes; Failed Reason** field will not activate.

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20191012001  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: Prepared

Preparation | **QUALITY CONTROL** | Submit

QC Sampling Time: 04:02 PM  
Ph Value: 11  
Strip Calculation Method: Biodex

**QUALITY CONTROL TEST & RESULTS**

Test Name	Background Reading	Portion 1	Portion 2	Portion 3	Free Pertechnetate	To Reduce Hydrolyse	Radiochemical Purity
First Strip	Measure	37	Measure		Measure		
	Actual	-83	Actual		Actual		
Second Strip	Measure	11	Measure		Measure		
	Actual	-12	Actual		Actual		

Quality Control Pass:   
Remarks:  **23**  
Quality Control Done By: System Administrator  
Date/Time: 12/10/2019 4:02 PM

Figure 3.8.2-16 Bulk Preparation

**STEP 23**

Enter **Remarks** if applicable

**Note**

**Quality Control Done By and Date/Time** will be automatically displayed and default by system.

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
Worksheet ID: HKL20191012001  
Created Date/Time: 12/10/2019 04:00:06 PM  
Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
Ordered Activity: 11 MBq  
Total Volume: 1 ml  
Status: Prepared

Preparation | **QUALITY CONTROL** | Submit

QC Sampling Time: 04:02 PM  
Ph Value: 11  
Strip Calculation Method: Biodex

**QUALITY CONTROL TEST & RESULTS**

Test Name	Background Reading	Portion 1	Portion 2	Portion 3	Free Pertechnetate	To Reduce Hydrolyse	Radiochemical Purity
First Strip	Measure	37	Measure		Measure		
	Actual	-83	Actual		Actual		
Second Strip	Measure	11	Measure		Measure		
	Actual	-12	Actual		Actual		

Quality Control Pass:   
Remarks:   
Quality Control Done By: System Administrator  
Date/Time: 12/10/2019 4:02 PM

Figure 3.8.2-17 Bulk Preparation

**Note**

- User is allowed to click on the **Strip Guideline** hyperlink and system will be displayed as per Figure 3.8.2-17.

Figure 3.8.2-17 displays multiple pages of Quality Control Test - Radiopharmaceutical Reconciliation tables. Each page includes a header with the title 'QUALITY CONTROL TEST - RADIOPHARMACEUTICAL RECONCILIATION' and a table with columns for 'Isotope', 'Radionuclide', 'Activity', 'Time', and 'Calculation'. The tables specify acceptance criteria and calculation formulas for various parameters such as % of total activity, % of radiochemical purity, and % of radiochemical impurity.

**Note**

Sample of Strip Guideline as per Figure 3.8.2-18.

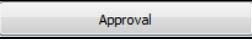
Figure 3.8.2-19 is a screenshot of the 'KIT PREPARATION' software interface. It shows a form for entering unit information (Unit Name: Kaunter Farmasi Nuklear, Worksheet ID: HKL20191012001, Created Date/Time: 12/10/2019 04:00:06 PM, Request Type: Internal) and product details (Product: Technetium-99m Mebrofenin, Ordered Activity: 11 MBq, Total Volume: 1 ml, Status: Prepared). A 'Submit' button is highlighted with a red circle and the number 24. Below the form is a 'QUALITY CONTROL TEST & RESULTS' table with columns for Test Name, Background Reading, Portion 1, Portion 2, Portion 3, Free Pertechnetate, To Reduce Hydrolyse, and Radiochemical Purity. The table shows data for 'First Strip' and 'Second Strip' with 'Measure' and 'Actual' values. A 'Quality Control Pass' checkbox is checked, and the 'Quality Control Done By' field is filled with 'System Administrator' and the date/time '12/10/2019 4:02 PM'.


Figure 3.8.2-19 Bulk Preparation

**STEP 24**

Click on the button to submit quality control results

**Note**

- The  button will disable after done the submission.
- The  button will be activated after the submission.

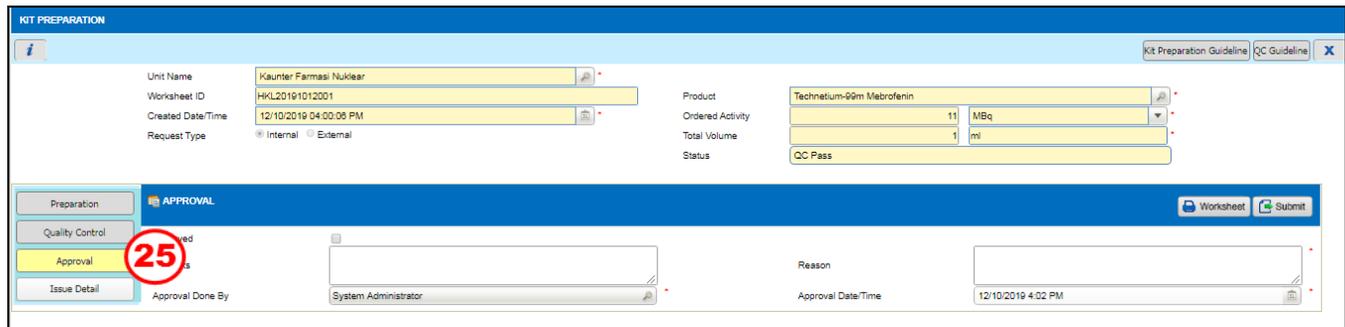

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
 Worksheet ID: HKL20191012001  
 Created Date/Time: 12/10/2019 04:00:08 PM  
 Request Type:  Internal  External


Product: Technetium-99m Mebrofenin  
 Ordered Activity: 11 MBq  
 Total Volume: 1 ml  
 Status: QC Pass

Preparation | **APPROVAL** | Quality Control | Issue Detail

Approved:   
 Reason:   
 Approval Done By: System Administrator  
 Approval Date/Time: 12/10/2019 4:02 PM


**Figure 3.8.2-20 Bulk Preparation**

**STEP 25**

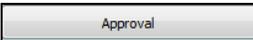Click on the  button and system will display screen as per Figure 3.8.2-20

**Note**

Status will be changed to QC Pass after done the submission.

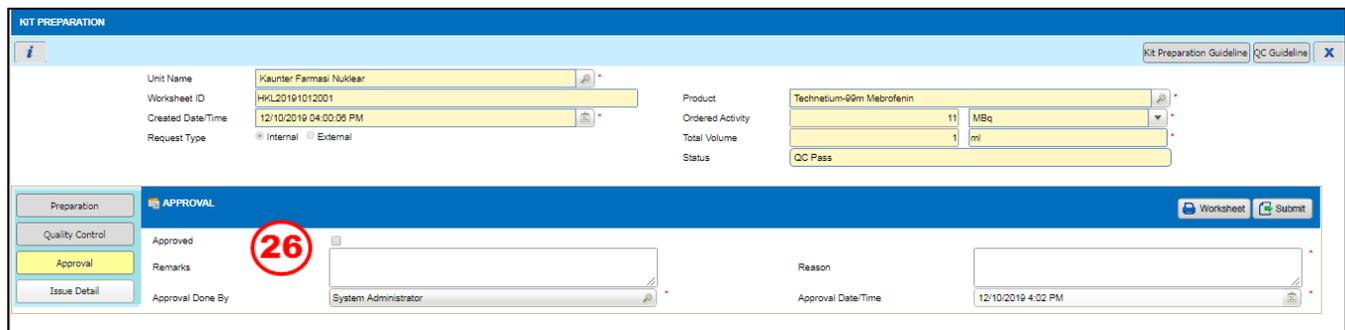

**KIT PREPARATION**

Unit Name: Kaunter Farmasi Nuklear  
 Worksheet ID: HKL20191012001  
 Created Date/Time: 12/10/2019 04:00:08 PM  
 Request Type:  Internal  External

Product: Technetium-99m Mebrofenin  
 Ordered Activity: 11 MBq  
 Total Volume: 1 ml  
 Status: QC Pass

Preparation | **APPROVAL** | Quality Control | Issue Detail

Approved:  **26**  
 Remarks:   
 Reason:   
 Approval Done By: System Administrator  
 Approval Date/Time: 12/10/2019 4:02 PM


**Figure 3.8.2-21 Bulk Preparation**


**STEP 26**

Select **Approved** either **Yes** or **No**

**Note**

- If the selection of **Approved, No**; **Reason** field will be activated and user have to enter the reason of not approve as per Figure 3.8.2-21.
- If the selection of **Approved, Yes**; **Reason** field will not activate as per Figure 3.8.2-21.
- Enter **Remarks** if applicable.
- The **Approval Done By** and **Date/Time** field will automatically default by system.

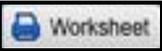- The  button is for user to reprint the worksheet.

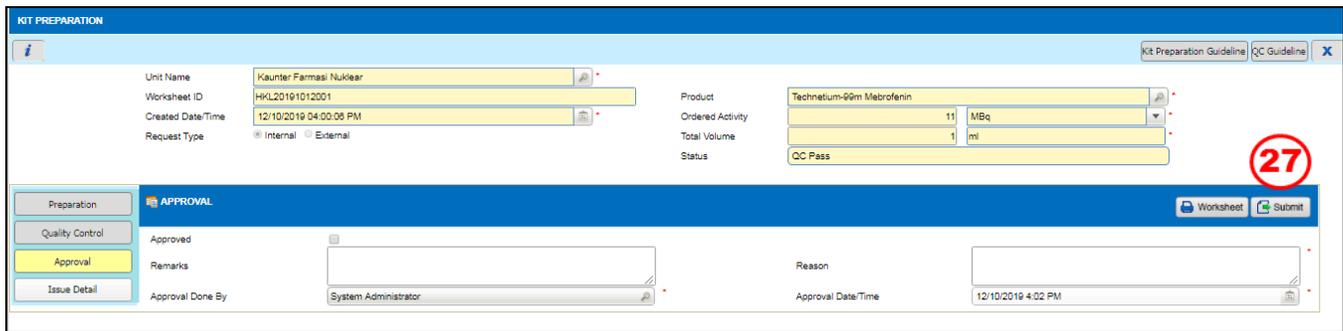


The screenshot shows the 'KIT PREPARATION' interface. The 'APPROVAL' section is active, displaying fields for 'Approved', 'Remarks', 'Reason', 'Approval Done By' (System Administrator), and 'Approval Date/Time' (12/10/2019 4:02 PM). A red circle with the number 27 highlights the 'Worksheet' button in the top right corner of the approval section.

Figure 3.8.2-22 Bulk Preparation

### STEP 27

Click on the  button to submit the approval details

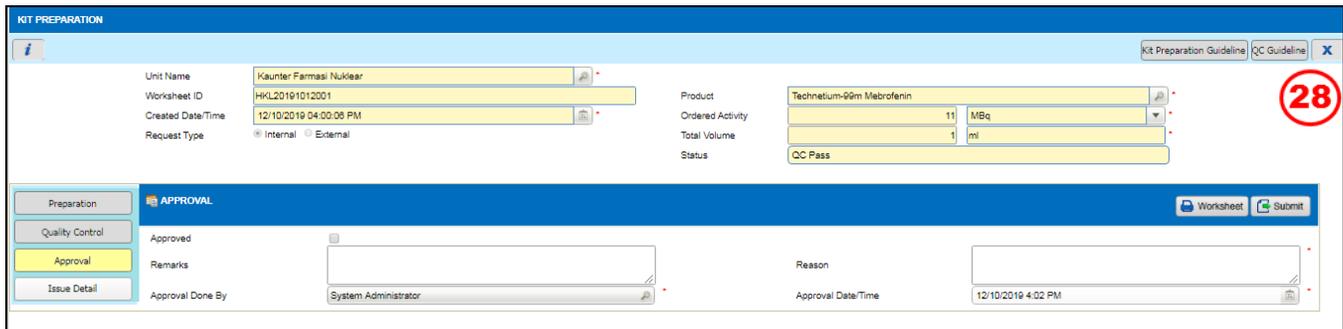

The screenshot shows the 'KIT PREPARATION' interface. The 'APPROVAL' section is active, displaying fields for 'Approved', 'Remarks', 'Reason', 'Approval Done By' (System Administrator), and 'Approval Date/Time' (12/10/2019 4:02 PM). A red circle with the number 28 highlights the 'Submit' button in the top right corner of the approval section.

Figure 3.8.2-23 Bulk Preparation

### Note


- **Status** will change to **Approved** after the submission had done.

### STEP 28

Click on the  button to close Bulk Preparation screen

### 3.8.3 Stock Verification

Stock Verification screen allow user to verify and measure radiopharmaceuticals stock received from pharmacy store.

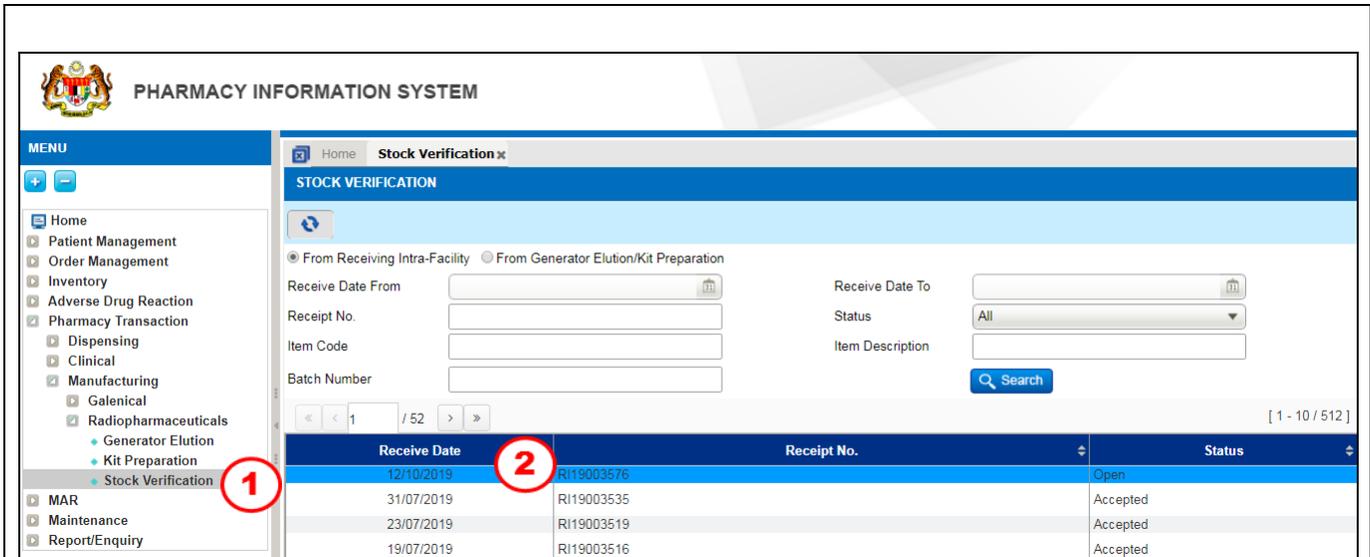

Receive Date	Receipt No.	Status
12/10/2019	RI19003576	Open
31/07/2019	RI19003535	Accepted
23/07/2019	RI19003519	Accepted
19/07/2019	RI19003516	Accepted

Figure 3.8.3-1 Stock Verification Listing Page

#### STEP 1


Click on 'Pharmacy Transaction' menu follow by 'Manufacturing' and click on 'Radiopharmaceuticals Stock Verification' sub-menu

#### Note

- Search criteria is provided for **Date** only.
- Click on the  button and system will be displayed related record.

#### STEP 2

Double click on the selected data and system will be displayed screen as per Figure 3.8.3-2

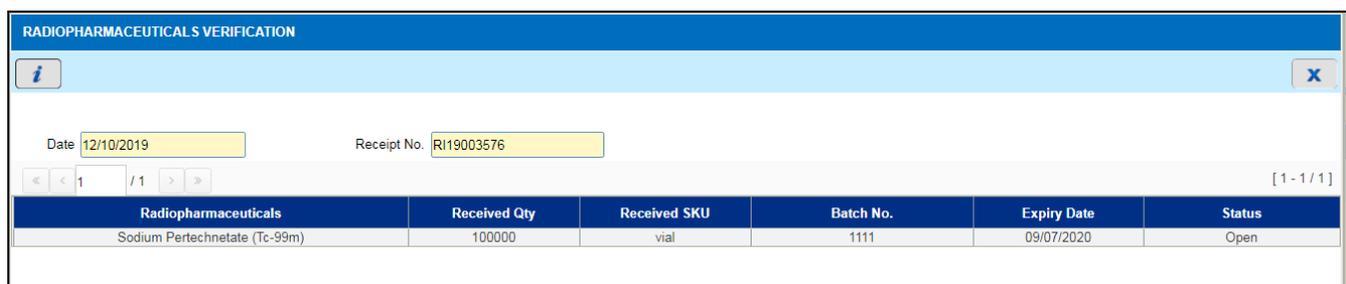

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	100000	vial	1111	09/07/2020	Open

Figure 3.8.3-2 Stock Verification

#### Note

- **Date** and **Receipt No** will be automatically displayed and default by system.
- By default, status is open.

**RADIOPHARMACEUTICALS VERIFICATION**

Date: 12/10/2019
Receipt No.: RI19003576

[ 1 - 1 / 1 ]

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	100000	vial	1111	09/07/2020	Open

**STOCK DETAILS**

Calibration Activity	<input type="text"/>	MBq *	<input type="text"/>	mCi *	Calibration Date/Time	<input type="text"/>	<input type="button" value="Calendar"/>	
Expected Activity	<input type="text"/>	MBq *	<input type="text"/>	mCi *	Expiry Date/Time	<input type="text"/>	<input type="button" value="Calendar"/>	
Measured Activity	<input type="text"/>	MBq *	<input type="text"/>	mCi *	Measured Date/Time	<input type="text"/>	<input type="button" value="Calendar"/>	
Activity Variation	<input type="text"/>	%	Volume					<input type="text"/>
						ml	<input type="button" value="Dropdown"/>	

**MEASUREMENT DETAILS**

Label Packaging Type	<input type="text"/>	<a href="#">Guideline</a>	
Survey on Surface	<input type="text"/>	<input type="text"/>	
Survey at 1 Meter Distance	<input type="text"/>	<input type="text"/>	
Wipe Test	<input type="text"/>	<input type="text"/>	
Remarks	<input style="height: 20px;" type="text"/>		

**VERIFICATION DETAILS**

Accept
 Reject

Reject Reason	<input style="height: 20px;" type="text"/>
Done By	<input type="text" value="System Administrator"/>
Updated By	<input type="text" value="System Administrator"/>

Date/Time	<input type="text" value="22/10/2019 12:26 PM"/>
Updated Date/Time	<input type="text" value="22/10/2019 12:26 PM"/>

Figure 3.8.3-3 Stock Verification

**STEP 3**

Click on the  **Accept** link to verify stock and system will be displayed screen as per Figure 3.8.3-3

**Note**

- Enter and select **Calibration Activity, Calibration Date/Time, Measured Activity, Measured Date/Time, Wipe Test, Survey at 1 Meter Distance, Expiry Date/Time, and Volume** field. All these fields are mandatory.
- Enter and select **Label Packaging Type, Survey on Surface, Survey at 1 Meter Distance, Wipe Test and Remarks**.
- **Done By, Date/Time, Updated By, and Updated Date/Time** will be automatically displayed by system.

PB\_U. MANUAL\_RADIOPHARMACEUTICALS-13<sup>th</sup> E

Page 69

**RADIOPHARMACEUTICALS VERIFICATION**

Date: 12/10/2019      Receipt No.: RI19003576

<< < 1 / 1 > >

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	100000	vial	1111	09/07/2020	Open

**STOCK DETAILS**

Calibration Activity:  MBq \*  mCi \*      Calibration Date/Time:  \* **4**

Expected Activity:  MBq \*  mCi \*      Expiry Date/Time:  \*

Measured Activity:  MBq \*  mCi \*      Measured Date/Time:  \*

Activity Variation:  %      Volume:  ml

**MEASUREMENT DETAILS**      **VERIFICATION DETAILS**

Label Packaging Type: White I [Guideline](#)

Survey on Surface:  mR/hr

Survey at 1 Meter Distance:  mR/hr

Wipe Test:  Bq per cm<sup>2</sup>

Remarks:

Accept       Reject

Reject Reason:

Done By: System Administrator      Date/Time:  \*

Updated By: System Administrator      Updated Date/Time:  \*

**5** 1 / 1

Figure 3.8.3-4 Stock Verification

**STEP 4**

Check on the  checkbox to accept or reject verification

**STEP 5**

Click on the  button to save record

**RADIOPHARMACEUTICALS VERIFICATION**

Date:  Receipt No.:

[ 1 - 1 / 1 ]

Radiopharmaceuticals	Received Qty	Received SKU	Batch No.	Expiry Date	Status
Sodium Pertechnetate (Tc-99m)	100000	vial	1111	09/07/2020	Open

**STOCK DETAILS**

Calibration Activity: <input type="text" value="1"/> MBq * <input type="text" value="0.03"/> mCi *	Calibration Date/Time: <input type="text" value="22/10/2019 12:26 PM"/>
Expected Activity: <input type="text" value="1"/> MBq * <input type="text" value="0.03"/> mCi *	Expiry Date/Time: <input type="text" value="31/10/2019 12:32 PM"/>
Measured Activity: <input type="text" value="1"/> MBq * <input type="text" value="0.03"/> mCi *	Measured Date/Time: <input type="text" value="22/10/2019 12:26 PM"/>
Activity Variation: <input type="text" value="0"/> %	Volume: <input type="text" value="1"/> ml

**MEASUREMENT DETAILS**

Label Packaging Type:  [Guideline](#)

Survey on Surface:  mR/hr

Survey at 1 Meter Distance:  mR/hr

Wipe Test:  Bq per cm<sup>2</sup>

Remarks:

**VERIFICATION DETAILS**

Accept  Reject

Reject Reason:

Done By:  Date/Time:

Updated By:  Updated Date/Time:

**Figure 3.8.3-5 Stock Verification**

**STEP 6**

Click on the  button to close Radiopharmaceuticals Verification screen

**Note**

- Status will change to Accepted as per Figure 3.8.3-6 at Stock Verification Listing Page.
- Add explanation for Rejected.



**PHARMACY INFORMATION SYSTEM**

**MENU**

- Home
- Patient Management
- Order Management
- Inventory
- Adverse Drug Reaction
- Pharmacy Transaction
- Dispensing
- Clinical
- Manufacturing
  - Galenical
  - Radiopharmaceuticals
    - Generator Elution
    - Kit Preparation
    - Stock Verification
- MAR
- Maintenance
- Report/Enquiry

**HOSPITAL KUALA LUMPUR** Logout

Home **Stock Verification**

**STOCK VERIFICATION**

From Receiving Intra-Facility  From Generator Elution/Kit Preparation

Receive Date From:

Receive Date To:

Receipt No.:

Status:

Item Code:

Item Description:

Batch Number:

[ 1 - 10 / 512 ]

Receive Date	Receipt No.	Status
12/10/2019	RI19003576	Accepted
31/07/2019	RI19003535	Accepted
23/07/2019	RI19003519	Accepted
19/07/2019	RI19003516	Accepted
19/07/2019	RI19003515	Accepted
11/07/2019	RI19003504	Accepted
10/07/2019	RI19003501	Accepted
28/05/2019	RI19003476	Accepted
28/05/2019	RI19003474	Accepted
27/05/2019	RI19003469	Accepted

**Figure 3.8.3-6 Stock Verification Listing Page**

## 4.0 Acronyms

Abbreviation	Definition
MOH	Ministry of Health
KKM	Kementerian Kesihatan Malaysia
PhIS	Pharmacy Information System
CPS	Clinical Pharmacy System
PKU	Packaging Keeping Unit
SKU	Store Keeping Unit
TDM	Therapeutic Drug Monitoring
CDR	Cytotoxic Drug Reconstitution
APPL	Approved Product Purchase List
RPL	Recommended Purchase List
EPO	Electronic Purchase Order

## 5.0 Links to Inventory Modules

No	Module	PDF Links	No	Module	PDF Links
1	<i>Finance</i>	<a href="#">Click Here</a>	15	<i>Internal Indent</i>	<a href="#">Click Here</a>
2	<i>Procurement Standard APPL</i>	<a href="#">Click Here</a>	16	<i>Issue</i>	<a href="#">Click Here</a>
3	<i>Procurement standard LP</i>	<a href="#">Click Here</a>	17	<i>Receive From Supplier</i>	<a href="#">Click Here</a>
4	<i>Procurement Standard Contract</i>	<a href="#">Click Here</a>	18	<i>Receive Inter Facility</i>	<a href="#">Click Here</a>
5	<i>Procurement Standard Quotation</i>	<a href="#">Click Here</a>	19	<i>Receive Intra Facility</i>	<a href="#">Click Here</a>
6	<i>Procurement Standard (RFQ)</i>	<a href="#">Click Here</a>	20	<i>Return to Supplier</i>	<a href="#">Click Here</a>
7	<i>Procurement Non Standard (Requisition Order)</i>	<a href="#">Click Here</a>	21	<i>Return to Supplying Unit</i>	<a href="#">Click Here</a>
8	<i>Quarantine</i>	<a href="#">Click Here</a>	22	<i>Slow Moving</i>	<a href="#">Click Here</a>
9	<i>Product Complaint</i>	<a href="#">Click Here</a>	23	<i>Stock Taking And Verification</i>	<a href="#">Click Here</a>
10	<i>Recalculate Buffer Level</i>	<a href="#">Click Here</a>	24	<i>Stock Transfer</i>	<a href="#">Click Here</a>
11	<i>Expiration And Condemn</i>	<a href="#">Click Here</a>	25	<i>Year End</i>	<a href="#">Click Here</a>
12	<i>Recall Product</i>	<a href="#">Click Here</a>	26	<i>Penalty</i>	<a href="#">Click Here</a>
13	<i>Payment</i>	<a href="#">Click Here</a>	27	<i>IWP Budget</i>	<a href="#">Click Here</a>
14	<i>External Indent</i>	<a href="#">Click Here</a>	28	<i>IWP Order Authorization</i>	<a href="#">Click Here</a>