



Pharmacy Information System (PhIS) and Clinic Pharmacy System (CPS)

Pharmacy Based

User Manual TDM Calculator

Version	: 12th Edition
Document ID	: PB_U. MANUAL_TDM CALCULATOR



CONFIDENTIAL COPYRIGHTED MATERIAL – *The information includes all concepts, comments, recommendations, and material, contained herein shall remain the property of Pharmacy Information System & Clinic Pharmacy System (PhIS & CPS) Project. No portion of this document shall be disclosed, duplicated or used in whole or in part of any purpose other than the purpose of the Pharmacy Information System & Clinic Pharmacy System (PhIS & CPS) Project execution only.*

Reference ID : PB_U. MANUAL_TDM CALCULATOR-12th E

Application reference: PhIS & CPS v2.5.1



Table of Contents

1.0 Introduction.....	1
1.1 Overview of PhIS	1
1.2 Purpose and Objectives	1
1.3 Organized Sections.....	1
2.0 Application Standard Features.....	2
2.1 PhIS Legend.....	2
3.0 Calculator for Vancomycin	4
4.0 Calculator for Gentamicin / Amikacin	7
5.0 Calculator for Valproic Acid.....	11
6.0 Calculator for Aminophylline / Theophylline	13
7.0 Calculator for Carbamazepine	15
8.0 Calculator for Digoxin.....	17
9.0 Calculator for Phenobarbitone	19
10.0 Calculator for Phenytoin.....	21
11.0 Acronyms.....	23
12.0 Links To Clinical Modules	23



1.0 Introduction

1.1 Overview of PhIS

Pharmacy Information System or better known as PhIS, is a complete and comprehensive system that integrates pharmacy related services geared towards pharmacy excellence. PhIS implementation would transform most of current manual process to electronic system would benefit facility end user in the health care sector.

There are 12 modules to assist services delivery by the health care sector which comprises of:

1. Order Management
2. Inpatient Pharmacy
3. Outpatient Pharmacy
4. Medication Counselling
5. Ward Pharmacy
6. Pharmacy Inventory
7. Manufacturing of Cytotoxic Drug Reconstitution, Parenteral Nutrition, IV Admixture & Eye Drop , Radiopharmaceuticals and Extemporaneous
8. Adverse Drug Reaction & Drug Allergic (ADR & DAC)
9. Clinical Pharmacokinetics Services (TDM)
10. Drug Information & Consumer Education (DICE)
11. Medication Therapy Adherence Clinic (MTAC)
12. Data Mining (PhARM)

1.2 Purpose and Objectives

This user manual outlines the TDM Calculator module and its key features and functionalities. The primary objective is to guide user through the process of completing PhIS application process.

User will understand the following activities in details:

- Calculator for Vancomycin
- Calculator for Gentamicin
- Calculator for Valproic Acid


1.3 Organized Sections

These are the sections within this document:

- Section 1: Introduction
- Section 2: Application Standard Features
- Section 3: Calculator for Vancomycin
- Section 4: Calculator for Gentamicin
- Section 5: Calculator for Valproic Acid
- Section 6: Acronyms
- Section 7: Links to Clinical Modules

2.0 Application Standard Features

2.1 PhIS Legend

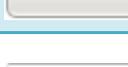Standard Legend			
	Login to PhIS		Logout from PhIS
	Change Login Password		Reset Login Screen
	Expand Menu		Collapse Menu
	Display Home Tab		Expand Module
	Collapse Module		Close All Open Tabs
	Refresh Screen		Search Record
	Add/Create New Record		Show Help
	Print	*	Mandatory Field
	Calendar Icon		Search Icon
	Close Window		Radio Button
<input type="checkbox"/>	Checkbox		Edit Record
	Cancel		Delete Record
	Save		Export and Open Report in Excel Format
	Add Item to the list		Delete Item from the list
	Request for Approval		Cancel the Request
	Approve Transaction		Reject Transaction
	Send for Approval		Dropdown Box
	Automatically Display/Retrieve Box		Empty Text Box



Therapeutic Drug Monitoring Module Legend			
Acknowledge	To Acknowledge	Drug Info	Drug Info
New	To Do New Order	Verify	To Verify
Cancel Order	Cancel the Order	+ Add New Sample	To Add New Sample
Send	To send request	CP2	Ward Pharmacy

Note

To learn more about Login Information, kindly click [Login Information](#) module for descriptive steps.

3.0 Calculator for Vancomycin

To view the calculator in the system, perform the steps below:

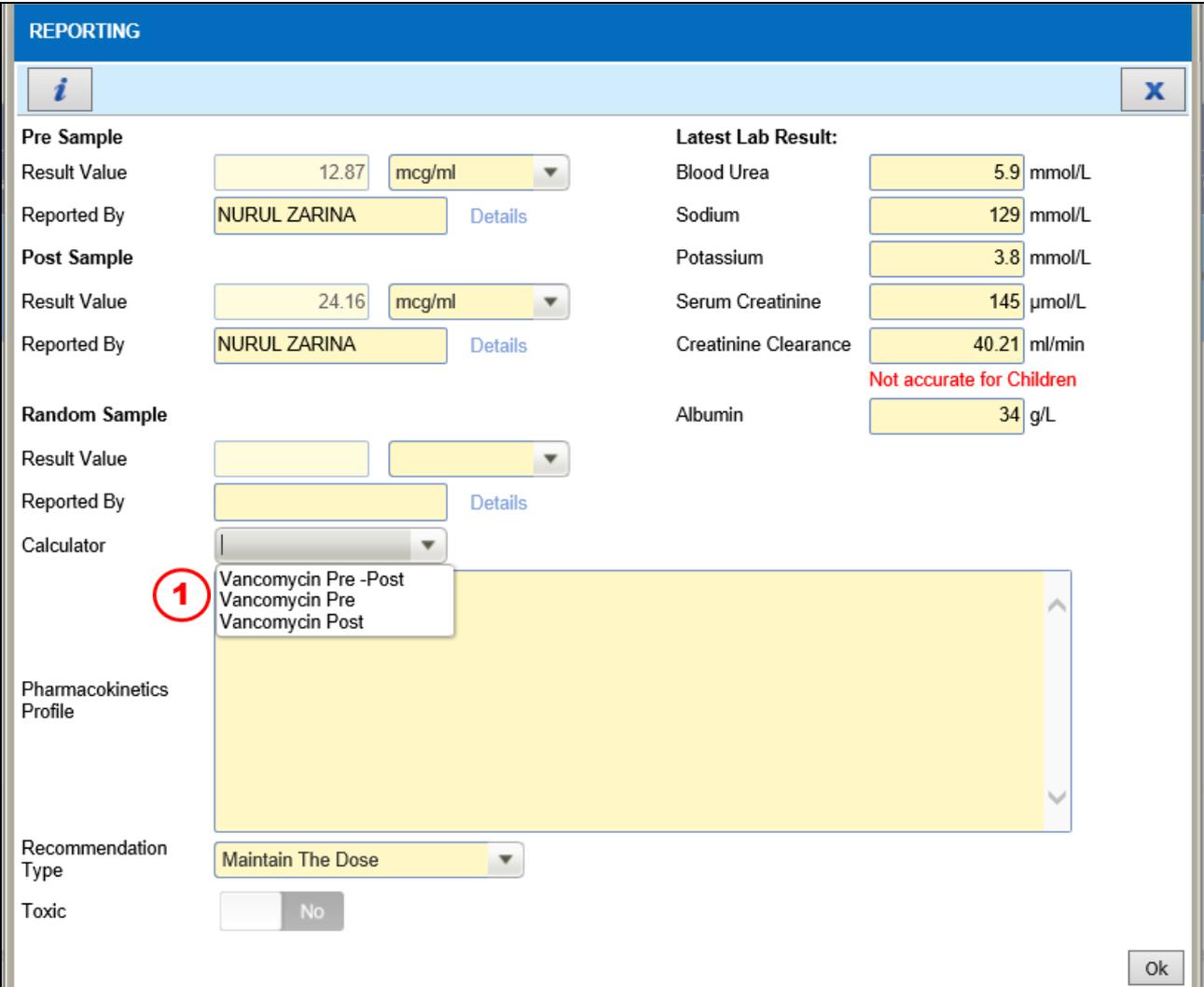


Figure 3.0-1 Pre Calculator for Vancomycin Detail

STEP 1

Select calculator from dropdown box example:

- Vancomycin in Pre
- Vancomycin in Post
- Vancomycin in Pre-Post

Note

- Choose the calculator and fill the result in the box as per Figure 3.0-1. Sample Details consists of:
 - Calculator for **Vancomycin Pre**
 - a. Dose medication in mg
 - b. Interval Time
 - c. Pre level result concentration
 - d. Serum Creatinine
 - e. Vd
 - f. CrCl - user can select option to use Adult or paediatrics calculation
 - g. Expected Cmax if assuming Expected Cmin = Pre level result

➤ Calculator for **Vancomycin Post**

- a. Dose medication in mg
- b. Interval Time
- c. Sampling Time: Duration between time after complete infusion and post sampling(t)
- d. Serum level: Post level result concentration, Serum Creatinine
- e. CrCl - user can select option to use Adult or paediatrics calculation
- f. New dose and Cmin based on desired Cmax
- g. Expected Cmin and Cmax based on new suggested dose

➤ Calculator for **Vancomycin Pre & Post**

- a. Dose medication in mg
- b. Interval Time
- c. Sampling Time: Duration between time after complete infusion and post sampling(t) and Duration between completion time of pre and post sampling (t2-t1)
- d. Serum level: Pre level result concentration, Post level result concentration and Serum Creatinine
- e. CrCl - user can select option to use Adult or paediatrics calculation
- f. Vd, Ke and t1/2
- g. Area Under Curve (AUC) Over 24hrs
- h. New dose and Cmin based on desired Cmax
- i. Expected Cmin and Cmax if Vd varies
- j. Expected Cmin and Cmax based on new suggested dose

- After save the calculator, the result will appear at the Pharmacokinetics Profile box
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 3.0-1

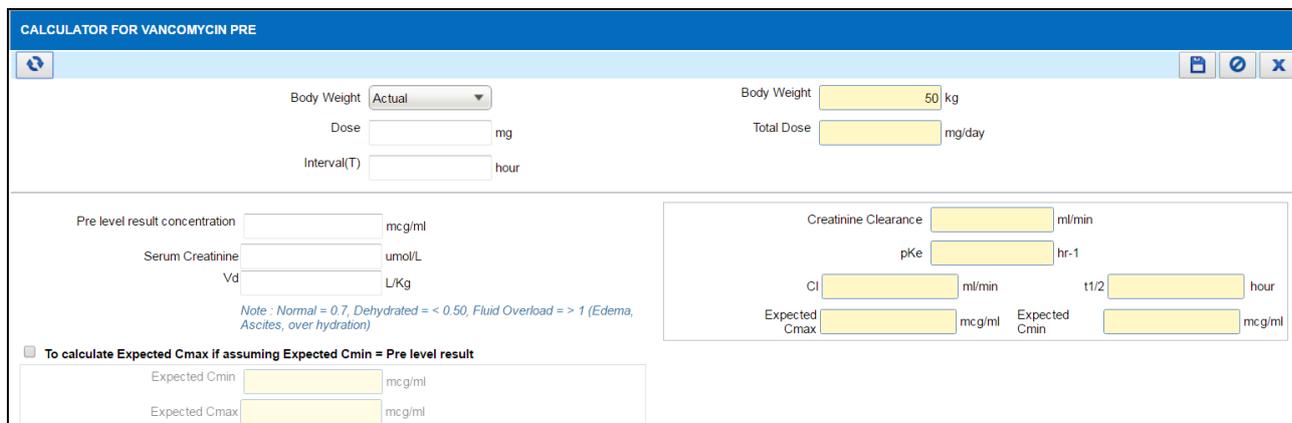


Figure 3.0-2 Pre Calculator for Vancomycin Detail

CALCULATOR FOR VANCOMYCIN POST

Body Weight: kg
 Dose: mg
 Interval(T): hour
 Total Dose: mg/day

Sampling Time
 Duration between time after complete infusion and post sampling(t): hour

Serum level
 Post level result concentration: mcg/ml
 Serum Creatinine: umol/L
 Creatinine Clearance: ml/min

To calculate expected Cmin and Cmax based on new suggested dose

pKa: hr-1
 Vd: L/kg
 Cl: ml/min
 t1/2: hour
 Cmax: mcg/ml
 Cmin: mcg/ml

New suggested dose: mg
 New Interval: hour
 Expected Cmax: mcg/ml
 Expected Cmin: mcg/ml

To calculate new dose and Cmin based on desired Cmax

Desired Cmax: mcg/ml
 Desired Interval: hour
 Expected Cmin: mcg/ml
 New Dose: mg

Figure 3.0-3 Post Calculator for Vancomycin Detail

CALCULATOR FOR VANCOMYCIN PRE & POST

Body Weight: kg
 Dose: mg
 Interval(T): hour
 Total Dose: mg/day

Sampling Time
 Duration between time after complete infusion and post sampling(t): hour
 Duration between completion time of pre and post sampling (t2-t1): hour

Serum level
 Pre level result concentration: mcg/ml
 Post level result concentration: mcg/ml
 Serum Creatinine: umol/L
 Creatinine Clearance: ml/min

To calculate Vd, Ke and t1/2

Ke: hr-1
 Vd: L/kg
 Cl: ml/min
 t1/2: hour
 Cmax: mcg/ml
 Cmin: mcg/ml

To calculate Area Under Curve (AUC) Over 24hrs

Vd: L/kg
 MIC: mg/L
 Note MIC = 2mg/L AUC = 400
 AUC:

To calculate expected Cmin and Cmax if Vd varies

Vd: L/kg
 Expected Cmax: mcg/ml
 Expected Cmin: mcg/ml

To calculate expected Cmin and Cmax based on new suggested dose

New suggested dose: mg
 Every Hour: hour
 Expected Cmax: mcg/ml
 Expected Cmin: mcg/ml

To calculate new dose and Cmin based on desired Cmax

Desired Cmax: mcg/ml
 Desired Interval: hour
 Expected Cmin: mcg/ml
 New Dose: mg

Figure 3.0-4 Pre & Post Calculator for Vancomycin Detail

4.0 Calculator for Gentamicin / Amikacin

To view the calculator in the system, perform the steps below:

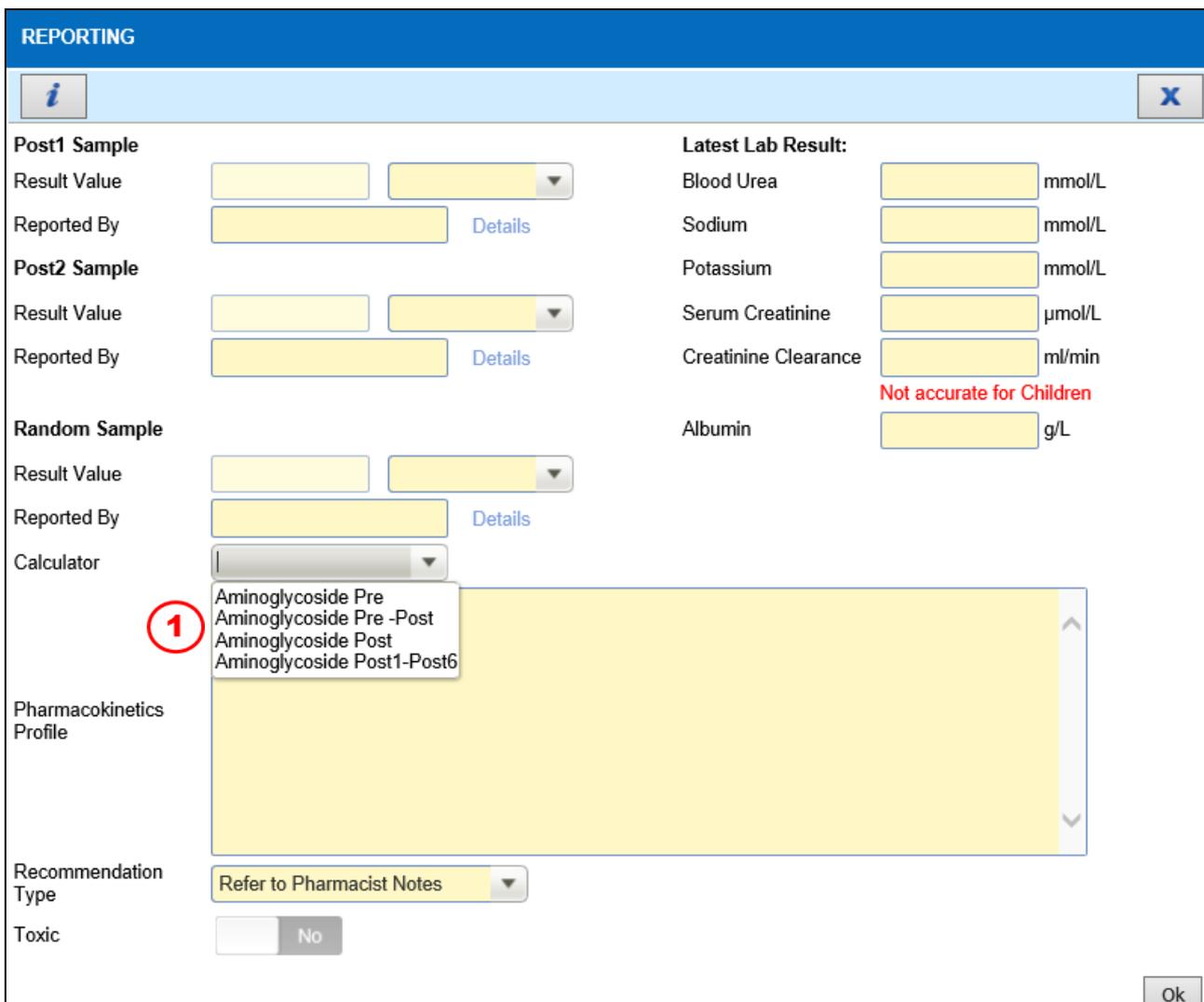


Figure 4.0-1 Calculator for Gentamicin / Amikacin

STEP 1

Select calculator from dropdown box example:

- Aminoglycoside Pre-Post
- Aminoglycoside Pre
- Aminoglycoside Post
- Aminoglycoside Post 1- Post 6

Note

- Choose the calculator and fill the result in the box as per Figure 4.0-1. Sample Details consists of:
 - Calculator for **Aminoglycoside Pre**
 - a. Dose medication in mg
 - b. Interval Time
 - c. Pre level result concentration
 - d. Serum Creatinine
 - e. CrCl - user can select option to use Adult or paediatrics calculation
 - f. Vd

g. Expected Cmax if assuming Expected Cmin = Pre level result

➤ Calculator for **Aminoglycoside Post**

- Dose medication in mg
- Interval Time
- Sampling Time: Duration between time after complete infusion and post sampling(t)
- Serum level: Post level result concentration, Serum Creatinine
- Vd, Ke and t1/2
- CrCl - user can select option to use Adult or paediatrics calculation
- New dose and Cmin based on desired Cmax
- Expected Cmin and Cmax based on new suggested dose

➤ Calculator for **Aminoglycoside Pre & Post**

- Dose medication in mg
- Interval Time
- Sampling Time: Duration between time after complete infusion and post sampling(t) and Duration between completion time of pre and post sampling (t2-t1)
- Serum level: Pre level result concentration, Post level result concentration and Serum Creatinine
- CrCl - user can select option to use Adult or paediatrics calculation
- Vd, Ke and t1/2
- New dose and Cmin based on desired Cmax
- Expected Cmin and Cmax based on new suggested dose

➤ Calculator for **Aminoglycoside Post 1 & Post 6**

- Dose medication in mg
- Interval Time
- Sampling Time: Duration between time after complete infusion and post 1 sampling(t) and Duration between completion time of post6 and post1 sampling (t6-t1)
- Serum level: Post 1 level result concentration, Post 6 level result concentration and Serum Creatinine
- CrCl - user can select option to use Adult or paediatrics calculation
- Vd, Ke and t1/2
- New dose and Cmin based on desired Cmax
- Expected Cmin and Cmax based on new suggested dose

- After save the calculator, the result will appear at the Pharmacokinetics Profile box
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 4.0-1

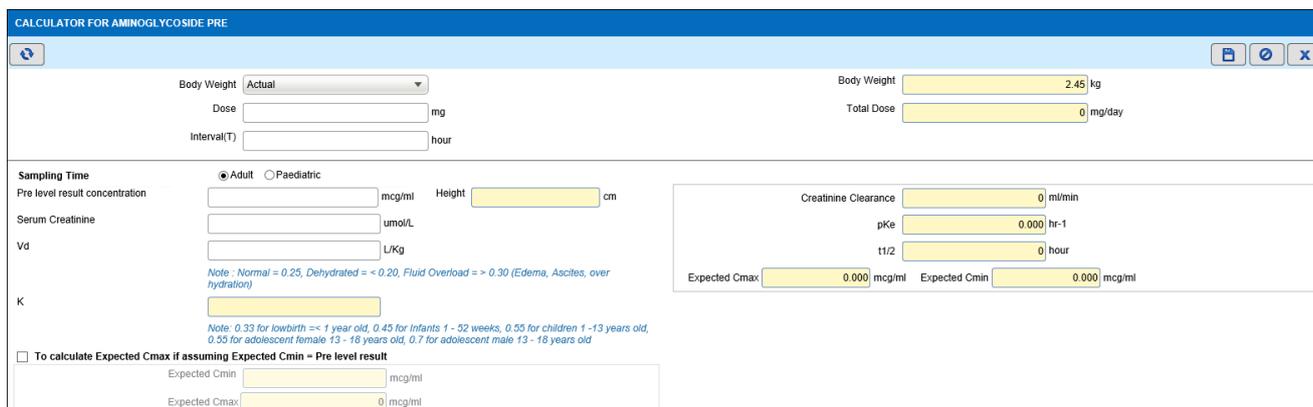


Figure 4.0-2 Calculator for Aminoglycoside Pre Detail (Adult)

Note

- If adult selected shall display calculation as shown in Figure 4.0-2
- K and height field will disable

- Serum Creatinine, system shall auto calculate value for Creatinine Clearance (CrCl)
- If paediatric selected shall display calculation as shown in Figure 4.0-3

Figure 4.0-3 Calculator for Aminoglycoside Pre Detail (Paediatric)

Note

- User able to insert value for:
 - Serum Creatinine
 - K
 - Height
- System shall auto calculate value for Creatinine Clearance (CrCl) based on formula as below:

$$\text{Formula} = \frac{K \times \text{height (cm)}}{88.4} \text{ Secr } \mu\text{mol/L}$$

Figure 4.0-4 Calculator for Aminoglycoside Post Detail

CALCULATOR FOR AMINOGLYCOSIDE PRE POST

Body Weight: | Body Weight: kg
 Dose: mg | Total Dose: mg/daily
 Interval(T): hour

Sampling Time
 Duration between time after complete infusion and post sampling(t): hour
 Duration between completion time of pre and post sampling (t2-t1): hour

Serum level
 Pre level result concentration: mcg/ml
 Post level result concentration: mcg/ml
 Serum Creatinine: umol/L
 Creatinine Clearance: ml/min

To calculate Vd, Ke and t1/2
 Ke: hr⁻¹
 Vd: L/Kg
 t1/2: hour
 Cmax: mcg/ml | Cmin: mcg/ml

To calculate new dose and Cmin based on desired Cmax
 Desired Cmax: mcg/ml
 Desired Interval: hour
 Expected Cmin: mcg/ml | New Dose: mg

To calculate Paediatric Creatinine Clearance
 Height: cm | Creatinine Clearance: ml/min/1.73m²
 K:
Note: 0.33 for lowbirth =<= 1 year old, 0.45 for Infants 1 - 52 weeks, 0.55 for children 1 -13 years old, 0.55 for adolescent female 13 - 18 years old, 0.7 for adolescent male 13 - 18 years old

To calculate expected Cmin and Cmax based on new suggested dose
 New suggested dose: mg
 New Interval: hour
 Expected Cmax: mcg/ml | Expected Cmin: mcg/ml

Note : mcg/ml = mg/L

Figure 4.0-5 Calculator for Aminoglycoside Pre Post Detail

CALCULATOR FOR AMINOGLYCOSIDE POST1 & POST6

Body Weight: | Body Weight: kg
 Dose: mg | Total Dose: mg/daily
 Interval(T): hour

Sampling Time
 Duration between time after complete infusion and post1 sampling(t): hour
 Duration between completion time of post6 and post1 sampling (t6-t1): hour

Serum level
 Post 1 level result concentration: mcg/ml
 Post 6 level result concentration: mcg/ml
 Serum Creatinine: umol/L
 Creatinine Clearance: ml/min

To calculate Vd, Ke and t1/2
 Ke: hr⁻¹
 Vd: L/Kg
 t1/2: hour
 Cmax: mcg/ml | Cmin: mcg/ml

To calculate new dose and Cmin based on desired Cmax
 Desired Cmax: mcg/ml
 Desired Interval: hour
 Expected Cmin: mcg/ml | New Dose: mg

To calculate Paediatric Creatinine Clearance
 Height: cm | Creatinine Clearance: ml/min/1.73m²
 K:
Note: 0.33 for lowbirth =<= 1 year old, 0.45 for Infants 1 - 52 weeks, 0.55 for children 1 -13 years old, 0.55 for adolescent female 13 - 18 years old, 0.7 for adolescent male 13 - 18 years old


To calculate expected Cmin and Cmax based on new suggested dose
 New suggested dose: mg
 New Interval: hour
 Expected Cmax: mcg/ml | Expected Cmin: mcg/ml

Note : mcg/ml = mg/L

Figure 4.0-6 Calculator for Aminoglycoside Pre 1-Post 6 Detail

5.0 Calculator for Valproic Acid

To view the calculator in the system, perform the steps below:

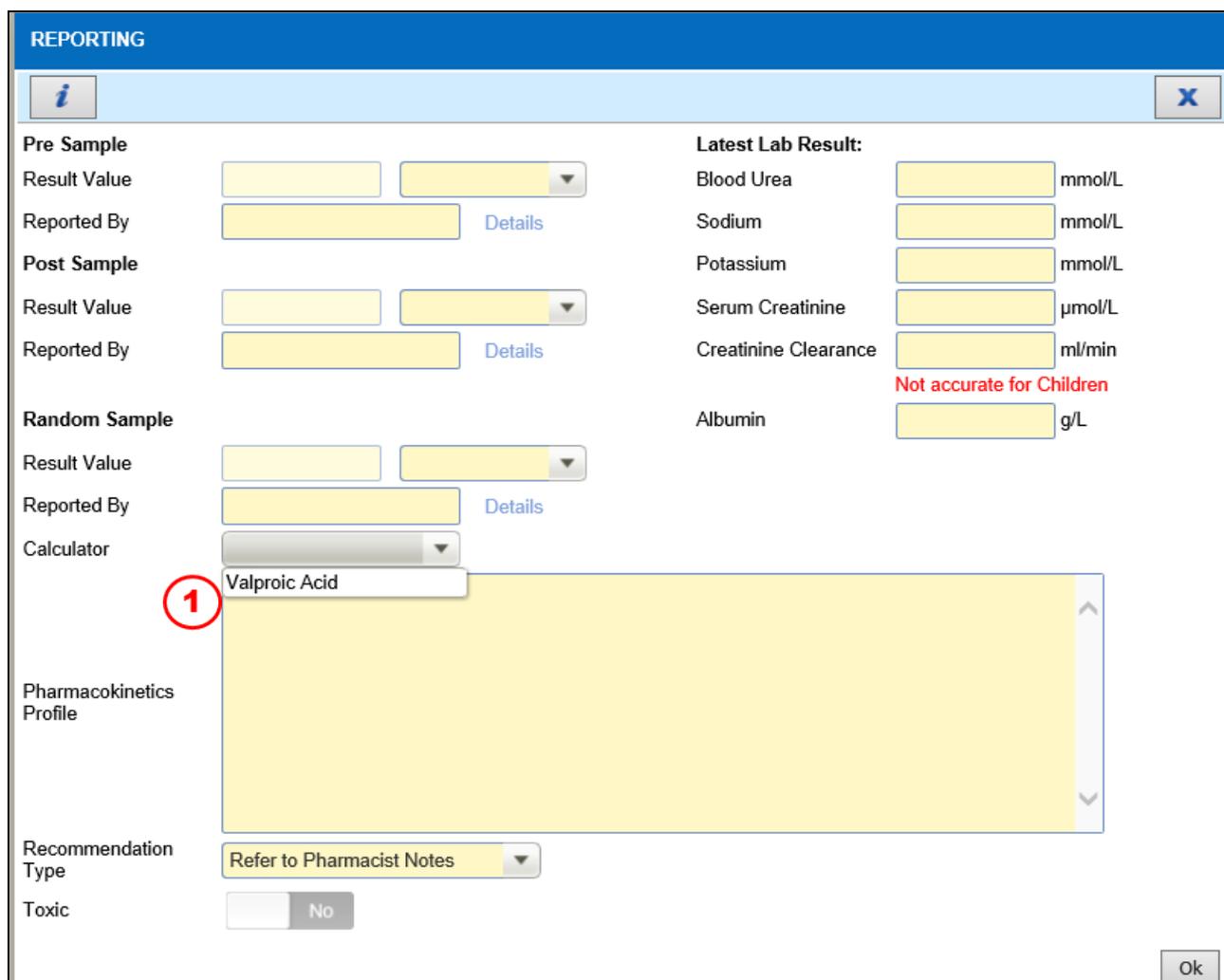

REPORTING

Pre Sample
 Result Value: [] []
 Reported By: [] Details

Post Sample
 Result Value: [] []
 Reported By: [] Details

Random Sample
 Result Value: [] []
 Reported By: [] Details
 Calculator: [Valproic Acid]

Latest Lab Result:
 Blood Urea: [] mmol/L
 Sodium: [] mmol/L
 Potassium: [] mmol/L
 Serum Creatinine: [] μmol/L
 Creatinine Clearance: [] ml/min
 Albumin: [] g/L
Not accurate for Children

Pharmacokinetics Profile

Recommendation Type: [Refer to Pharmacist Notes]
 Toxic: [No]

Ok

Figure 5.0-1 Calculator for Valproic Acid

STEP 1

Select calculator from dropdown box example:

- Valproic Acid

Note


- Choose the calculator and fill the result in the box as per Figure 5.0-1. Sample Details consists of:
 - **Valproic Acid** calculator
 - Dose medication in mg
 - Interval (T) in Hours
 - Total Dose in mg/day
 - Choose between monotherapy or polytherapy
 - ✓ **Monotherapy**- single or combination therapies of different enzyme activities
 - ✓ **Polytherapy** - combination therapies of similar enzyme activities
 - Cp Measured in mcg/ml
 - Constant CI in mg/kg/hr
 - Vd , Ke and t1/2 based on CI
 - New Maintenance Dose (NMD) based on I, Cp measured and Cp desired
 - Cpss level based on current dose and CI

- **Expected Cp based on CI and New Suggested Dose (NSD)**
- After save the calculator, the result will appear at the Pharmacokinetics Profile box
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 5.0-1

Figure 5.0-2 Calculator for Valproic Acid Detail

6.0 Calculator for Aminophylline / Theophylline

To view the calculator in the system, perform the steps below:

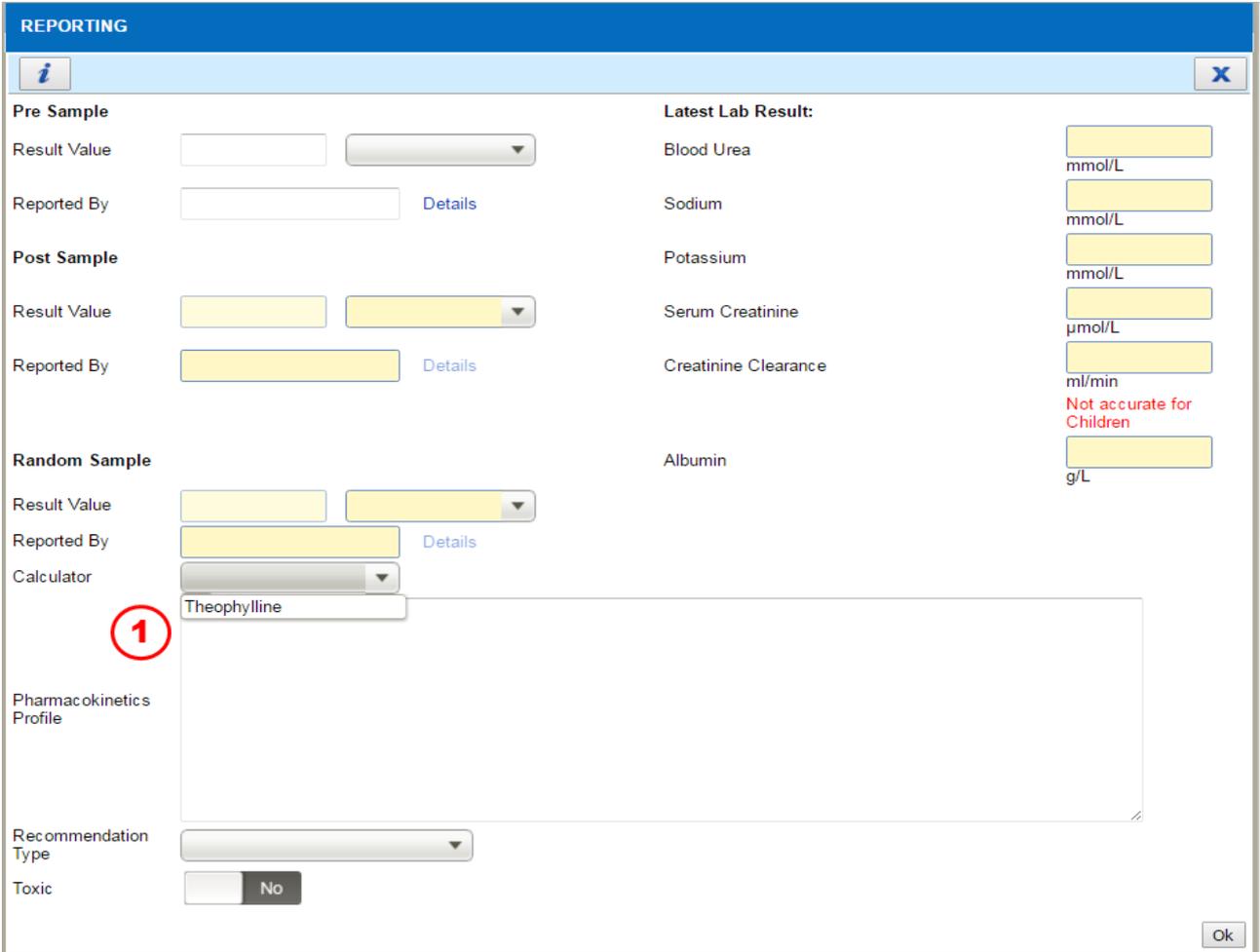

The screenshot shows a software window titled "REPORTING" with a blue header. It contains several input fields and sections:

- Pre Sample:** Includes "Result Value" (text input), "Reported By" (text input with "Details" link), and "Latest Lab Result:" section with "Blood Urea" (mmol/L), "Sodium" (mmol/L), "Potassium" (mmol/L), "Serum Creatinine" (µmol/L), and "Creatinine Clearance" (ml/min, with a note "Not accurate for Children").
- Post Sample:** Includes "Result Value" (text input), "Reported By" (text input with "Details" link), and "Albumin" (g/L).
- Random Sample:** Includes "Result Value" (text input), "Reported By" (text input with "Details" link), and "Calculator" (dropdown menu). A red circle with the number "1" highlights the "Theophylline" option in this dropdown.
- Pharmacokinetics Profile:** A large empty text area.
- Recommendation Type:** A dropdown menu.
- Toxic:** A toggle switch currently set to "No".

Figure 6.0-1 Calculator for Aminophylline / Theophylline

STEP 1

Select calculator from dropdown box example:

- Aminophylline / Theophylline

Note

- Choose the calculator and fill the result in the box as per Figure 6.0-1. Sample Details consists of:

- **Theophylline** calculator
- Choose between oral or iv

Oral

- Dose medication in mg/hr
- Cp measured in mcg/ml
- Clearance based on patient's factor(s)
- Ke and t1/2 based on Cl
- Loading Dose (LD) if NO Theophylline given within 24 hours
- Loading Dose (LD) if NO Theophylline given within 24 hours
- Without Therapy Period (T) if level is SUPRA-therapeutic based on Cp desired and Ke
- IV Aminophylline to Oral
- Cpss level based on current dose and Cl

- j. Loading Dose (LD) if Theophylline given within 24 hours
- k. Incremental Loading Dose (ILD) if level is SUB-therapeutic based on Cp desired and Vd
- l. Infusion rate (Ko) based on Cp desired and Cl

Intravenous (iv)

- a. Dose medication in mg
- b. Interval (T) in hours
- c. Cp measured in mcg/ml
- d. Clearance based on patient's factor(s)
- e. Ke and t1/2 based on Cl
- f. Loading Dose (LD) based on Cp desired and Vd
- g. Incremental Loading Dose (ILD) if level is SUB-therapeutic based on Cp desired and Vd
- h. New Maintenance Dose (IMD) based on Cp desired and Cl
- i. Cpss level based on current dose and Cl
- j. New Dose (Ko) based on Cp and Cl when convert Oral to IV Aminophylline

- After save the calculator, the result will appear at the Pharmacokinetics Profile box.
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 6.0-1

Figure 6.0-2 Calculator for Theophylline for Drug Aminophylline

7.0 Calculator for Carbamazepine

To view the calculator in the system, perform the steps below:

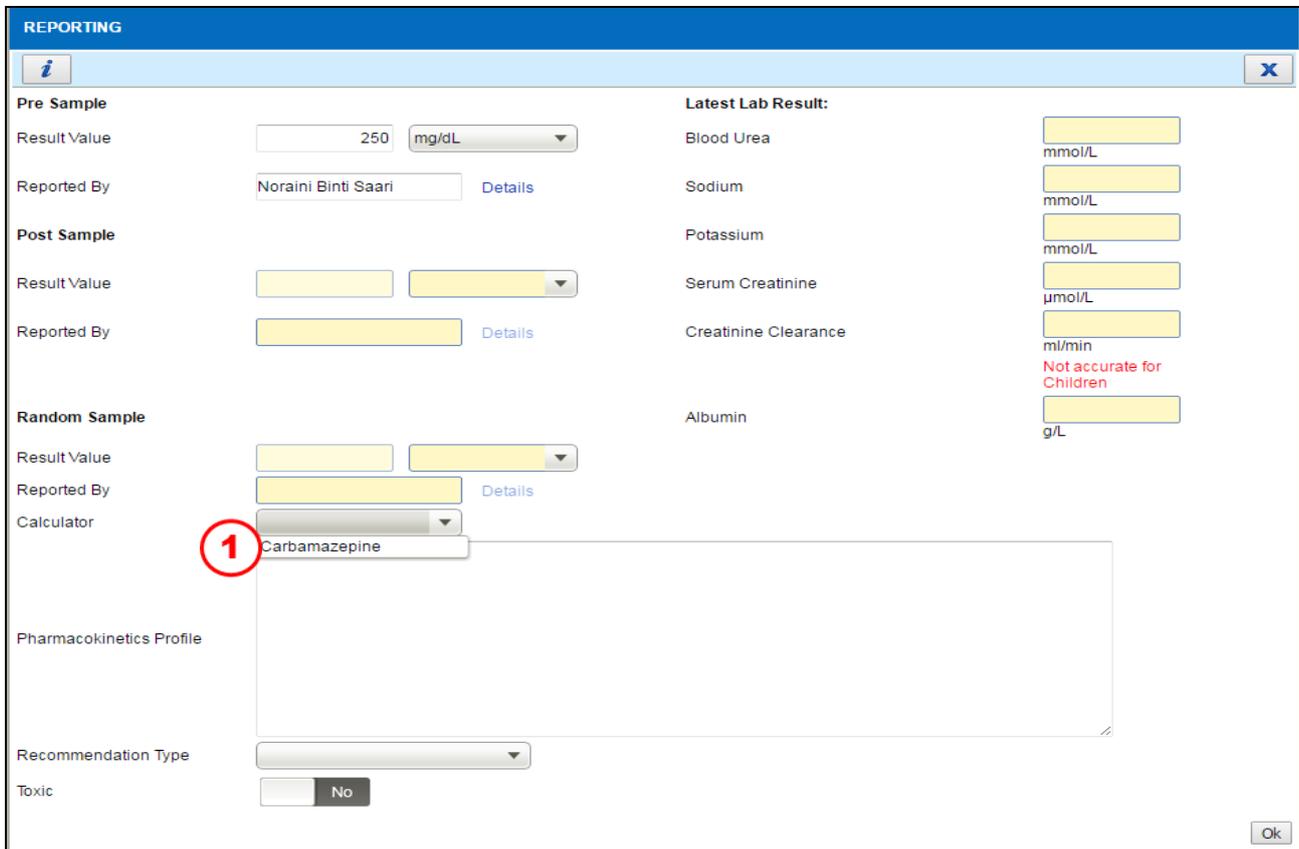


Figure 7.0-1 Calculator for Carbamazepine

STEP 1

Select calculator from dropdown box example:

- Carbamazepine

Note

- Choose the calculator and fill the result in the box as per Figure 7.0-1. Sample Details consists of:
 - **Carbamazepine** calculator
 - Dose medication in mg
 - Interval Time in Hours
 - Total Dose in mg/day
 - Choose between monotherapy or polytherapy
 - ✓ **Monotherapy**- single or combination therapies of different enzyme activities
 - ✓ **Polytherapy** - combination therapies of similar enzyme activities
 - Cp Measured in mcg/ml
 - Vd , Ke and t1/2
 - New Maintenance Dose (NMD) based on Cp measured and Cp desired
 - Cpss level based on current dose and Cl
 - Expected Cp based on New Suggested Dose

- After save the calculator, the result will appear at the Pharmacokinetics Profile box.
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 7.0-1


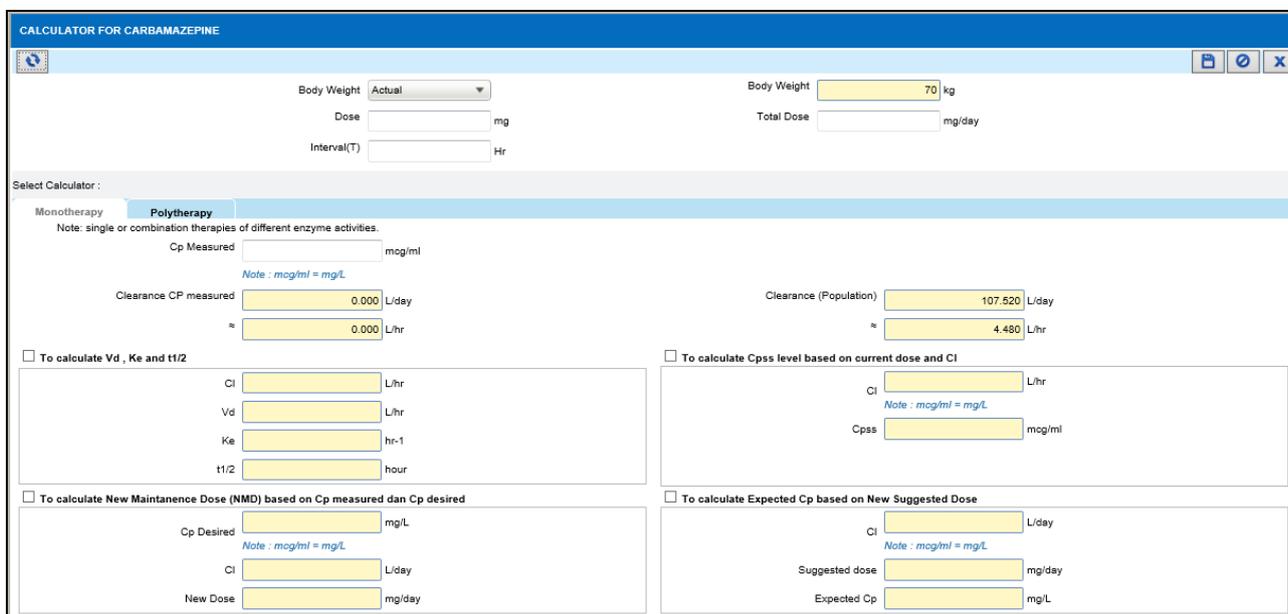


Figure 7.0-2 Calculator for Carbamazepine for Drug Carbamazepine

8.0 Calculator for Digoxin

To view the calculator in the system, perform the steps below:

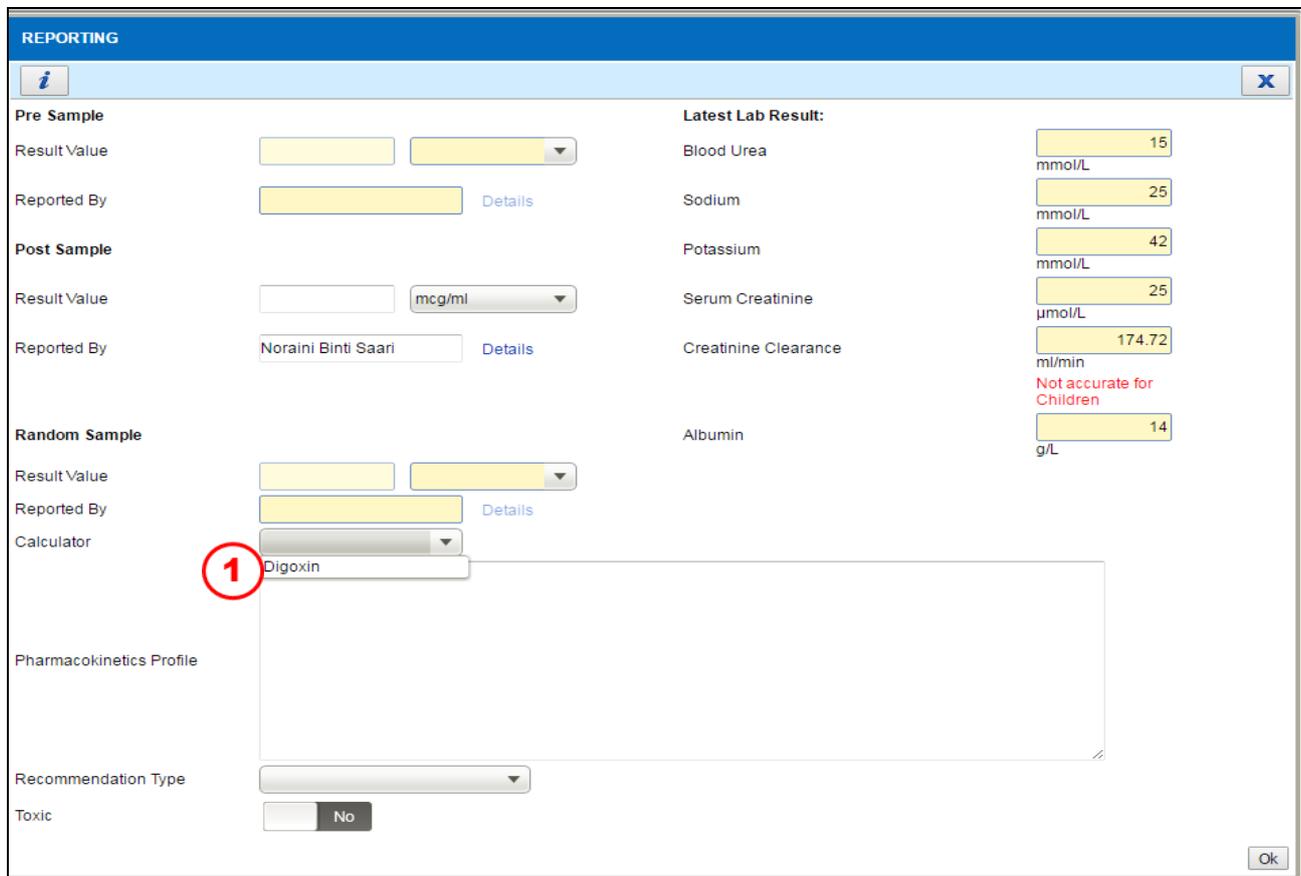

The screenshot shows a 'REPORTING' window with the following sections:

- Pre Sample:** Result Value (input), Reported By (input with 'Details' link).
- Post Sample:** Result Value (input, unit 'mcg/ml'), Reported By (input 'Noraini Binti Saari' with 'Details' link).
- Random Sample:** Result Value (input), Reported By (input with 'Details' link), Calculator (dropdown menu with '1' circled around 'Digoxin').
- Pharmacokinetics Profile:** (empty text area).
- Recommendation Type:** (dropdown menu).
- Toxic:** (checkbox 'No').
- Latest Lab Result:**
 - Blood Urea: 15 mmol/L
 - Sodium: 25 mmol/L
 - Potassium: 42 mmol/L
 - Serum Creatinine: 25 µmol/L
 - Creatinine Clearance: 174.72 ml/min (Note: Not accurate for Children)
 - Albumin: 14 g/L

Figure 8.0-1 Calculator for Digoxin

STEP 1

Select calculator from dropdown box example:

- Digoxin

Note

- Choose the calculator and fill the result in the box as per Figure 8.0-1. Sample Details consists of:
 - **Digoxin calculator**
 - Total Dose in mg/day
 - Serum Creatinine in umol/L
 - CrCl - user can select option to use Adult or paediatrics calculation
 - Choose patient condition:
 - ✓ **Without CHF and not renally impaired**
 - ✓ **Without CHF and renally impaired**
 - ✓ **With CHF and not renally impaired**
 - ✓ **With CHF and renally impaired**
 - Cp Measured in mcg/ml
 - Ke and t1/2 based on Cl
 - Vd based on Patient's Factor(s)
 - Cpss level based on current dose and Cl
 - Expected Cp based on Cl and New Suggested Dose
 - Loading Dose(LD) based on Vd and Cp desired
 - New Maintenance Dose(NMD) based on Cl and Cp desired

- After save the calculator, the result will appear at the Pharmacokinetics Profile box.
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 8.0-1

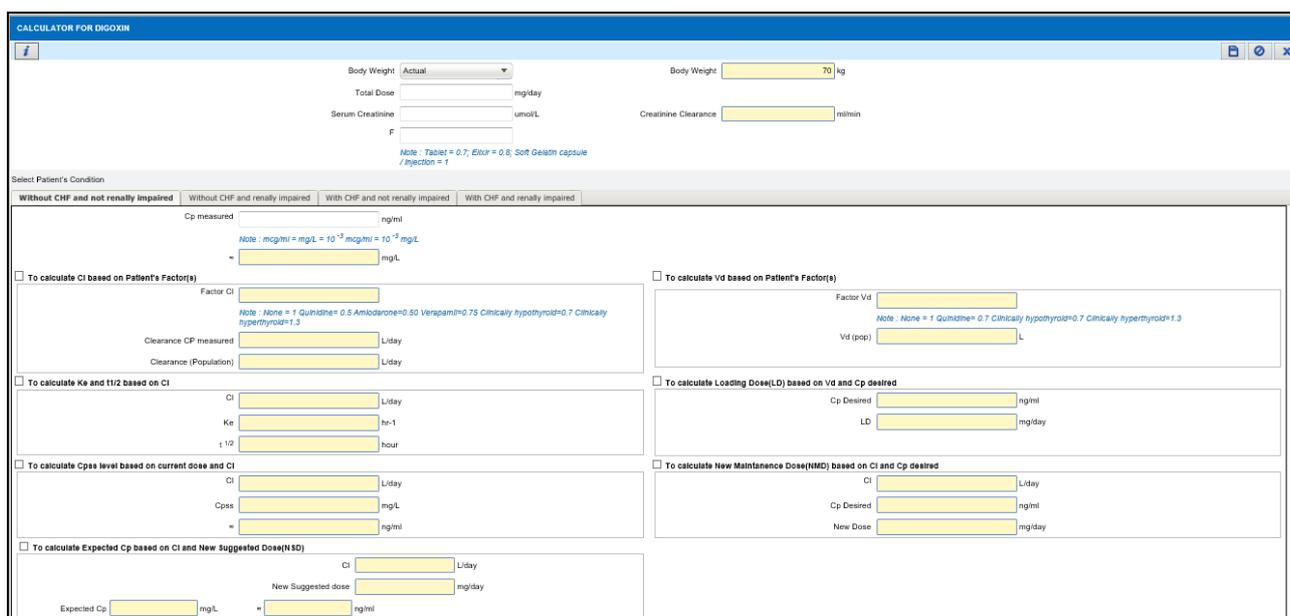


Figure 8.0-2 Calculator for Digoxin for Drug Digoxin

9.0 Calculator for Phenobarbitone

To view the calculator in the system, perform the steps below:

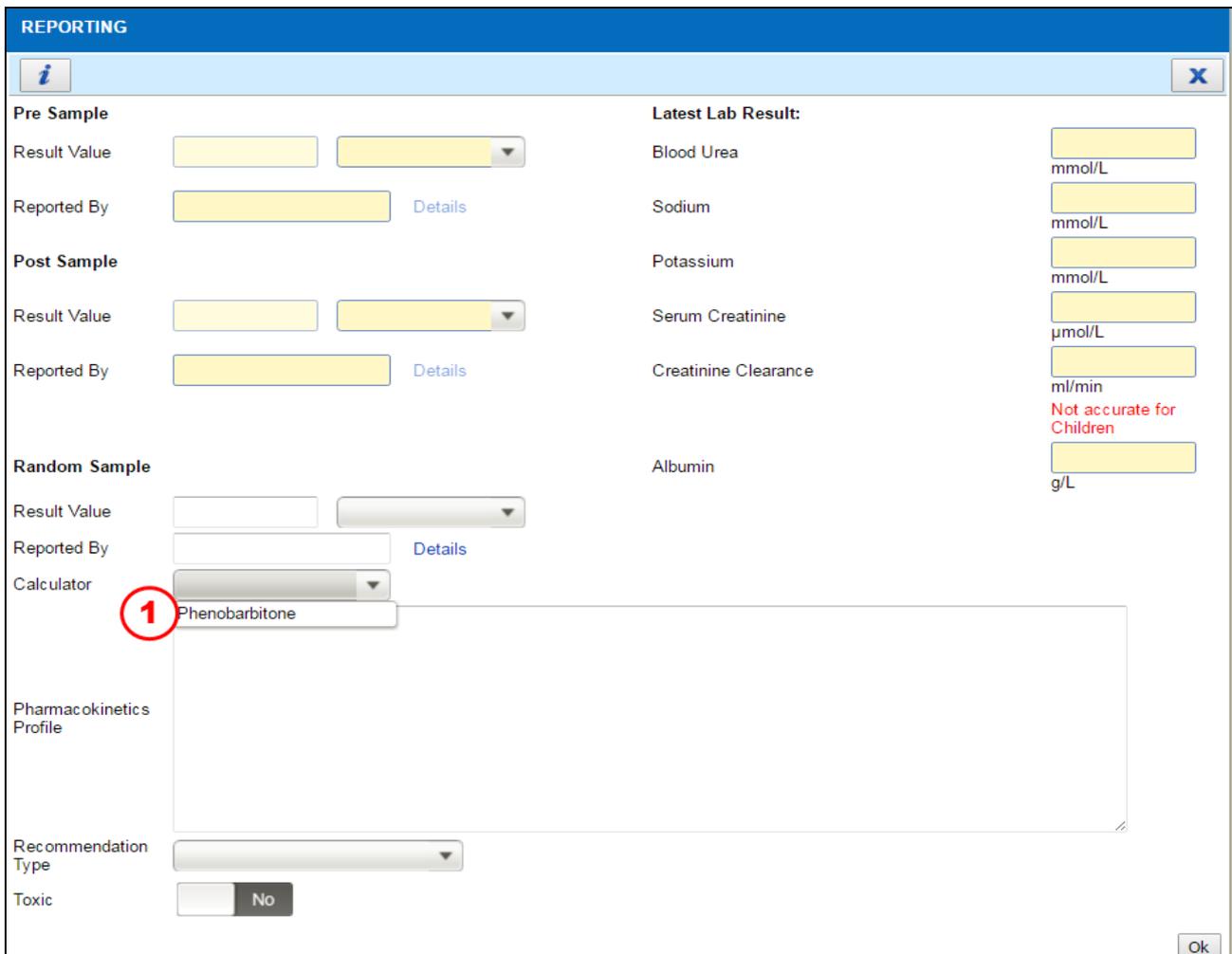


Figure 9.0-1 Calculator for Phenobarbitone

STEP 1

Select calculator from dropdown box example:

- Phenobarbitone

Note

- Choose the calculator and fill the result in the box as per Figure 9.0-1. Sample Details consists of:
 - **Phenobarbitone calculator**
 - Dose medication in mg
 - Interval Time in Hours
 - Total Dose in mg/day
 - Cp Measured in mcg/ml
 - Constant Cl in ml/kg/hr
 - Vd, Ke and t1/2 based on Cl
 - New Maintenance Dose (NMD) based on Cl, Cp measured and Cp desired
 - Cpss level based on current dose and Cl
 - Expected Cp based on Cl and New Suggested Dose (NSD)

- After save the calculator, the result will appear at the Pharmacokinetics Profile box.
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 9.0-1

CALCULATOR FOR PHENOBARBITONE

Body Weight: Actual (dropdown) | Interval(T): hour

Body Weight: kg | S:

Dose: mg | Total Dose: mg/day

Note: Oral / IV = 1

Cp measured: mcg/ml | Constant CI: ml/kg/hr

Note: mcg/ml = mg/L | Note: Adult = 4 ml/kg/hr Children (≤ 12years 11month) = 8 ml/kg/hr Neonate = 4 ml/kg/hr

Clearance CP measured: L/day | Clearance (Population): L/day

≈ L/hr | ≈ L/hr

To calculate Vd, Ke and t1/2 based on CI

CI: L/hr | Constant Vd: L/kg

Note: Adult = 0.6 L/kg, Neonate = 0.8 - 1 L/kg

Vd: L | Ke: hr⁻¹ | t 1/2: hour

To calculate Cpss level based on current dose and CI

CI: L/hr | Cpss: mcg/ml

Note: mcg/ml = mg/L

To calculate New Maintenance Dose (NMD) based on CI, Cp measured and Cp desired

Cp desired: mcg/ml | CI: L/day

Note: mcg/ml = mg/L

New Dose: mg/day

To calculate Expected Cp based on CI and New Suggested Dose(NSD)

CI: L/day | NSD: mg/day

Expected Cp: mcg/ml

Note: mcg/ml = mg/L

Figure 9.0-2 Calculator for Phenobarbitone for Drug Phenobarbitone

10.0 Calculator for Phenytoin

To view the calculator in the system, perform the steps below:

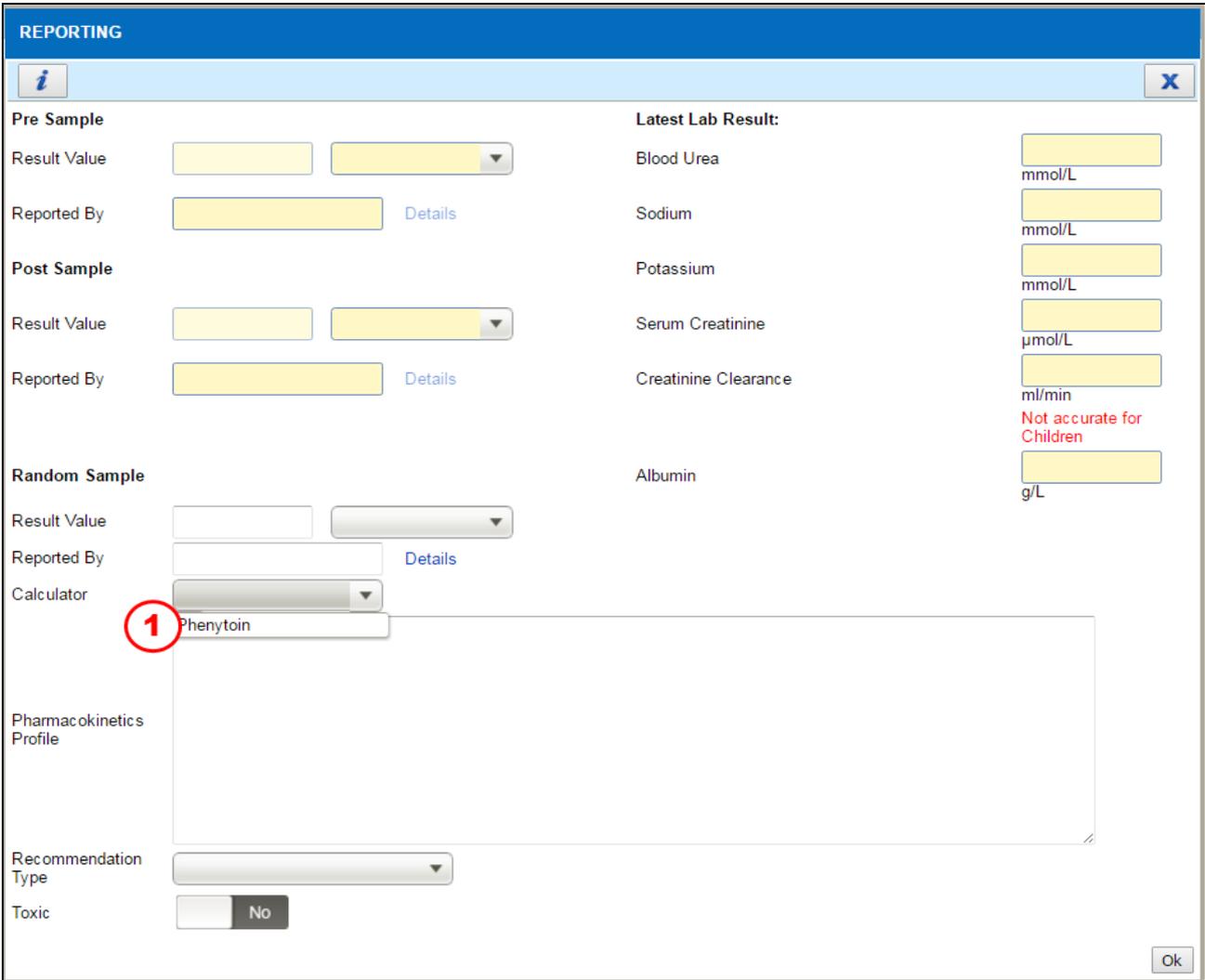


Figure 10.0-1 Calculator for Phenytoin

STEP 1

Select calculator from dropdown box example:

- Phenytoin

Note

- Choose the calculator and fill the result in the box as per Figure 10.0-1. Sample Details consists of:

➤ Phenytoin calculator

- Total Dose in mg/day
- Albumin g/dl
- Choose patient condition between **Without Hypoalbuminemia, Hypoalbuminemia or Hypoalbuminemia with ESRF Without Hypoalbuminemia**
 - a. Cp Measured in mcg/ml
 - b. Vmax (population)
 - c. Vd, Ke, t1/2 and Cl based on Vmax
 - d. Cpss level based on Vmax
 - e. Incremental Loading Dose (ILD) based on Cp desired and Cp measured

- f. Loading Dose (LD) based on C_p desired
- g. therapy withhold period (T) when level is toxic
- h. New Maintenance Dose (NMD) based on V_{max} and C_p desired

Hypoalbuminemia/Hypoalbuminemia with ESRF

- a. C_p Measured in mcg/ml
- b. Corrected C_p Measured in mcg/ml
- c. V_{max} (population)
- d. V_d , K_e , $t_{1/2}$ and Cl based on V_{max}
- e. C_{pss} level based on V_{max}
- f. Incremental Loading Dose (ILD) based on C_p desired and C_p measured
- g. Loading Dose (LD) based on C_p desired
- h. therapy withhold period (T) when level is toxic
- i. New Maintenance Dose (NMD) based on V_{max} and C_p desired

- After save the calculator, the result will appear at the Pharmacokinetics Profile box.
- Calculator is optional only user can direct type at the Pharmacokinetics Profile box without using the calculator as per Figure 10.0-1

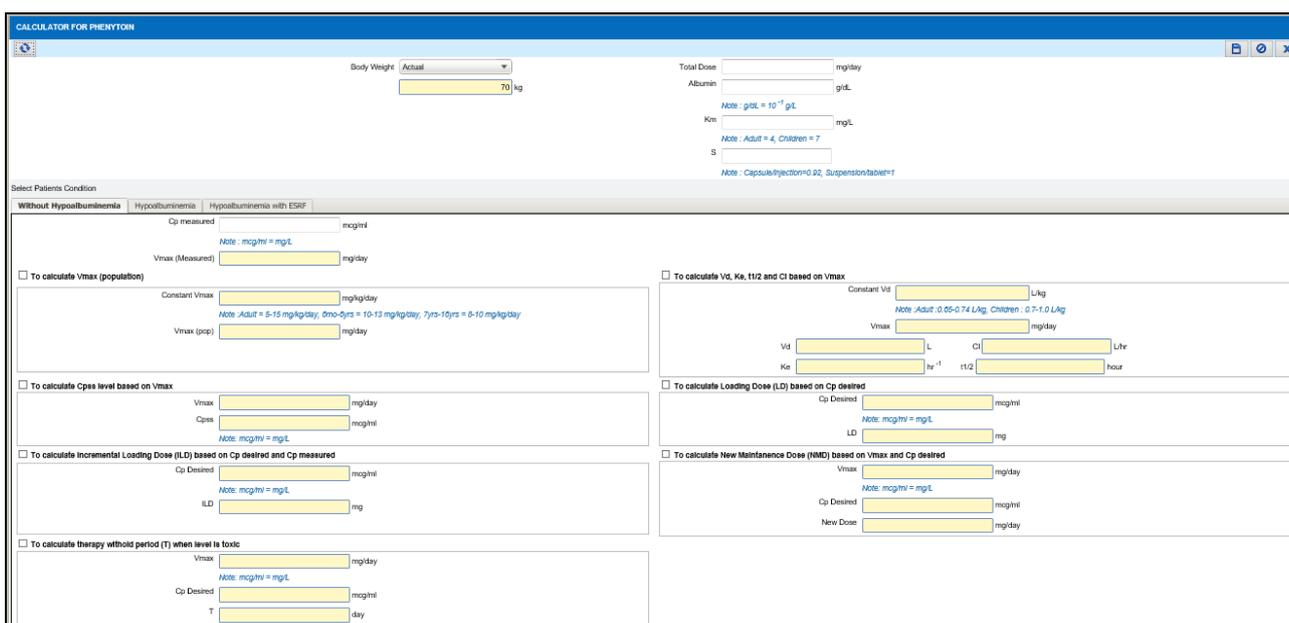


Figure 10.0-2 Calculator for Phenytoin for Drug Phenytoin

11.0 Acronyms

Abbreviation	Definition
PhIS	Pharmacy Information System
CPS	Clinical Pharmacy System
PM	Patient Management
eGL	Electronic Guarantee Letter
HRMIS	Human Resource Management Informasi System
MRN	Medical Record Number
MOH	Ministry Of Health

12.0 Links To Clinical Modules

No	Module	PDF Links	No	Module	PDF Links
1	<i>Inpatient</i>	Click Here	12	<i>CDR Dispensing</i>	Click Here
2	<i>CDR Order</i>	Click Here	13	<i>Methadone Dispensing</i>	Click Here
3	<i>TDM Order</i>	Click Here	14	<i>PN Dispensing</i>	Click Here
4	<i>PN Order</i>	Click Here	15	<i>Order Management</i>	Click Here
5	<i>IV Order</i>	Click Here	16	<i>Patient Management</i>	Click Here
6	<i>Prepacking</i>	Click Here	17	<i>Radiopharmaceuticals</i>	Click Here
7	<i>Galenical</i>	Click Here	18	<i>Outpatient</i>	Click Here
8	<i>MTAC</i>	Click Here	19	<i>Special Drug Request</i>	Click Here
9	<i>ADR & DAC</i>	Click Here	20	<i>MAR</i>	Click Here
10	<i>Medication Counselling</i>	Click Here	21	<i>DICE</i>	Click Here
11	<i>Ward Pharmacy</i>	Click Here	22		