
Pharmacy Information System (PhIS) and Clinic Pharmacy System (CPS)

PhIS & CPS THP & Zebra Configuration Manual

Version : 1.0

Document ID : PhIS_CPS/THP_&_Zebra_Configuration_Manual_v1.0



Table of Contents

Revision History	iv
1. THP Printer Configuration.....	5
2. Zebra Printer Configurations	8



© 2011-2023 Pharmaniaga Logistics Sdn Bhd (PLSB)

CONFIDENTIAL COPYRIGHTED MATERIAL – *The information includes all concepts, comments, recommendations, and material, contained here in shall remain the property of Pharmaniaga Logistics Sdn Bhd (PLSB). No portion of this document shall be disclosed, duplicated or used in whole or in part of any purpose other than the purpose of the Pharmacy Information System & Clinic Pharmacy System (PhIS & CPS) Project execution only.*



Revision History

Version No.	Date of Release	Consolidated by	Verified by	List of changes from Previous Version
1.0	10-8-2023	Arif Syahrir Mohd Isah	Sophian Ezudean Mustapha	Initial document

Verified and Acknowledge by

Name	Role / Designation	Organization	Signature & Date
Sophian Ezudean	Manager	Pharmaniaga Logistic S/B	
		Program Perkhidmatan Farmasi (PPF), KKM	
		Bahagian Pengurusan Maklumat (BPM), KKM	


Peer Review

Version No.	Reviewed / Verified by	List of changes from Previous Version
1.0	(QA) Juridah	- Overall review


1. THP Printer Configuration

Printing Configuration

- Right-click on the printer icon and select the **Printing Preferences**.

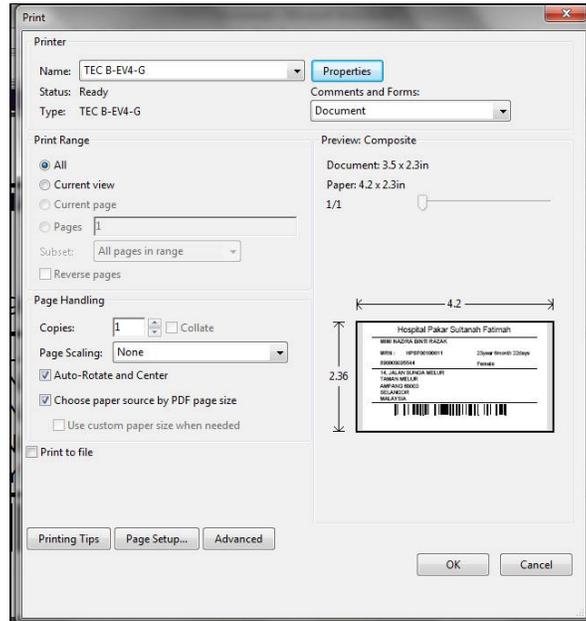

- Click on the **Page Setup** tab to allow specifying the default media size, orientation, and other properties of the media.

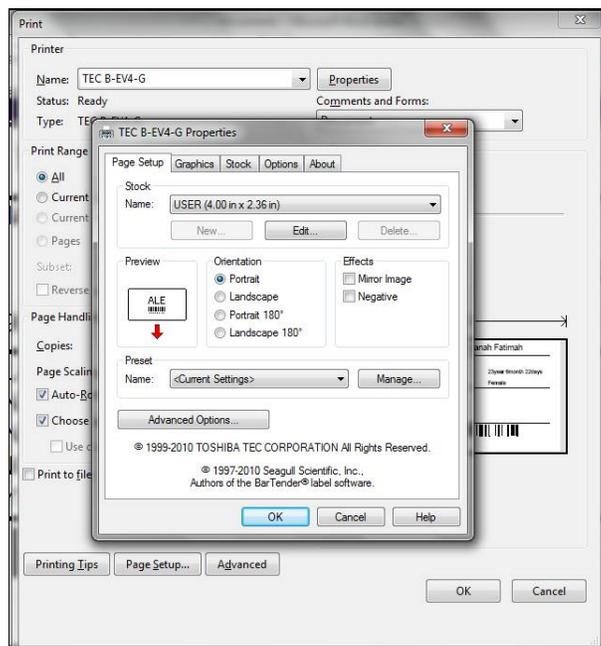


- **Name**

Allow specifying a previously loaded stock by name to use, modify, or delete.


- **New, Edit, Delete**

Allow adding a new label stock specification to the list, or editing / deleting a label stock already on the list.

- Set Width (90mm) and Height (60mm) on the Page Setup Tab. Then Click OK

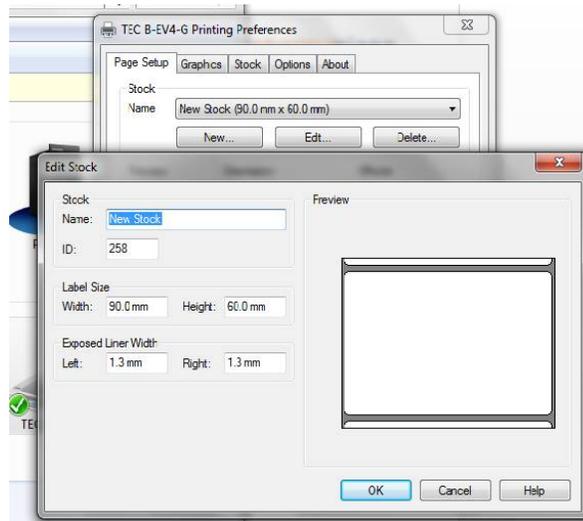


- Click on the Stock tab to allow specifying the default Print Method, Sensor, Label Gap on Media Settings, and Print Speed on Issue Settings.

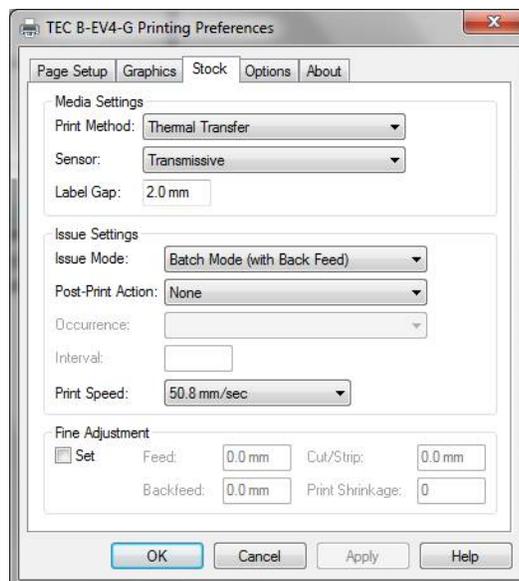



- On the Media Settings, set print Method to “Thermal Transfer” and Sensor to “Transmissive”. Select the “50.8mm/sec” for Print Speed, then Click OK to save configuration.

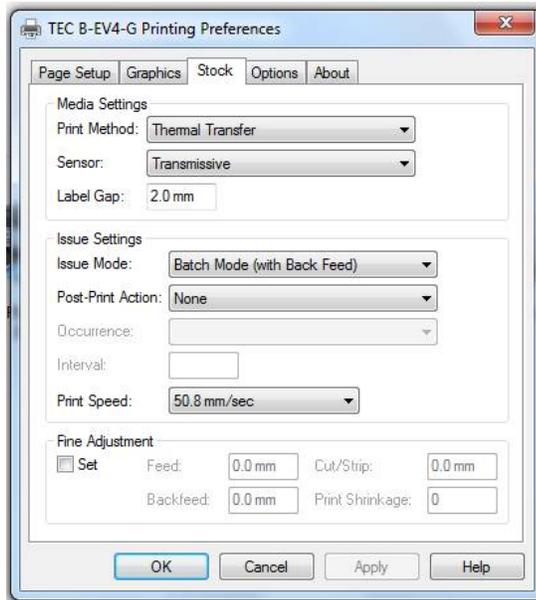



2. Zebra Printer Configurations

A. Printer Configuration

Printing Preference

- Right-click on the printer icon and select the **Printing Preferences**

- **Option >**

Setting:

No of Copies: **1**

Speed: **101.6 mm/s**

Darkness: **15**

Stock: **User defined**

Paper Format:

Select **mm** and **portrait**

Size:

Width: **93.70**

Height: **63.00**


Unprintable area:

Left: **1.00**

Right: **0.00**


Top: **0.00**

Bottom: **0.00**

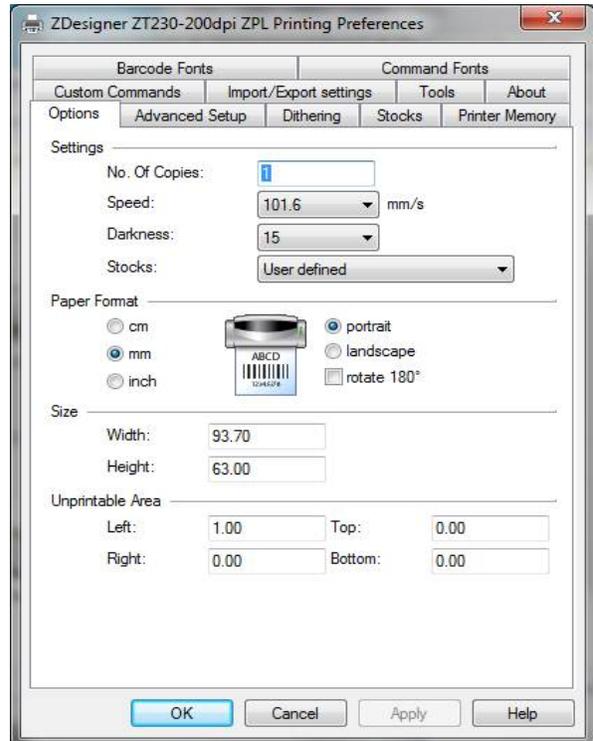




- **Advance Setup >**
Operation Mode:
Select ***Tear off***
Media Type:
Select ***Thermal Transfer***
Tracking Mode:
Select ***Continuous***
Adjustment:
Top: 0.00 mm
Tear Off: 0.00 mm
Left Position: 0.00 mm

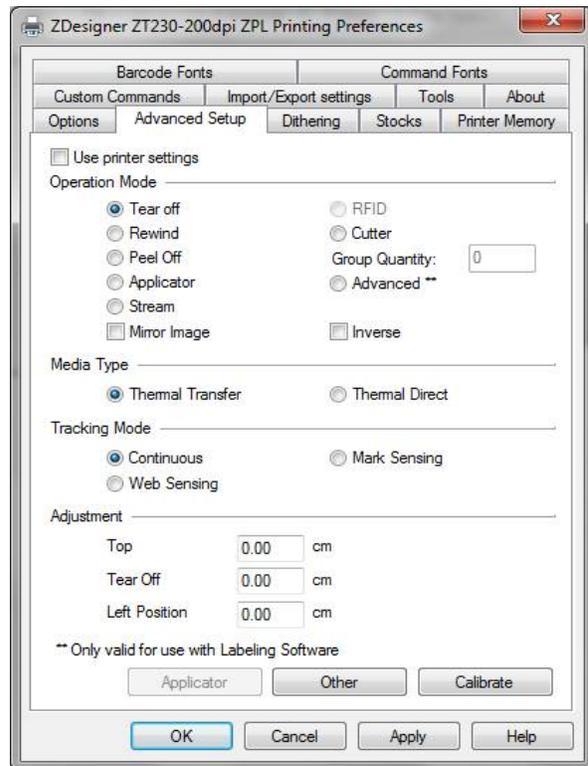



- **Dithering >**
Dithering type:
Select ***Smooth***

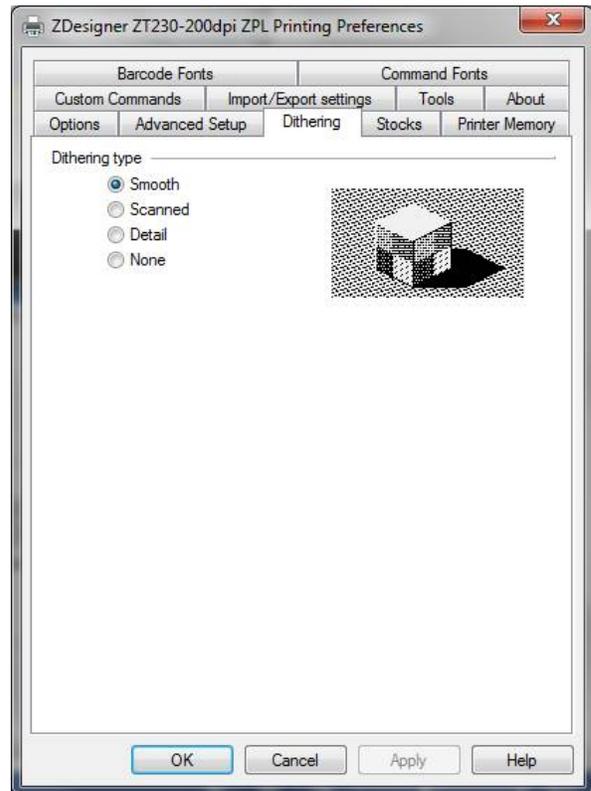



- **Barcode Fonts >**
Defined barcode fonts:
Select **Code128**

